

PHYSICAL THERAPY PRODUCT CATALOG

- Orthopedics
- Mobilization
- Examination
- Treatment
- Functional Training
- Vestibular Rehabilitation
- Sports Rehabilitation

Translating evidence into practice via ergonomically efficient products.



 **HENRY SCHEIN[®]**
www.henryschein.com

made in the USA
with US & imported parts



ISO 13485:2003
FM 584901



TABLE OF CONTENTS

PT400	2
New! PT400M	3
PT300	4
PT250	5
New! PT150	6
PT100	6
POWERLINE™	7
STORABLE MAT	8
PORTABLE MANIPULATION TABLE	9
STABILIZATION STRAPS	10
ACCESSORIES	11
New! BOIANCE™ FLOAT	12
PRONE PILLOW	12
STOOLS	13
FABRIC CHOICES	14



“Oakworks® PT Series table’s low 17” height range is ideal for lower patient positions often necessary for manipulations, NMRs and ADLs that require the therapist to establish treatment planes with significant leverage and it makes wheelchair transfers safe and easy. Its sturdy base enhances the therapist’s ability to establish effective manipulation vectors while being ergonomically adjustable to ensure safety for both the patient and the therapist.”


Dr. Robert DuVall,
PT, DHS, ATC, OCS, SCS, FAAOMPT CSCS

OAKWORKS PT SERIES TABLES

The perfect platform for PT work!


Table designs are PT facilitated by PT practitioners.

FIRM RESPONSE FOAM

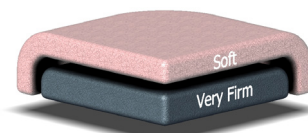

This padding has been designed to deliver a firm yet comfortable surface that maximizes Physical Therapists’ ergonomic efficiency during examination and intervention procedures. Additionally, this surface minimizes Manipulation Dry Hits with its super high density, twin layer design. The Firm Response foam enables increased reliability of patient examinations.

STIFF TWIN BEAM CONSTRUCTION

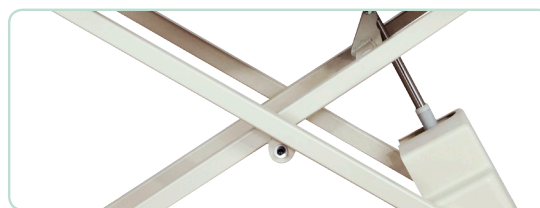

Stiffer than cantilever designs and with fewer moving parts, our design supports all your work with a 17”- 34” height range and a powerful ETL listed motor.

PT SERIES SAFETY LISTINGS

FDA listed to regulation #878.4960 - Class 1; CE Declaration of Conformity to Annex I of the EC Directives 90/385/EEC and 93/42EEC; cETLus Test Reports to IEC/EN 60601-1 3rd Edition, IEC/EN 60601-1-2 3rd Edition, IEC/EN 62366 1st Edition; Authorization to Mark (ATM) to IEC 60601-1-6 3rd Edition, ATM to ANSI/AAMI ES60601-1 3rd Edition, ATM to CAN/CSA-C22.2 No. 60601-1

OAKWORKS®
MEDICAL

made in the USA
with US & imported parts



PT400

550 LB. LIFTING CAPACITY

- Electric controlled center break
- Multiple angles & configurations
- Uniform stability maintained in 17" - 34" electric height range



Shown with optional
Locking Twin Casters


Foot Control
included!



**Watch
the video!**

SPECIFICATIONS	
Lifting Capacity	550 lbs.
Height Range	17"- 35"
Width	27", 29", 31"
Length	73"
Padding	1.75" Firm Response, 2.5" Comfort Foam™
Fabric	TerraTouch™
Controls	1 foot control (multi-function)
Options	<i>Hand Control, Locking Twin Casters</i>
Safety Listings	FDA, cETLus, CE Marked, CB Certificate <i>Details on page 1</i>


Bold=Standard, Italic=Additional Charge (except padding)

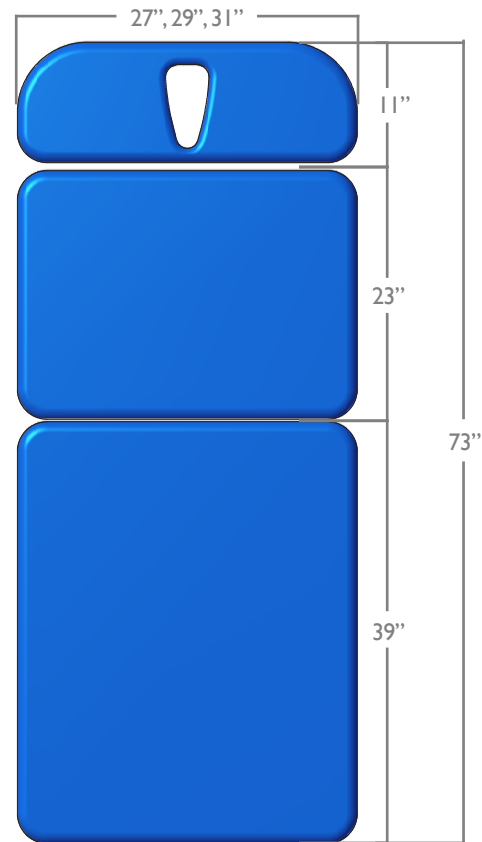


Chair
Position



FEATURES

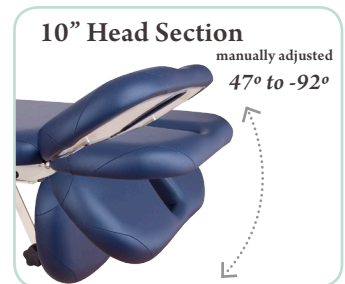


New PT400M

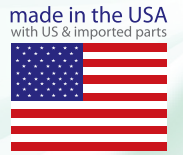


550 LB. LIFTING CAPACITY

- Powered center break & height
- Optional arm rests & table extenders
- Contoured ergonomic design
- Advanced head support included

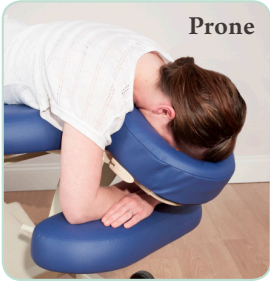

Prone

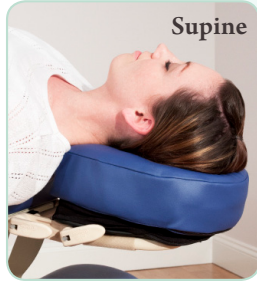

Supine




Watch the video!

SPECIFICATIONS

Lifting Capacity	550 lbs.
Height Range	18"- 36"
Width	29.5"
Length	77" (including head rest)
Padding	2.5" Comfort Foam™, 1.75" Firm Response
Fabric	TerraTouch™
Controls	1 foot control (multi-function)
Options	<i>Hand Control, Second Foot Control, Locking Twin Casters</i>
Accessories	<i>PT400M Rotating Arm Rest (set of 2), PT400M Extension Pad, PT400M Shelf Bracket</i>
Safety Listings	FDA, cETLus, CE Marked, CB Certificate <i>Details on page 1</i>

Bold=Standard, Italic=Additional Charge (except padding)

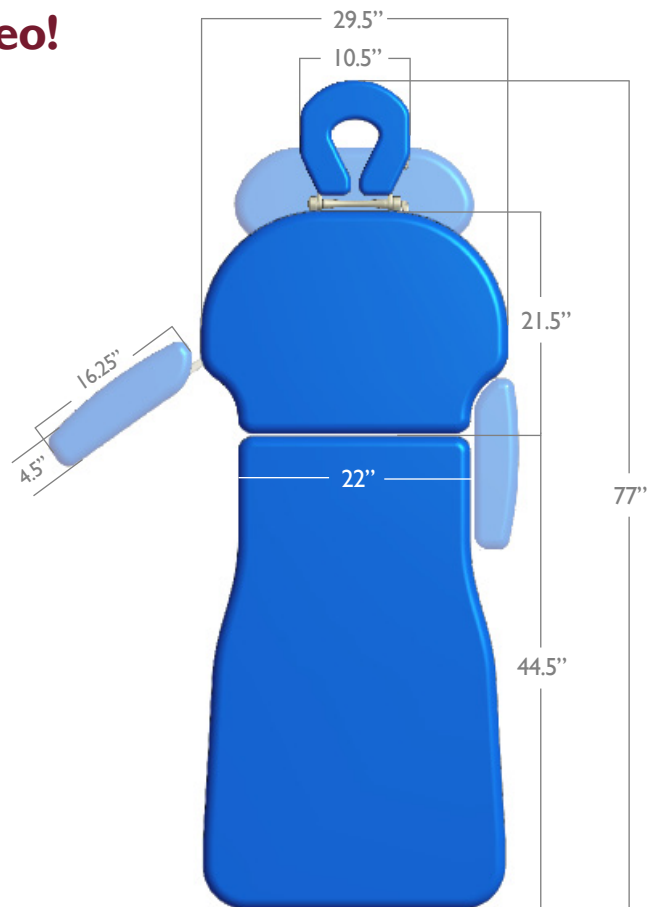

Powered Center Break
0° to 45°



Rotating Arm Rests
removable accessory



Removable Extension Pad
bracket not included

PT300

550 LB. LIFTING CAPACITY

- Super low 16" height
- Smooth High-Low excursion range
- Patient examination ergonomic efficiency
- Dix-Hallpike maneuver ergonomic efficiency
- Wide enough for Brandt-Daroff exercises for treatment of BPPV
- Firm Response foam minimizes Dry Hits
- Ease of VAT with drop head section
- Efficient/Dry needling patient positioning

Shown with optional
Locking Twin Casters



Foot Control included!

made in the USA
with US & imported parts



FEATURES

0° - 82° LIFT-ASSIST BACKREST

This easy to use system makes it possible to adjust the table effortlessly even when your patient is on the table. The powerful gas cylinders are controlled by easy to reach finger controls that lock the table securely in place after movement.



SPECIFICATIONS

Lifting Capacity	550 lbs.
Height Range	16"- 34"
Width	29", 31"
Length	73"
Padding	1.75" Firm Response, 2.5" Comfort Foam™
Fabric	TerraTouch™
Controls	1 foot control
Options	<i>Hand Control, Locking Twin Casters, Retractable Side Rail (31" width only)</i>
Safety Listings	FDA, cETLus, CE Marked, CB Certificate <i>Details on page 1</i>

Bold=Standard, Italic=Additional Charge (except padding)



**Watch
the video!**

PT250

550 LB. LIFTING CAPACITY

- Super low 16" height
- Vertebral Artery Testing efficiently and safely
- Optimize Examination/ Dry needling

made in the USA
with US & imported parts



Foot Control included!

SPECIFICATIONS

Lifting Capacity	550 lbs.
Height Range	16"- 34"
Width	29", 31"
Length	72"
Padding	1.75" Firm Response, 2.5" Comfort Foam™
Fabric	TerraTouch™
Controls	1 foot control
Options	<i>Hand Control, Locking Twin Casters</i>
Safety Listings	FDA, cETLus, CE Marked, CB Certificate <i>Details on page 1</i>


Bold=Standard, Italic=Additional Charge (except padding)

"The drop down head section feature enables patient safe and Physical Therapist efficient cervical examination and treatment."

Dr. Robert DuVall,
PT, DHS, ATC, OCS, SCS, FAAOMPT CSCS

FEATURES & OPTIONS

11" Head Section +40° to -90°

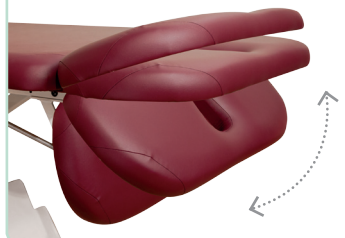

Optional Locking
Twin Casters

Enhances maneuverability in small spaces and allows PT to treat patients from either side of the table.



Watch the video!

New PT150

550 LB. LIFTING CAPACITY



Foot Control included!

- Super low 16" height
- Lift-Assist backrest 0°- 90°
- Stiff twin beam construction
- Mobilization
- Neuromuscular Re-education

"Multi-functional table for ergonomically efficient manual therapy, neuromuscular re-education both conveniently applied to therapeutic activities necessary to achieve optimal patient outcomes."

Dr. Robert DuVall,
PT, DHSc, ATC, OCS, SCS, FAAOMPT CSCS

SPECIFICATIONS for PT150 & PT100

Lifting Capacity	550 lbs.
Height Range	16"- 34"
Width	29", 31"
Length	73" (PT100), 75" (PT150)
Padding	1.75" Firm Response, 2.5" Comfort Foam™
Fabric	TerraTouch™
Controls	1 foot control
Options	<i>Hand Control, Locking Twin Casters</i>
Safety Listings	FDA, cETLus, CE Marked, CB Certificate <i>Details on page 1</i>

Bold=Standard, Italic=Additional Charge (except padding)

OPTION



Optional Locking Twin Casters
Enhances maneuverability in small spaces and allows PT to treat patients from either side of the table.


PT100

550 LB. LIFTING CAPACITY



Foot Control included!

- Super low 16" height
- Stiff twin beam construction
- Mobilization
- Neuromuscular Re-education
- Therapeutic Activities

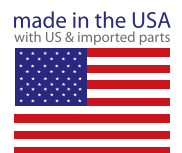

POWERLINE™

500 LB. WEIGHT CAPACITY

- Heavy-duty construction
- Strongest table in its class
- Easy setup
- Ships in one box with hardware included!
- 3 year guarantee



shown with Backrest top & optional shelf



Watch the video!

SPECIFICATIONS

Weight Capacity	500 lbs.
Height	30"
Width	27" or 30"
Length	72"
Top	Flat or Backrest
Padding	2.5" Semi-Firm, 1.75" Firm Response
Fabric	TerraTouch™
Accessories	<i>Stabilization Strap</i>
Options	<i>Storage Shelf</i>
Safety Listings	FDA, CE Marked

Bold=Standard, *Italic=Additional Charge*

FEATURES



HEAVY DUTY CONSTRUCTION

Super strong hardware used throughout, built to last!

OPTIONS *(additional charge)*



INSET SHELF

The optional shelf is notched at the corners so it stays in place, while providing added support and stability to the table. Offers easy access to storage & frequently used accessories.




shown with Flat top

made in the USA
with US & imported parts



STORABLE MAT TABLE

550 LB. WEIGHT CAPACITY

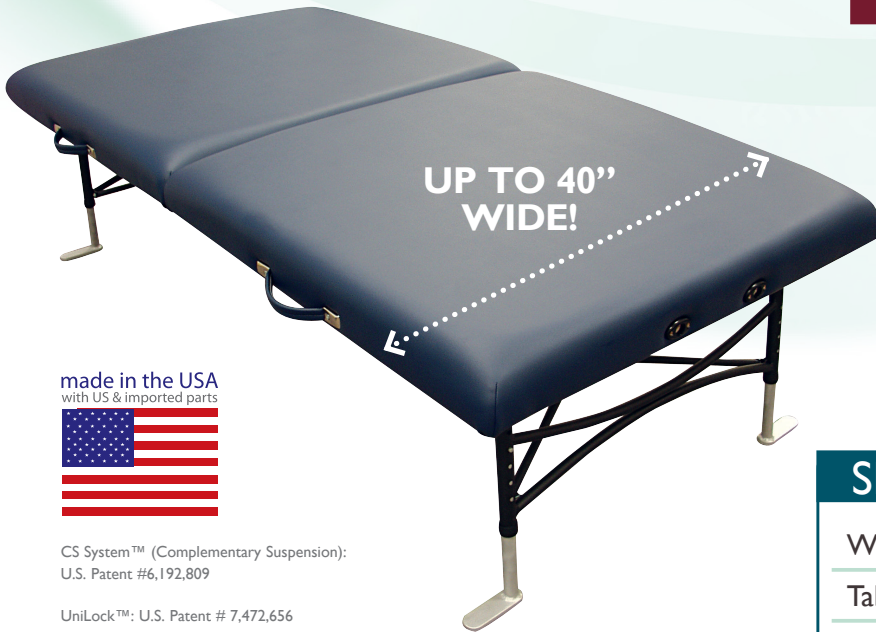

- 16"-25" height range
- Ultra stiff leg bracing
- Space saving design
- 40" width option

made in the USA
with US & imported parts



CS System™ (Complementary Suspension):
U.S. Patent #6,192,809

UniLock™: U.S. Patent # 7,472,656

SPECIFICATIONS

Weight Capacity	550 lbs.
Table Weight	40-55 lbs. <i>(depending on width)</i>
Height	16"- 25"
Width	33", 35", 40"
Length	78"
Top	Flat
Padding	2.5" Semi-Firm, <i>1.5" Firm Response</i>
Fabric	TerraTouch™
Accessories	<i>Stabilization Strap (33" & 35" only)</i>
Options	<i>Stabilizing Feet (4 included with 40")</i>
Safety Listings	FDA, CE Marked

The **ONLY** portable mat table on the market with a **550 lb. Weight Capacity!**

The Oakworks® storable mat is as flexible as it is stable. Perfect for space constrained facilities where mat work is required, the portable mat table provides up to a 40" by 78" working surface while easily stores out of the way in a super-compact suitcase sized space.



QUICK PORTABILITY

With its space saving design, the Storable Mat folds up easily for storage and is ready to follow wherever you go.



STABILIZING FEET

For maximum stability each 40" Storable Mat comes standard with stabilizing feet, which easily fold out of the way for storage.

Bold=Standard, Italic=Additional Charge (except padding & 40" width)


Ultra stiff leg bracing holds up under heavy lateral force



PORTABLE MANIPULATION TABLE

550 LB. WEIGHT CAPACITY

- Adjustable Heights
- Easily portable at only 25 lbs.
- Optimal body mechanics for any therapist
- Firm Response foam minimizes Dry Hits and increases reliability of examinations.
- Designed for Manual Manipulation, Range of Motion & Manual Muscle Testing

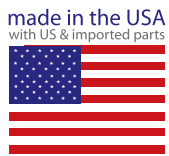

CS System™ (Complementary Suspension):
U.S. Patent #6,192,809
UniLock™: U.S. Patent # 7,472,656



“The Portable Manipulation Table from Oakworks® offers an outstanding platform for performing thrust and non-thrust manipulation. The team at Oakworks® has listened to the feedback of practicing manual physical therapists and created a high value and excellent performing product. You can’t ask anything more out of a portable manipulation table!”

Timothy W. Flynn, PT, PhD

Desirable portability for Physical Therapists in school, internship, residency or fellowship training environments.



SPECIFICATIONS

Weight Capacity	550 lbs.
Table Weight	25 lbs.
Height	24" - 34"
Width	27"
Length	73"
Top	Flat
Padding	1.5" Firm Response
Fabric	TerraTouch™
Accessories	<i>Carry Case, Stabilization Strap</i>
Options	<i>Sports Stripe</i>
Safety Listings	FDA, CE Marked

Bold=Standard, Italic=Additional Charge



ADJUSTABLE

24"-34" height range is ergonomically correct for any therapist.



DURABLE

TerraTouch™ fabric & Firm Response foam matched with our strong aluminum frame provides endless stability.



STABILIZATION STRAPS


GREAT FOR RANGE OF MOTION, MOBILIZATION & ISOLATION

FDA listed

- Adjustable in length over a wide range for stabilization of any body part
- Designed to work with Oakworks® PowerLine™ table, Oakworks® Portable Manipulation table and Storable Mat (excluding 40”).
- Secure and strong, heavy duty design will last for years
- Wide straps offer more surface area for client comfort
- Great for range of motion, mobilization, isolation



Watch the video!

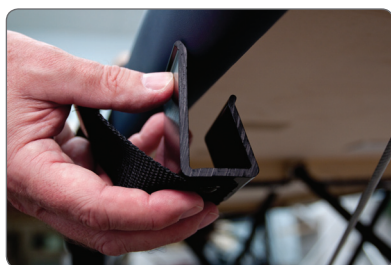

ACCESSORIES

made in the USA
with US & imported parts



EVERYTHING YOU NEED TO POSITION YOUR PATIENTS

AERO•CEL™ BOLSTERS *Available in all TerraTouch™ colors!*

Oakworks® Aero•Cel™ bolsters are 100% Polyurethane fabric - no PVC - with pull handles at both ends.



6" Half-round
3"H x 26"L

8" Half-round
4"H x 26"L

6" Semi-round
4.5"H x 26"L
most popular

8" Semi-round
6.5"H x 26"L

6" Round
6"D x 26"L

8" Round
8"D x 26"L



PILLOW COVERS

Easy to clean, the pillow cover is a hygienic barrier between your pillow and client. Choose from standard or king size.



RECTANGULAR BOLSTER

Designed specifically to support a patient's leg when in a side lying position, this rectangular bolster offers two levels of height and ample surface area. 7"W x 27"H x 10"L



WEDGE BOLSTER

Easily used as a backrest or leg support. 24" W x 10" H x 24" L with a 36° incline.



TABLE EXTENDERS

Extend the length of your table. Insert into the table's accessory outlets or store out of the way when not in use.



BOIANCE™ FACE REST CRESCENT

Water spheres within the multi-layer foam support the face and weight of the head and delivers comfort while avoiding sinus pressure.



QUICKLOCK™ FACE REST PLATFORM & AERO•CEL™ FACE REST CRESCENT


The perfect combination for accurate positioning of the cervical spine. Platform available in Black or Vanilla.

These products are not FDA Listed.

Introducing THE BOIANCE™ FLOAT

PATENT PENDING WATER SPHERE TECHNOLOGY

FDA listed

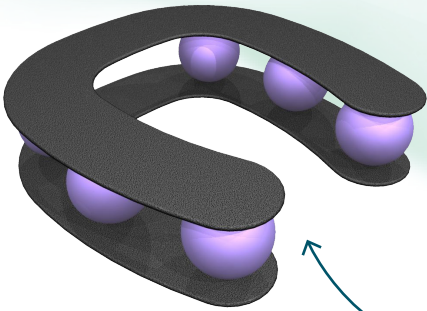

shown with the Boiance™ face rest crescent

- Incredible freedom of movement in all directions
- Eliminate compression of the cervical spine
- Just insert between cradle and pad!

6 Water Spheres, the same ones used in the famous Boiance Face Cradle Pillow, float the pillow and isolate it from the platform and table. Just place the Float between your face rest pad and platform.



Watch the video!

PRONE PILLOW

FOR PERFECT PRONE & SUPINE POSITIONING




- Increased ventilation & improved comfort
- Low profile reduces cervical extension/flexion
- Superior positioning for any size patient
- Removable, easy to clean pads



20" L x 20.5" W x 6.25" H

US patent # D555,251 S
(Made in China)

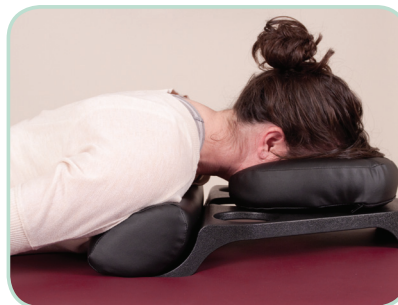

PROPRIETARY PAD SYSTEM

- Removable pads are easy to clean & less costly to replace
- High quality covered soft cushions last much longer and are more comfortable than any vinyl coated rigid foam, and they are replaceable!

SPECIFICATIONS	
Weight	2.6 lbs.
Dimensions	20" L x 20.5" W x 6.25" H
Removable Pads	Face Rest, Bolster
Fabric	TerraTouch™ Coal
Safety Listings	FDA, CE Marked

STOOLS

300 LB. WEIGHT CAPACITY

- Super high density, long lasting foam
- Professional quality gas-assist system
- High quality TerraTouch™ fabric

Assembled in the USA
with imported parts



PREMIUM STOOL with BACKREST

Colors: All TerraTouch™ Colors
Low Height: 18" - 24"
High Height: 20" - 28"



PREMIUM STOOL

Colors: All TerraTouch™ Colors
Low Height: 18" - 24"
High Height: 20" - 28"



PREMIUM STOOL with SADDLE SEAT

Colors: All TerraTouch™ Colors
Low Height: 19" - 25"
High Height: 22" - 29"



PROFESSIONAL STOOL

All TerraTouch™ Colors
Low Height: 18" - 24"
High Height: 21" - 28"



PROFESSIONAL STOOL with SADDLE SEAT

All TerraTouch™ Colors
Low Height: 19" - 25"
High Height: 22" - 29"



BASIC STOOL

Colors: Black, Heron or Ruby ONLY
Height Range: 18" - 24"
(Made in China)

These products are not FDA Listed.


FABRIC

DURABLE & COMFORTABLE

TERRATOUGH™ PVC & Latex FREE!

TerraTouch™ fabric is designed to meet our high standards for longevity and suppleness. An Oakworks' exclusive, this PVC-free premium polyurethane has excellent resistance to oil, abrasion, and tears. TerraTouch™ fabric will meet the table requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables and will pass testing for Cytotoxicity, Skin Irritation and Sensitization as well as CA TB117 Fire Retardance testing. For disinfecting purposes you may use Protex, MadaCide, MadaWipe, Accell TB, Virox® (Canada), PDI Duper Sani-Cloth Germicidal wipe or a diluted 1/10 bleach solution.

Please refer to operating manuals for cleaning guidelines.

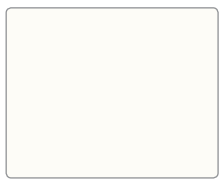


Pure White-T26



Butter-T09



Saffron-T21

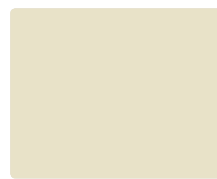

Opal-T12



Clay-T18



Taupe-T28



Earth-T07



Espresso-T27



Ruby-T08



Orchid-T19



Heron-T05




Ocean-T20



Sapphire-T17



Sky Blue-T04

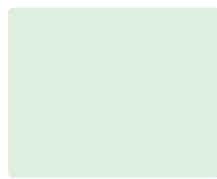

Seafoam-T30



Sage-T22



Blue Grass-T24




Forest-T14



Stone-T13



Pewter-T29

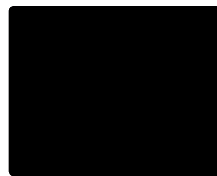

Coal-T01

We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate. In order to receive an accurate representation of color, please contact us for a sample.

OAKWORKS[®] MEDICAL

QUALITY POLICY

It is the policy of Oakworks[®] to provide world class products and services that consistently meet or exceed our customers' expectations. We pledge to diligently maintain the integrity of our Quality Management Systems in our pursuit of total customer satisfaction.



POSITIONING PATIENTS
for POSITIVE OUTCOMES

OAKWORKS[®] MEDICAL

Serving many specialities...



Imaging Tables



Procedure Chairs



Ultrasound Tables

 **HENRY SCHEIN[®]**
www.henryschein.com


GLOBAL COOLING PARTNER



ISO 13485:2003
FM 584901

made in the USA
with US & imported parts

