


CAPROSYN™ SUTURE

CAPROSYN™ SUTURE





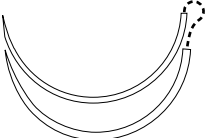
Description	<p>CAPROSYN monofilament synthetic absorbable sutures are prepared from POLYGLYTONE™ 6211 synthetic polyester which is composed of glycolide, caprolactone, trimethylene carbonate, and lactide.</p>
Characteristics	<p>The advanced extrusion process of the molecule of POLYGLYTONE™ 6211 gives the suture:</p> <ul style="list-style-type: none">■ Strong and secure wound approximation during the critical wound healing period■ Predictable absorption profile■ Excellent handling■ Excellent knot security
Indications	<p>CAPROSYN synthetic absorbable sutures are indicated for use in general soft tissue approximation and/or ligation, but not for use in cardiovascular or neurological surgery, microsurgery, or ophthalmic surgery.</p>
Contraindications	<p>This suture, being absorbable, should not be used where extended approximation of tissue is required.</p>
Actions	<p>CAPROSYN synthetic absorbable sutures elicit a minimal acute inflammatory reaction in tissues, which is followed by gradual encapsulation of the suture by fibrous connective tissue. Progressive loss of tensile strength and eventual absorption of CAPROSYN synthetic absorbable sutures occurs by means of hydrolysis. Absorption begins as a loss of tensile strength followed by a loss of mass.</p>
Tensile Strength	<p>Subcutaneous tissue implantation studies in animals indicate that CAPROSYN synthetic absorbable sutures retain a minimum of 50–60% USP and EP knot strength at 5 days post-implantation, and a minimum of 20–30% of USP and EP knot strength at 10 days post-implantation. All of the original tensile strength is essentially lost by 21 days post-implantation.</p>
Absorption	<p>Absorption begins as a loss of tensile strength without appreciable loss of mass. Absorption of CAPROSYN synthetic absorbable sutures is essentially complete by 56 days.</p>
Sterilization	<p>Ethylene oxide</p>
Packaging	<p>CAPROSYN synthetic absorbable sutures are available in USP sizes 1 (4 metric) through 6-0 (0.7 metric). They are available undyed (natural) and dyed (purple). The sutures are supplied sterile, in pre-cut lengths or ligating reels, non-needled or affixed to various needle types: in one, two, or three dozen boxes.</p>










NON-NEEDLE PRODUCTS	Length	Color	Units Per Box	M USP	6-0	5-0	4-0	3-0	2-0	0	1
				0.7	1	1.5	2	3	3.5	4	
PRE-CUT SUTURE	12x18	Undyed	24					C-2	C-13	C-14	C-15
	3x18	Violet	36							C-74	
REEL	60	Undyed	24					LC-112	LC-113	LC-114	
STANDARD LENGTH	60	Undyed	24					C-212			C-215
	60	Undyed	36						C-213	C-214	
REVERSE CUTTING ▼											
3/8 Circle 12mm  C-1	18	Undyed	12	C-3810							
3/8 Circle 12mm  C-1	18	Undyed	12		C-3711						



REVERSE CUTTING ▼	Length	Color	Units Per Box	USP								
				M	6-0 0.7	5-0 1	4-0 1.5	3-0 2	2-0 3	0 3.5	1 4	
3/8 Circle 19mm  C-13	18	Undyed	36			SC-3634						
	30	Undyed	36				SC-691	SC-692	SC-693			
3/8 Circle 24mm  C-14	30	Undyed	12					SC-693-G				
1/2 Circle 13mm  C-21	18	Undyed	12				C-1744					
1/2 Circle 37mm  GS-11	30	Undyed	36						CC-543			
1/2 Circle 37mm  GS-11, GS-21	60	Undyed	36					CC-982	CC-983			









REVERSE CUTTING ▼	Length	Color	Units Per Box	M USP	6-0	5-0	4-0	3-0	2-0	0	1
					0.7	1	1.5	2	3	3.5	4
1/2 Circle 40mm  GS-12	36	Undyed	36								CC-641
3/8 Circle 11.0mm Curvature 135° Radius 4.67mm Diameter 0.017"  HE-1	18	Undyed	12	✦ C-1780-GK							
3/8 Circle 11.0mm Curvature 135° Radius 4.67mm Diameter 0.017"  HE-1	18	Undyed	12	✦ C-1790-GK							
3/8 Circle 13.0mm Curvature 155° Radius 4.80mm Diameter 0.021"  HE-3	18	Undyed	12		✦ C-1792-GK	✦ C-1793-GK					
3/8 Circle 7.4mm Curvature 110° Radius 4.11mm Diameter 0.012"  HE-6	12	Undyed	12	✦ C-3794-GK							



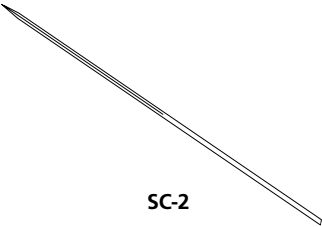




PREMIUM REVERSE CUTTING ▼	Length	Color	Units Per Box	USP M	6-0	5-0	4-0	3-0	2-0	0	1
					0.7	1	1.5	2	3	3.5	4
3/8 Circle 11mm P-10	18	Undyed	12	✦	SC-3625-G	SC-3626-G					
	18	Violet	12	✦	SC-3635-G						
3/8 Circle 13mm P-13	18	Undyed	12	✦	SC-5688-G	SC-5689-G	SC-5690-G				
	18	Violet	12	✦	SC-5588-G	SC-5587-G	SC-5590-G				
3/8 Circle 16mm P-11	18	Undyed	12	✦	SC-5635-G	SC-5626-G	SC-5686-G				
	18	Undyed	36			SC-5626					
3/8 Circle 19mm P-12	18	Undyed	12			SC-5616-G	SC-5617-G	SC-5618-G			
	18	Violet	12				SC-5627-G				
	30	Undyed	12				SC-5637-G	SC-5638-G			
3/8 Circle 24mm P-14	18	Undyed	36				SC-5679	SC-5640			
	30	Undyed	12					SC-5780-G			
	30	Undyed	36				SC-5779				
1/2 Circle 16mm P-24	18	Undyed	12	✦	SC-5641-G	SC-5642-G	SC-5643-G				
	18	Undyed	36			SC-5642	SC-5643				






PREMIUM X-CUTTING 	Length	Color	Units Per Box	M USP	6-0	5-0	4-0	3-0	2-0	0	1
					0.7	1	1.5	2	3	3.5	4
3/8 Circle 11mm  DX-11	18	Undyed	12			✦ XC-490					
3/8 Circle 13mm  DX-13	18	Undyed	12			✦ XC-834	✦ XC-835				
3/8 Circle 16mm  DX-16	18	Undyed	12			✦ XC-844	✦ XC-845				
3/8 Circle 19mm  DX-19	18	Undyed	12			✦ XC-822	✦ XC-823	✦ XC-497			
1/2 Circle 16mm  DXH-16	18	Undyed	12				✦ XC-814				



CUTTING ▼	Length	Color	Units Per Box	USP							
				M	6-0	5-0	4-0	3-0	2-0	0	1
				0.7	1	1.5	2	3	3.5	4	
Straight 11mm  SC-11	18	Undyed	12			SC-604					
Straight 11mm  SC-11	18	Undyed	12			SC-828					
Straight 60mm  SC-2	30	Undyed	36				SC-643				
TAPER POINT ⊕											
1/2 Circle 17mm  CV-23	30	Undyed	36	UC-201	UC-202	UC-203	UC-204				
	30	Violet	36		UC-215	UC-214	UC-213				
1/2 Circle 22mm  CV-25	30	Undyed	36			UC-403	UC-404				

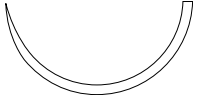

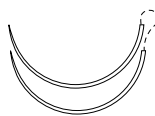



TAPER POINT	Length	Color	Units Per Box	USP							
				6-0	5-0	4-0	3-0	2-0	0	1	
				0.7	1	1.5	2	3	3.5	4	
1/2 Circle 22mm  GS-23	30	Undyed	36					UC-863	UC-864		
1/2 Circle 27mm  GS-22	30	Undyed	36					CC-883	CC-884	CC-885	
	30	Violet	36					CC-830	CC-831	CC-832	
1/2 Circle 37mm  GS-21	5X18	Undyed	12						CC-10-M	CC-102-M	
	5X18	Violet	12					CC-9-MG	CC-8-MG		
	5x30	Undyed	12							CC-6-M	
	5X30	Violet	12						CC-12-MG	CC-11-MG	
	3X30	Undyed	24						3CC-812	3CC-813	
	30	Undyed	36					CC-810	CC-811	CC-812	CC-813
	30	Violet	36						CC-843	CC-844	CC-845
	36	Undyed	36					CC-922	CC-923	CC-924	CC-925
	36	Violet	36						CC-953	CC-824	




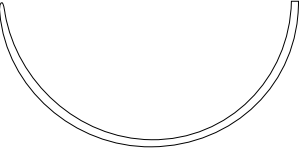


TAPER POINT	Length	Color	Units Per Box	M USP	6-0	5-0	4-0	3-0	2-0	0	1
					0.7	1	1.5	2	3	3.5	4
1/2 Circle 40mm GS-24	30	Undyed	36					CC-852	CC-801	CC-802	CC-803
	36	Undyed	36					CC-912	CC-915	CC-913	CC-914
	36	Violet	36					CC-966	CC-965	CC-964	CC-963
1/2 Circle 48mm GS-25	30	Undyed	36					CC-873	CC-872		
	36	Undyed	36							CC-904	CC-905
	36	Violet	36						CC-973	CC-974	CC-975
1/2 Circle 65mm GS-26	30	Undyed	24					CC-62	CC-63		
1/2 Circle 65mm GS-26	60 ³⁰	Undyed	24								CC-824-L
1/2 Circle 34mm HGS-20	D-TACH™ 5X18	Undyed	12							CC-100-M	CC-975-M
	D-TACH™ 5X18	Violet	12							CC-90-MG	

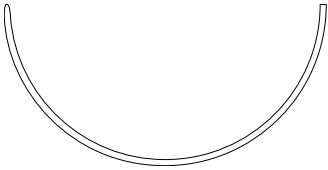
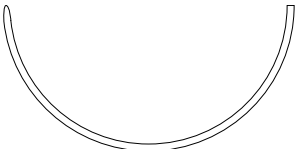



TAPER POINT	☉	Length	Color	Units Per Box	M USP						
					6-0	5-0	4-0	3-0	2-0	0	1
					0.7	1	1.5	2	3	3.5	4
1/2 Circle 37mm  HGS-21	▲ D-TACH®	Undyed	12							CC-103-M	CC-101-M
	5X18										
	5X18	Undyed	12								CC-201-M
	▲ D-TACH®	Violet	12							CC-113-MG	
	5X18										
	▲ D-TACH®	Undyed	12								CC-38-M
	5X30										
	▲ D-TACH®	Violet	12							CC-200-MG	
5X30											
30	Undyed	36							CC-942	CC-943	
30	Violet	36							CC-952		
1/2 Circle 26mm  V-20	▲ D-TACH®	Undyed	12				GC-69-M				
	5x18										
	30	Undyed	36			GC-121	GC-122	GC-123	GC-124		
	30	Violet	36			GC-321	GC-322	GC-323			
1/2 Circle 26mm  V-20	30	Undyed	36				GC-522				
1/2 Circle 37mm  V-26	30	Undyed	36					GC-127			



TAPER POINT	Length	Color	Units Per Box	M USP	6-0	5-0	4-0	3-0	2-0	0	1
					0.7	1	1.5	2	3	3.5	4
5/8 Circle 27mm  GU-46	30	Undyed	36					UC-979	UC-978		
	36	Undyed	36							UC-967	
	30	Violet	36					UC-879	UC-878	UC-877	
5/8 Circle 1/2 Circle 27mm 26mm  GU-46 V-20	30	Undyed	36						CC-558		
	30	Violet	36						CC-658		
5/8 Circle 37mm  GU-45	30	Undyed	36						UC-245	UC-246	
	30	Violet	36						UC-345		
BLUNT POINT - PROTECT POINT™											
1/2 Circle 64mm  BP-27	30	Undyed	24							CC-47	



BLUNT TAPER POINT - PROTECT POINT™	Length	Color	Units Per Box	USP							
				6-0 0.7	5-0 1	4-0 1.5	3-0 2	2-0 3	0 3.5	1 4	
1/2 Circle 85mm  BGS-28	36	Undyed	24								CC-49
1/2 Circle 65mm  BGS-29	30	Undyed	24							CC-50	
TAPERCUTTING											
1/2 Circle 37mm  KV-34	36	Undyed	12						CC-1945-G		
	36	Undyed	36				CC-962			CC-195	CC-196