



EasyOne® Products

ndd Medical Technologies is committed to setting new standards in pulmonary function testing by offering innovative, easy to use products and superior customer support.

- **Enhances patient satisfaction with portable testing**
- **Accommodates changing healthcare landscape**
- **Flexible EMR integration and export capabilities**



Your choice is an EasyOne

Point of care testing for fast, accurate and easy diagnosis.


Improve efficiency and quality of care by providing complete lung function testing on site

- Complete lung function in 1 square foot
- Automatic calibration
- No expensive service contracts
- No warm up time
- Challenge testing



EasyOne Pro

EasyOne Pro LAB

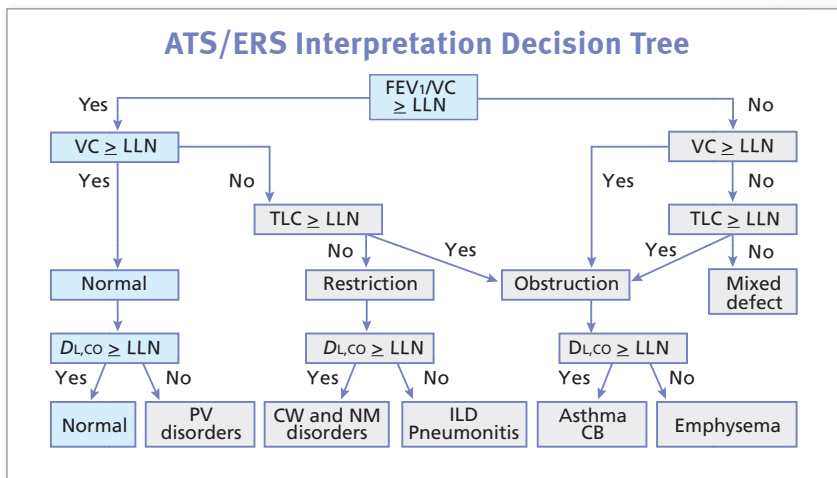

EasyOne Pro®

Portable single breath DLCO device allowing physicians to provide their patients with prompt, accurate diagnosis and treatment.


EasyOne Pro® LAB

All the benefits of the EasyOne Pro plus multiple-breath nitrogen washout to obtain complete lung volumes, FRC and LCI measurements.



Cart sold separately for use with EasyOne Pro or EasyOne Pro LAB


Spirometry is essential for the diagnosis and management of asthma and COPD

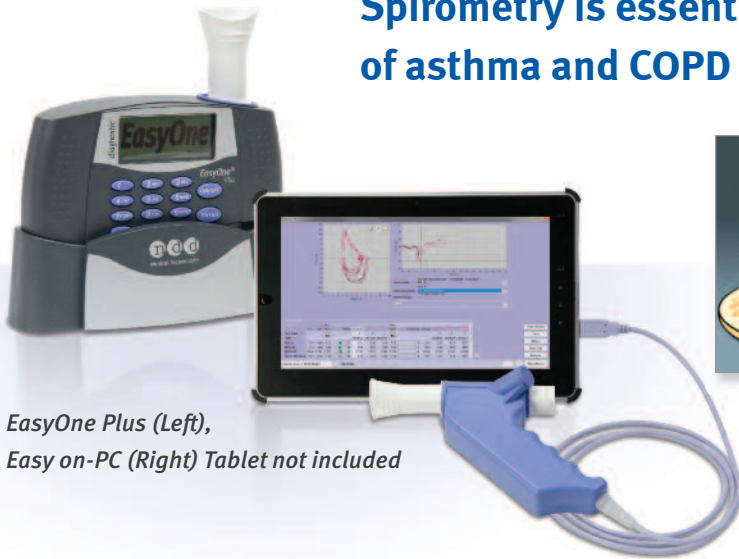

EasyOne Plus (Left),
Easy on-PC (Right) Tablet not included



Incentive Screens

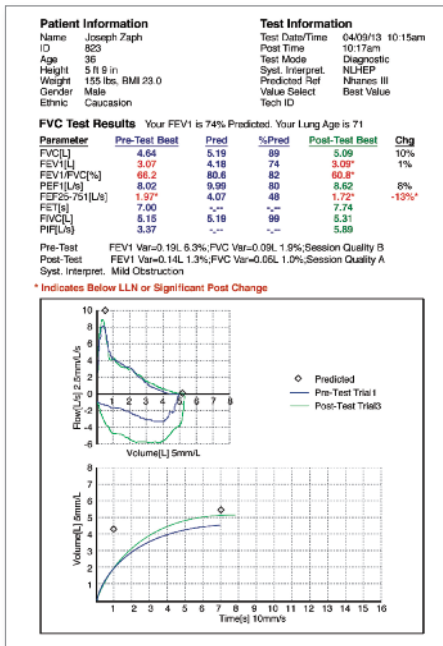

EasyOne® Plus

True portability for testing anytime, anywhere.

- Easy to use
- Instant results
- Powered by 2 AA batteries

Easy on-PC

This easy-to-operate spirometer uses the power of your PC, laptop, or tablet, leveraging premium ultrasound technology for a complete spirometry solution.

- Challenge testing
- Pediatric incentives
- Real time curves

Think ndd-Technology



No moving parts, no codes to enter, no screens to catch sputum, and no disposables to calibrate. The ultrasonic flow measurement is independent of gas composition, pressure, temperature, and humidity thus eliminating errors due to these variables.

Disposables

barriette®

A virtual filter validated for optimum hygienic separation of spirette/patient from the analyzer.



spirette®

The modern consumable for hygienic, precise, maintenance free measurement.

Systems

Part Numbers	Description	Includes
2001-2P	EasyOne Plus Diagnostic Spirometer	Diagnostic spirometer, USB cradle, printer & software
2001-2NP	EasyOne Plus Diagnostic Spirometer	Diagnostic spirometer, USB cradle & software
2700-3	Easy on-PC Spirometer	Spirometry sensor, software & carrycase
3000-00	EasyOne Pro	Single breath DLCO system with flow sensor, gas supply tube, one-way valve, DLCO valve, pressure-tube, temperature/humidity sensor and touch stylus; DLCO regulator and 2-year warranty
3100-00	EasyOne Pro LAB	Multiple breath nitrogen washout system with flow sensor, gas supply tube, one-way valve, DLCO valve, pressure-tube, temperature/humidity sensor and touch stylus; DLCO regulator and 2-year warranty

Disposables

Part Numbers	Description
2050-1	Case of 50 spirettes
2050-5	Case of 200 spirettes
3050-1	Case of 50 barriettes
3050-2	Case of 100 barriettes
3150-1	Case of 40 FRC barriettes (for use with EasyOne Pro LAB only)
3150-2	Case of 80 FRC barriettes (for use with EasyOne Pro LAB only)

Please contact your local distributor or ndd Medical Technologies
at 877-904-0090 for ordering information.

Proven Worldwide

- NASA, COPDGene, BOLD
- Meets ATS/ERS, NIOSH/OSHA and Disability standards



nddmedical



nddmedical



nddmed



nddmedical

