



MEDICAL PRODUCTS

· SPECIMEN TRANSPORT & SHIPPING PRODUCTS ·
· MEDICAL/SURGICAL DISPOSABLES ·

COM-PAC INTERNATIONAL
2006 PRODUCT CATALOG



2006 CATALOG



PRINTED IN THE USA

Com-Pac International, Inc.
800 Industrial Park Road
P.O. Box 2707
Carbondale, IL 62901
sales@com-pac.com
www.com-pac.com

TABLE OF CONTENTS

PACKAGING KITS

• Designing Solutions	3
• Infecon 2000	4
• Infecon 2050	5
• Infecon 3000	6
• Infecon 3050 FX-11	7
• Infecon Diagnostic Mailer	8
• Infecon Diagnostic Specimen Case	8
• Bitran Storage & Transport Box	9
• Diagnostic Foam Rack Kit	9

INSULATED PACKAGING

• Infecon 5000	10
• Infecon 5500	11
• Infecon 6000	12

SPECIMEN BAGS

• Bitran Series "S" Bags	13
• Bitran Series "PE" Bags	13
• Bitran Liquid-Tight Formalin Bags	14
• Infecon Specimen Transport Bags	14
• Inspex Liquid Tight Adhesive Specimen Bag	15
• Infecon 24 Hour Urine Collection Container Bags	15
• Saf-T-Zip Bags	16
• Pneumatic Tube Bags	16

ACCESSORIES

• Compartmental Bubble Pouches	17
• Custom Products	17
• Bitran Data Tags	18
• Absorbent	18
• Speci-six - 6 Compartment Zipper Bags	18
• Shipping Labels & Documents	19

MEDICAL/SURGICAL PRODUCTS

• Sponge Counter Bags	22
• Heart Monitor Bags	22
• Ice Bags	23
• Tray Covers	23
• Custom Products	24

DESIGNING SOLUTIONS

The time honored tradition of Com-Pac International innovation continues...for over 30 years, we have provided laboratories with high integrity packaging for specimen transport. We thank our customers for their loyalty and spreading the word - our Bitran bags remain the only guaranteed liquid tight zipper bags available. Good news travels and over the years, many healthcare practitioners from a variety of disciplines have come to us to design solutions for their specific needs.


With over 100 years of combined experience, our packaging engineers meet the challenge with your needs in mind. Merging compliance, packaging integrity and convenience, we provide cost effective solutions to your everyday challenges.

The first section of this catalog is dedicated to the Laboratory Specimen Transport Products, the stock line of products labs have depended on for many years. A special thanks to our customers who call and request the additional product offerings. Take a look at the new products in this issue:

- Inspex bags – liquid tight adhesive bags
- Pneumatic Tube bags
- Speci-Six - six compartment zipper bags
- Saf-T-Zip bags with absorbents
- Infecon bags - available with and without absorbent – for packaging 24 hour urine collection bottles.

We are pleased to provide an additional section to this publication highlighting a few of the Medical / Surgical Disposable products we manufacture. Feel free to call us with any questions or product ideas.

So share the news – tell your friends.

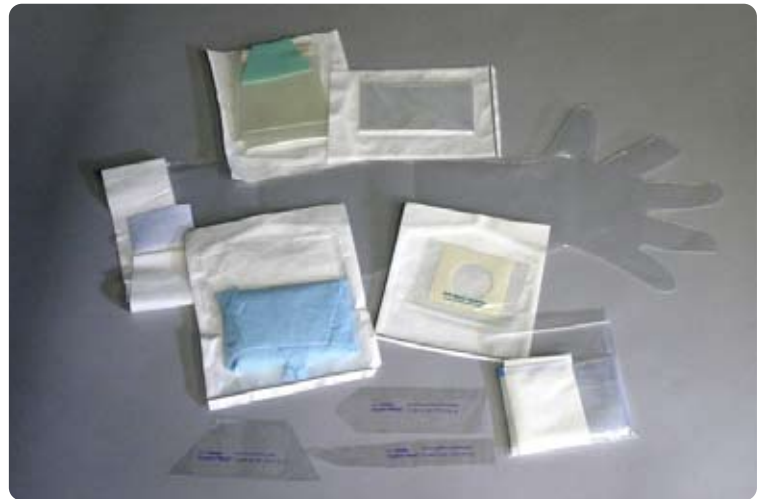

INFECON® 2000

Catalog No: INF-2000

Infectious Substance/Diagnostic Specimen Shipper - 6.2 U.N. Certified

The economical, dual purpose Infecon 2000 is the perfect container to use when shipping test tubes containing infectious substances or diagnostic specimens. UN 6.2 certified shipper includes:

- A liquid-tight .5 liter 95 kPa pressure vessel
- A perforated six compartmental bubble to provide protection
- 50ml capacity absorbent sheet
- An outer box with U.N. certification markings
- Infectious Substance and Diagnostic Specimen labels
- Shipper Declarations
- Easy-to-follow instructions
- Pull Tab/Locking Box

Count: 6/case

Empty Unit Weight: .38 lbs.

Dimensions: 4" x 4" x 8" w/.5 liter vessel

- Meets 49 CFR regulations for both air and ground shipments
- Meets ICAO/IATA Packing Instruction 602 and 650
- Certification available upon request

REPLACEMENT PARTS FOR INF - 2000

PART #	DESCRIPTION	QTY. PER UNIT
INF - 2000	Shipper	6/case
INF - 2003	50ml Capacity Absorbent Sheet	100/pack
INF - 2004	.5 liter Vessels (2" diameter 6" depth)	6/case
INF - 2005	Replacement Shipper Boxes	8/case
INF - 2011	Outer Shipping Box	10/case
INF - 9800	Shipper Declarations	150/pack
INF - 9830	2x2 Infectious Substance Label	50/pack
INF - 9847	2x2 Diagnostic Specimen Label	50/pack
BLP - 1000	Six Compartmental Bubble Pouch	25/case

INFECON 2050

Catalog No: INF-2050

Infectious Substance/Diagnostic Specimen Shipper - 6.2 U.N. Certified



Infecon 2050 utilizes all the popular features of the Infecon 2000 but has a larger outer box as preferred by some overnight carriers.

Count: 6/Case
Empty Unit Weight: .66 lbs.
Inside Dimensions: 6" x 6" x 8" w/.5 liter vessel

- Meets 49 CFR regulations for both air and ground shipments
- Meets ICAO/IATA Packing Instruction 602 and 650
- Certification available upon request

REPLACEMENT PARTS FOR INF - 2050

PART #	DESCRIPTION	QTY. PER UNIT
INF - 2050	Shipper	6/case
INF - 2003	50ml Capacity Absorbent Sheet	100/pack
INF - 2004	.5 liter Vessels (2" diameter 6" depth)	6/case
INF - 9800	Shipper Declarations	150/pack
INF - 9830	2x2 Infectious Substance Label	50/pack
INF - 9847	2x2 Diagnostic Specimen Label	50/pack
BLP - 1000	Six Compartment Bubble Pouch	25/case

INFECON 3000

Catalog No: INF-3000

Infectious Substance Shipper - 6.2 U.N. Certified



Count: 12/Case
Empty Unit Weight: 1.2 lbs.
Inside Dimensions: 6" x 6" x 8" w/1 liter vessel

- Meets 49 CFR regulations for both air & ground shipments
- Meets ICAO/IATA Packing Instruction 602 and 650
- Certification available upon request

The number one choice for shipping infectious substances, the Infecon 3000 is an easy-to-use packaging system. Each complete shipper includes:

- A liquid-tight 1 liter 95kPa pressure vessel
 - Pressure vessel is transparent for easy viewing of the contents
 - Pressure vessel is autoclavable for reuse
- Foam Stabilizer
- A perforated six compartmental bubble pouch to provide protection
- 50 ml capacity absorbent sheet
- An outer box with U.N. certification markings
- Mini Super Label which eliminates the need for individual labels
- Shipper Declarations
- Complete Instruction Guide

REPLACEMENT PARTS FOR INF - 3000

PART #	DESCRIPTION	QTY. PER UNIT
INF - 3000	Shipper	12/case
INF - 3001	Complete Refurbishment Kit Shipping boxes with required marking "O" rings; Foam stabilizers; Bubble Pouch; 50ml absorbent sheet; Shipper Declarations; Labels; Instructions	12/case
INF - 3002	Fifteen 6"x10" bubble wrap; Fifteen 50ml absorbent sheet	15/case
INF - 3003	36 compartmental bubble pouch; 36 50ml absorbent sheets 36 250ml absorbent sheets	36/case
INF - 3004	One liter vessels (4" diameter 6" depth)	6/case
INF - 3005	"O" Rings	36/pack
INF - 3011	Boxes and Labels for the INF-3000	12/bundle
INF - 9835	6x7 Infectious Substance Label	12/pack
BLP - 1000	Perforated six compartment bubble pouch	25/case
INF - 2003	50ml Capacity Absorbent Sheet	100/pack

INFECON 3050

Catalog No: INF-3050

FX-11 Infectious Substance Shipper - 6.2 U.N. Certified



Each complete shipper includes:

- A liquid-tight 1 liter 95kPa pressure vessel
 - Pressure vessel is transparent for easy viewing of contents
 - Pressure vessel is autoclavable for reuse
- A perforated six compartmental bubble pouch to provide protection
- 50 ml capacity absorbent sheet
- Shipper Declarations
- Infectious Label
- Foam Stabilizer
- An outer box with U.N. Certification markings
- FedEx Compliant

Count: 6/Case
Empty Unit Weight: 1.46 lbs.
Inside Dimensions: 12.75" x 6.5" x 8" w/1 liter vessel

- Meets 49 CFR regulations for both air and ground shipments
- Meets ICAO/IATA Packing Instruction 602 and 650
- Certification available upon request
- Meets IATA Operator Variation for Federal Express

REPLACEMENT PARTS FOR INF - 3050

PART #	DESCRIPTION	QTY. PER UNIT
INF - 3050	Shipper	12/case
INF - 3004	One liter vessels (4" diameter 6" depth)	6/case
INF - 9800	Shipper Declarations	150/pack
INF - 9835	6 x 7 Infectious Substance Label	12/pack
BLP - 1000	Perforated six compartment bubble pouch	25/case
INF - 2003	50ml Capacity Absorbent Sheet	100/pack

INFECON DIAGNOSTIC MAILER

Catalog No: DIA-1000



The Infecon Diagnostic Mailer contains:

- 9 x12 liquid-tight Infecon bag
- 250 ml absorbent
- Six compartment bubble pouch
- Sturdy outer box

Count: 36/Case
Empty Unit Weight: 25 lbs.
Unit Dimensions: 9.25" x 4.5" x 2.25"

INFECON DIAGNOSTIC SPECIMEN CASE

Catalog No: INF-3600



The Infecon Diagnostic Specimen Case is durable, leakproof, and convenient to use. It is great to use for transporting diagnostic specimens from one location to another within the lab. Infecon Diagnostic Case includes:

- Two Pick-n-Pluck foam inserts
- Two 250ml capacity absorbent sheets
- Holds 77 tubes per insert/2 inserts per case

Count: 1/Case
Empty Unit Weight: 4 lbs.
Inside Dimensions: 13" x 7.75" x 7.25"
Outside Dimensions: 15.5" x 8.5" x 9"

BITRAN® STORAGE & TRANSPORT BOXES

Catalog No: 4759-L



Convenient packaging designed to keep bagged specimens organized and oriented. This simple solution offers twenty 2.25" x 2.25" cells in its standard configuration and its adjustable grid system provides great flexibility, accommodating a variety of tubes, vials, bottles, and jars. The polyethylene coated grid is collapsible for space saving storage. Boxes and grids shipped flat.

Ideal for use with Bitran and Infecon liquid-tight specimen bags (pages 14-15). Versatile grid system also fits inside the Infecon Speci-Freez Insulated Shipper (page 12). Simply slide grid into the Speci-Freez shipper, and pack specimens with gel packs or dry ice for shipping temperature controlled specimens.

Count: 10/Case
Empty Unit Weight: 1.12 lbs.
Inside Dimensions: 10.25" x 12.25" x 8.5"

Use grid system in INF-6000

DIAGNOSTIC FOAM RACK KIT

Catalog No: FRB-2000



This kit includes:

- Pick-n-Pluck foam
- 14 x 12 gusseted liquid-tight Bitran bag with recloseable zipper
- Requisition pouch
- Two 250ml absorbent sheets
- OSHA-required biohazard symbol
- Holds 45 tubes

Count: 6/Case
Unit Dim: 9" x 5" x 4"

INFECON 5000

Catalog No: INF-5000

Insulated Infectious Shipper - 6.2 U.N. Certified



For a complete, integrated packaging system for shipping infectious substances on dry ice, the Infecon 5000 is a superior choice. Each complete shipper includes:

- Urethane insulation provides temperature protection exceeding 60 hours with only 5.6 pounds of dry ice
- An outer box with U.N. certification markings
- A liquid-tight .5 liter 95kPa pressure vessel
 - Pressure vessel is transparent for easy viewing of its contents
 - Sterilize with alcohol for reuse
- A perforated six compartment bubble pouch to provide protection
- 50ml capacity absorbent sheet
- Shipper declarations and compliant labels
- Packing instructions printed on shipper box flap
- Ideal for international shipping

Count: 6/Case
Empty Unit Weight: 3.1 lbs.
Inside Dimensions: 7.5" x 4.25" x 8.5" w/.5 liter vessel
Outside Dimensions: 11" x 8" x 12.25"

- Meets 49 CFR regulations for both air and ground shipments
- Meets ICAO/IATA Packing Instruction 602 & 904
- Certification available upon request

REPLACEMENT PARTS FOR INF - 5000

PART #	DESCRIPTION	QTY. PER UNIT
INF - 5000	Shipper	6/case
INF - 5001	Complete Refurb	6/case
INF - 2003	50ml absorbent sheet	100/pack
INF - 2004	.5 liter vessels (4" diameter 6" depth)	6/case
INF - 5005	Shipping Boxes	8/bundle
INF - 9800	Shipper Declarations	150/pack
INF - 9855	1S/01 Super Label	10/pack
INF - 9860	Dry Ice Super Label	10/pack
INF - 9875	Insulated Infectious Substance/Dry Ice Return	10/pack
BLP - 1000	Six compartment bubble pouch	25/case

INFECON 5500

Catalog No: INF-5500

Insulated Infectious Substance Shipper - 6.2 U.N. Certified



For maximum performance, order the Infecon 5500 insulated shipper. Each complete shipper includes:

- Temperature protection exceeding 100 hours with only 16 pounds of dry ice
- An outer box with U.N. certification markings
- A liquid-tight 1 liter 95kPa pressure vessel
 - Pressure vessel is transparent for easy viewing of its contents
 - Pressure vessel is autoclavable for reuse
- 100ml capacity absorbent sheet
- A six compartment bubble pouch to provide protection
- Shipper declarations and compliant labels
- Packing instructions printed on shipper box flap
- Ideal for international shipping

Count: 4/Case
Empty Unit Weight: 5.12 lbs.
Inside Dimensions: 9.5" x 7.5" x 9.5" w/1 liter vessel
Outside Dimensions: 13" x 10.75" x 13"

- Meets 49 CFR regulations for both air and ground shipments
- Meets ICAO/IATA Packing Instruction 602 and 904
- Certification available upon request

REPLACEMENT PARTS FOR INF - 5500

PART #	DESCRIPTION	QTY. PER UNIT
INF - 5500	Shipper	4/case
INF - 5501	Complete Refurb	4/case
INF - 2003	50ml absorbent sheet	100/pack
INF - 3004	One liter vessels (4" diameter 6" depth)	6/case
INF - 5505	Shipping boxes	4/case
INF - 9800	Shipper Declarations	150/pack
INF - 9855	Infectious Substance Dry Ice Super Label	10/pack
INF - 9860	Dry Ice Super Label	10/pack
INF - 9875	Insulated Infectious Substance/Dry Ice Return	10/pack
BLP - 1000	Six compartmental bubble pouch	25/case

INFECON 6000

Catalog No: INF-6000

Speci-Freez® Insulated Overpack Shipper



Count: 4/Case
Empty Unit Weight: 2.8 lbs.
Inside Dimensions: 12.25" x 10.25" x 8"
Outside Dimensions: 15.5" x 14" x 12.5"

This shipper is the perfect overpack insulated container for shipping infectious substances or diagnostic specimens which must be kept frozen or refrigerated. Specifreez shipper includes:

- Temperature protection for 76 hours which holds approximately 22 pounds of dry ice
- Expanded Polystyrene Insulated Tub
- Infectious Substance and Dry Ice Super Label
- Sturdy Outer Box
- Dry ice super label
- Infectious substance and dry ice return label which encourages the return of the package for reuse
- Instruction sheet
- Corresponding products to use with the Specifreez are: **INF-2000, INF-2050, INF-3000, 4759-L**

Important Note:

Air Eligibility when using our eligible inner package

Ideal for use with Grid System



REPLACEMENT PARTS FOR INF - 6000

PART #	DESCRIPTION	QTY. PER UNIT
INF - 6000	Shipper	4/case
INF - 6001	Outer boxes	10/bundle
4757-L	Versatile grid system	10/case
INF - 9800	Air Eligibility Labels	50/pack
INF - 9850	Infectious Substance/Dry Ice Super Label	150/pack
INF - 9860	Dry Ice Super Label	10/pack
INF - 9875	Infectious Substance/Dry Ice Return Label	10/pack

Specifreez is a registered trademark of Com-Pac International, Inc.

BITRAN SERIES “S” BAGS



Independently tested and certified liquid tight. Reports available upon request.

- The Series “S” 4740-S through 4745-S are packaged in packs of 250 each, four packs per case. The 16”x 16” and 24”x 24” bags are packaged in quantities of 100 per case.

- For long term handling and storage of specimens:
- Patented manufacturing process eliminates leaks and offers unprecedented safety;
 - Double-track zipper is airtight and easy to close;
 - Appropriate for use as a primary container for most diagnostic specimens;
 - Manufactured from durable 3-mil Saranex® material which provides an exceptional vapor barrier and odor control; custom sizes and printing available.
 - Custom sizes and printing available

CATALOG #:	DIMENSIONS:
4740 - S	2” x 4”
4741 - S	3” x 6”
4742 - S	6” x 6”
4743 - S	7” x 8”
4744 - S	9” x 12”
4745 - S	12” x 12”
4746 - S	16” x 16”
4747 - S	24” X 24”

BITRAN SERIES “PE” BAGS



The Series “PE” 4752-PE are packaged in packs of 250 each, four packs per case. The 16”x 16” and 24”x 24” bags are packaged in quantities of 100 per case.

Independently tested and certified liquid tight. Reports available upon request.

Bitran is a registered trademark of Com-Pac International, Inc.
Saranex is a registered trademark of Dow Chemical Corp.

- For short-term handling and storage of specimens:
- Patented manufacturing process eliminates leaks and offers unprecedented safety;
 - Double-track zipper is airtight and easy to close;
 - Appropriate for use as a primary container for most specimens;
 - Manufactured from durable 2.8-mil polyethylene that remains flexible when frozen;
 - Custom sizes and printing available.

CATALOG #:	DIMENSIONS:
4752 - PE	3” x 6”
4753 - PE	6” x 6”
4754 - PE	7” x 8”
4755 - PE	9” x 12”
4756 - PE	12” x 12”
4758 - PE	16” x 16”
4759 - PE	24” x 24”

BITRAN LIQUID-TIGHT FORMALIN BAGS



Safely handling and storing specimens in formalin solution is easy with Bitran Formalin Bags. This liquid-tight bag with vapor barrier effectively prevents the transmission of harmful fumes associated with formalin. This unique bag is made by a patented manufacturing process which eliminates leaks. Made from strong 3-mil Saranex, the Bitran Formalin Bag provides long-term protection. For added convenience, the bag is preprinted with a biohazard symbol, formalin warnings, and emergency response block as required by OSHA. The double-track zipper provides easy closure and airtight security. Packaged in packs of 250 each, four packs per case.

CATALOG #: **DIMENSIONS:**

6741 - S	3" x 6"
6743 - S	7" x 8"
6744 - S	9" x 12"
6745 - S	12" x 12"

INFECON SPECIMEN TRANSPORT BAGS



For the secure transport and handling of diagnostic specimens, the Infecon Specimen Transport Bag is the perfect choice. Infecon Bags are specially manufactured using a patented process which guarantees the elimination of leaks. The handy document pouch prevents misplacing vital paperwork. Made of durable 2.8-mil polyethylene, the Infecon Bag remains flexible when frozen. Preprinted with the OSHA required biohazard symbol, the bag also comes with a strong, easy-to-close orange double-track Saf-T-Zip closure. The larger sizes meet the need for transporting oversized diagnostic specimens. These bags are available in custom sizes and with custom print. The 5020-IP through the 5024-IP bags are packaged in units of 250 per pack, four packs per case; the 5025-IP, 5026-IP and 5028-IP bags are packaged in units of 100 per case.

CATALOG #:	DIMENSIONS:
5020 - IP	6" x 6"
5021 - IP	6" x 9"
5023 - IP	8" x 10"
5022 - IP	9" x 12"
5024 - IP	12" x 12"
5025 - IP	16" x 16"
5026 - IP	24" x 24"
5028 - IP	13" x 18"

*Independently tested and certified liquid tight.
Reports available upon request.*

INSPEX® LIQUID TIGHT ADHESIVE SPECIMEN BAG



Designed for Tamper Evident specimen protection during transport, the InspeX bags provide liquid tight containment and security. Attempts at tampering are evident - you must destroy the bag to remove its contents. The OSHA required biohazard symbol and simple closing instructions are preprinted on the bag. The document pouch keeps requisitions separate and dry. A recloseable adhesive strip on the document pouch secures paperwork during transport. InspeX bags are packaged 1,000 bags per case – 10 bundles of 100 bags each for inventory control.

CATALOG #:	DIMENSIONS:
8021 - XP	6.5" x 10" ID
8022 - XP	6.625" x 12" ID

INFECON 24 HR. URINE COLLECTION CONTAINER BAGS



State of the art - Infecon Liquid Tight bags that reduce the risk of danger associated with the spill of specimens and contamination. The challenge of moving a 24 hour urine collection is removed when you combine the integrity of Infecon with a 300 ml absorbent sheet. When it comes in contact with liquid, the absorbent turns to a gel-like substance. A convenient document pouch provides a safe and separate location for requisitions. Packaged 100/case.

CATALOG #:	DIMENSIONS:	ABSORBENT:
5028 - IP - 300	13" x 18"	YES
5028 - IP	13" x 18"	NO

SAF-T-ZIP® BAGS



Saf-T-Zip is a registered trademark of Com-Pac International, Inc.

Com-Pac's exclusive Saf-T-Zip bags are a cost-effective alternative for specimen transportation.

- Color coded for easy sorting
- OSHA required biohazard symbol
- Document pouch
- Packaged 100 per bundle
- 10 bundles per case

For added convenience, order the Saf-T-Zip with 50ml absorbents all ready in the bag.

Size: 6" x 9" with document pouch

CATALOG #	COLOR	ABSORBENT
7250-KT	Assortment	50 ml
7250-RD	Red	50 ml
7250-OR	Orange	50 ml
7250-YE	Yellow	50 ml
7250-GR	Green	50 ml
7250-BL	Blue	50 ml

Size: 6" x 9" with document pouch

CATALOG #	COLOR	ABSORBENT
7210-KT	Assortment	N/A
7210-RD	Red	N/A
7210-OR	Orange	N/A
7210-YE	Yellow	N/A
7210-GR	Green	N/A
7210-BL	Blue	N/A

PNEUMATIC TUBE BAGS




All in one protective packaging for transporting vials, tube and bottle through pneumatic tube systems.

Each kit contains:

- Compartment bubble pouch
- Liquid-tight Bitran zipper bag with pocket
- Instructions for use
- Packaged 100 bags/case

CATALOG#:	DIMENSIONS:
DA-1718	7" x 8"
DA-1912	9" x 12"

COMPARTMENTAL BUBBLE POUCHES

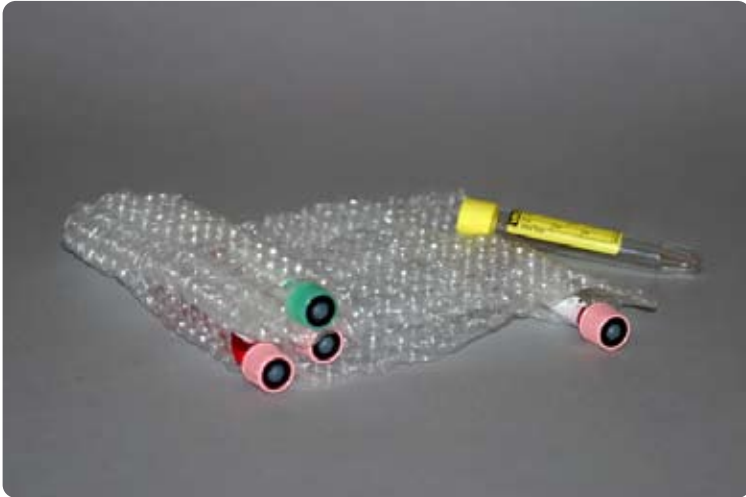

When transporting specimen tubes, use the economical BLP-1000 Bubble Pouches. These Bubble Pouches feature six perforated 2" compartments and are designed to keep specimens separated. The protective cushioning material reduces risk of breakage.

Perforated compartments can be torn apart into six single 2" pouches. Packaged 25 pouches/case.

CATALOG#: DIMENSIONS:
BLP-1000 13.75" x 6"

CUSTOM PRODUCTS



Our custom products meet each customer's high performance requirements. Bags can be imprinted with up to six colors for facility logs, special instructions, or handling information. Choose features such as document pockets, gussets, tape, and zipper colors.

Custom specimen shippers are designed for a facility's specific needs with a variety of containment systems. Our innovative packaging can be your advantage.

SPECI-SIX - 6 COMPARTMENT ZIPPER BAGS



Revolutionary design for transporting 6 test tubes together while keeping them separated from each other. Saf-T-Zip double track zipper closure provides easy access and convenience. Manufactured with 3 mil low density polyethylene printed with numbered compartments. Designed for multiple blood draws from a single patient and ideal for test kits.

CATALOG#: 6100 - ST
COUNT: 1000/case
DIMENSIONS: 10" x 6"

BITRAN DATA TAGS



These tags feature a two-part specimen information form for attachment to all specimen bags. These waterproof tags prevent ink from smearing and the pressure sensitive tape holds one part to the bag while the other part remains in patient records.

CATALOG#: 4757 - L
COUNT: 250/pack, 4 packs/case
DIMENSIONS: 3.5" x 6"

ABSORBENT



INFECON 50-ML CAPACITY ABSORBENT

Catalog No: INF - 2003
Count: 100/pack
Dimensions: 3" x 3" sheet

INFECON 100-ML CAPACITY ABSORBENT

Catalog No: INF - 1003
Count: 100/pack
Dimensions: 3" x 6" sheet

INFECON 250-ML CAPACITY ABSORBENT

Catalog No: INF - 1004
Count: 100/pack
Dimensions: 4" x 8" sheet

Infecon absorbent fast and reliable for that extra protection you require.

SHIPPING LABELS & DOCUMENTS



INFECON 9825 INFECTIOUS SUBSTANCE SUPER LABEL

Catalog #: INF-9825
Count: 10/pack
Dimensions: 10.5 x 6"

When shipping infectious substances at ambient temperatures, choose the one-piece Infecon Super Label 9825. By eliminating the need for ordering and placing individual labels on shipments, this all-inclusive label saves time and eliminates frustration. The Super Label has the appropriate spaces to identify the shipment of infectious substances, the proper U.N. markings, spaces for the shipper and recipient addresses, and the required emergency phone number. Just fill out the information and affix the label to the outer box of the package.

INFECON 9850 INFECTIOUS SUBSTANCE / DRY ICE SUPER LABEL

Catalog #: INF-9850
Count: 10/pack
Dimensions: 10.5 x 8.25"

When shipping infectious substances on dry ice, eliminate the need for placing individual labels on shipments by choosing the Infecon 9850 Super Label. The one sheet design contains all necessary labels for both infectious substances and dry ice, and the required U.N. marking. Spaces are provided for the shipper and receiver addresses, the infectious substance identification, and the dry ice weight.

INFECON 9800 SHIPPER'S DECLARATIONS

Catalog #: INF-9800
Count: 150 Declarations

Dangerous Goods Declaration form for shipping infectious substances and/or for shipping specimens requiring dry ice.

INFECON 9860 DRY ICE SUPER LABEL

Catalog #: INF-9860
Count: 10/pack
Dimensions: 10.5" x 5.75"

When shipping diagnostic specimens on dry ice, the DI Super Label simplifies the task of labeling. Just fill in the shipper and receiver addresses, and the weight of the dry ice, then apply the Super Label to the outer box—easy, quick, and economical.

INFECON 9875 INSULATED INFECTIOUS SUBSTANCE/ DRY ICE SHIPPERS RETURN LABEL

Catalog #: INF-9875
Count: 10/pack
Dimensions: 10.5" x 8"

For returning empty insulated shippers which previously contained substance and dry ice, order the Infecon 9875 Return Label. Place the label over the IS/DI label used to ship the package. The return label allows the empty shipper to be returned as an undeclared package. Order this return label to place over the INF-9850 and the INF-9825 labels.

LABELS FOR SHIPPING INFECTIOUS SUBSTANCES & DIAGNOSTIC SPECIMENS



1. Infecon 9830 2 x 2 Infectious Substance Label

Catalog #: INF-9830

Count: 50/Pack



2. Infecon 9835 6.2 Infectious Substance Mini Super Label

Catalog #: INF-9835

Count: 12/Pack

Dimensions: 6" x 7"



3. Infecon 9836 Class 9 Dry Ice Mini Super Label

Catalog #: 12/Pack

Dimensions: 6" x 6"



4. Infecon 9838 Inner Packaging Complies Label

Catalog #: INF-9838

Count: 12/Pack

Dimensions: 4" x 2.75"



5. Infecon 9839 Cargo Aircraft Only Label

Catalog #: INF-9839

Count: 12/pack

Dimensions: 4.75" x 4.5"



6. Infecon 9847 Diagnostic UN 3373 Specimen Label-Diamond

Catalog #: INF-9847

Count: 50/Pack

Dimensions: 2" x 2"




7. Infecon 9844 Diagnostic UN 3373 Specimen Label

Catalog #: INF-9844

Count: 12/Pack

Dimensions: 4" x 4"

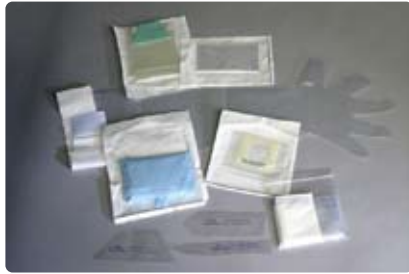


MEDICAL/SURGICAL DISPOSABLE PRODUCTS

Com-Pac International, Inc. is proud to announce the addition of our Medical/Surgical products. Because you, our customer, have requested it. Com-Pac International is committed to the development and improvement of your specialized medical products for the operating room, Patient Care, Cardiology, and more. We achieve this by seeking solutions with state-of-the-art equipment and our talented engineering staff. Our products are designed and manufactured to consistently meet high quality and performance standards. We are committed to doing our best every day to provide the products and services needed to assist you in delivering the best patient care possible.

SPONGE COUNTER BAGS

Catalog #: MSD-1034-SC

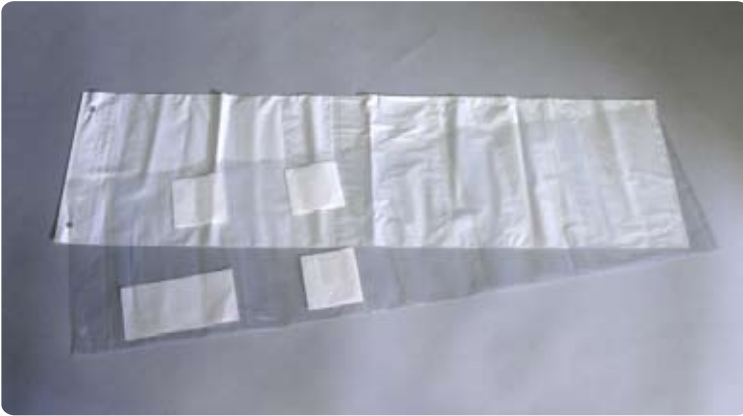

Count: Packaged 50/Dispenser Box/Case; 250/Case
Dim: 10 1/4" X 34" + 2 1/4" HDR

Prevention is the key in today's operating environment; it is all about reducing risk for patients and healthcare workers. What better way then to use the sponge counter bags?

- **Safe:** Minimizes exposure to blood borne pathogens, reduces mess and hazards and help meet OSHA requirements.
- **Convenient:** Counting sponges is easy and efficient. Simply roll up for easy disposal.
- **Economical:** Low cost per procedure. The sponge is easy to see, aiding in monitoring blood loss in patients
- Available with clear or opaque background.

HEART MONITOR BAGS

Catalog #: MSD-0507-HM



Count: Packaged 250 / Case.
Unit Dim: 5 X 7 3/4

Embossed protective heart monitor bags. Each of these bags are constructed of FDA approved 4 mil embossed material.

- Complete system with zipper closure for easy access
- One 2 inch slit in bag for feeding through monitor wires.
- One half inch slit in 8 mil extension from zipper closure to secure gator clip and two quarter inch hang holes for securing bag to patient.
- Prevents damage to monitor from moisture and spills.

COM-PAC ICE BAGS



Guaranteed liquid tight zipper bags are ideal for making ice packs. The double track zipper provides a leakproof seal and the blended polyethylene material remains flexible when cold. Latex free.

Follow these steps to make an ice pack:

- Fill approximately one half of the bag with crushed or cubed ice.
- Add enough water to barely cover the ice.
- Squeeze the air out of the bag.
- Close the bag by sliding your thumb and forefinger across the zippers from one edge to the other.
- Repeat to assure the bag is sealed.
- Wrap the bag in a damp towel, and apply

Part Number	Size	Qty / Pack	Qty / Case
4752-PE	3" x 6"	250	1000
4753-PE	6" x 6"	250	1000
4754-PE	7" x 8"	250	1000
4755-PE	9" x 12"	250	1000
4756-PE	12" x 12"	250	1000
4758-PE	16" x 16"	N/A	100
4759-PE	24" x 24"	N/A	100

TRAY COVERS

Catalog #: MSD-2036-6



Clear disposable plastic sleeves for instrument trays. These sleeves are constructed of 2.2mil FDA approved material. Designed with a five inch Cuff, labeled for hand placement to allow easy application and disposal.


- Prevents cross contamination.
- Simply remove after each patient and dispose.
- Packaged for cleanliness and inventory control.

Count: Packaged 50/inner pack; 200 bags/case

Dim: 20" x 36" bag

CUSTOM PRODUCTS

Com-Pac Custom Medical / Surgical Disposables

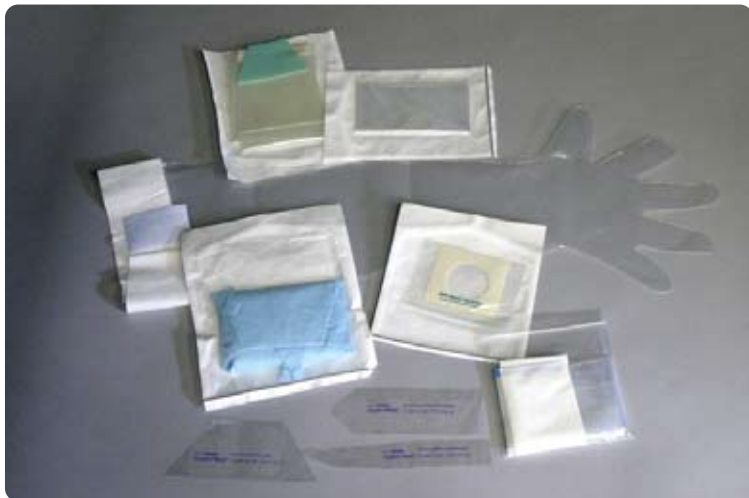

Com-Pac International takes pride in our ability to help our customer design, manufacture, and distribute custom flexible packaging. Our R&D team with over 100 years of experience specializes in design, manufacturing techniques and testing methods to meet the demanding performance requirements for each of our customers unique product in the medical, food and industrial markets. We empower our work force with the resources to exceed our customer expectations for service, product quality and added value.

- Chair Covers
- Laproscope Covers
- X-ray Equipment Covers
- Instrument Barrier Shields
- Telescoping Drapes
- Light Sleeves
- Drainage Bags
- Liquid Tight Oral Feeding Apparatus
- Kit Assembly
- Chemotherapy Bags
- Disposal Systems
- Adhesive Document Protective Covers
- Protective Boots and Gloves
- and much more....

A decorative border with a wavy, ribbon-like appearance. It features horizontal stripes of red, white, and blue, with white stars scattered throughout. The border is set against a solid dark red background.

MADE IN THE USA

WWW.COM-PAC.COM

PH. 800-824-0817

FAX. 618-529-2234