

Bard Access Systems Hickman® Catheters

The most requested name, the most complete line.



BARD

Original Hickman® excellence for all your long-term vascular access needs.


Bard understands that your patients have a variety of vascular access needs. To help you meet these needs, we offer the most complete line of vascular access devices for adult and pediatric care. Our Hickman, Broviac® and

Leonard® Catheters are available with percutaneous introducer kits to facilitate placement. And all these catheters feature the VitaCuff® Antimicrobial Cuff, which helps provide increased protection against infection.

Hickman Triple-Lumen Catheters provide versatile long-term vascular access capabilities.

Hickman and Broviac Single-Lumen Catheters are available in sizes from 2.7 Fr to 9.6 Fr.

Time-saving repair kits for Hickman, Broviac and Leonard Catheters include the external catheter replacement segment, splice sleeve, splice connector, medical adhesive, syringe, 18-gauge blunt needle and instructions.

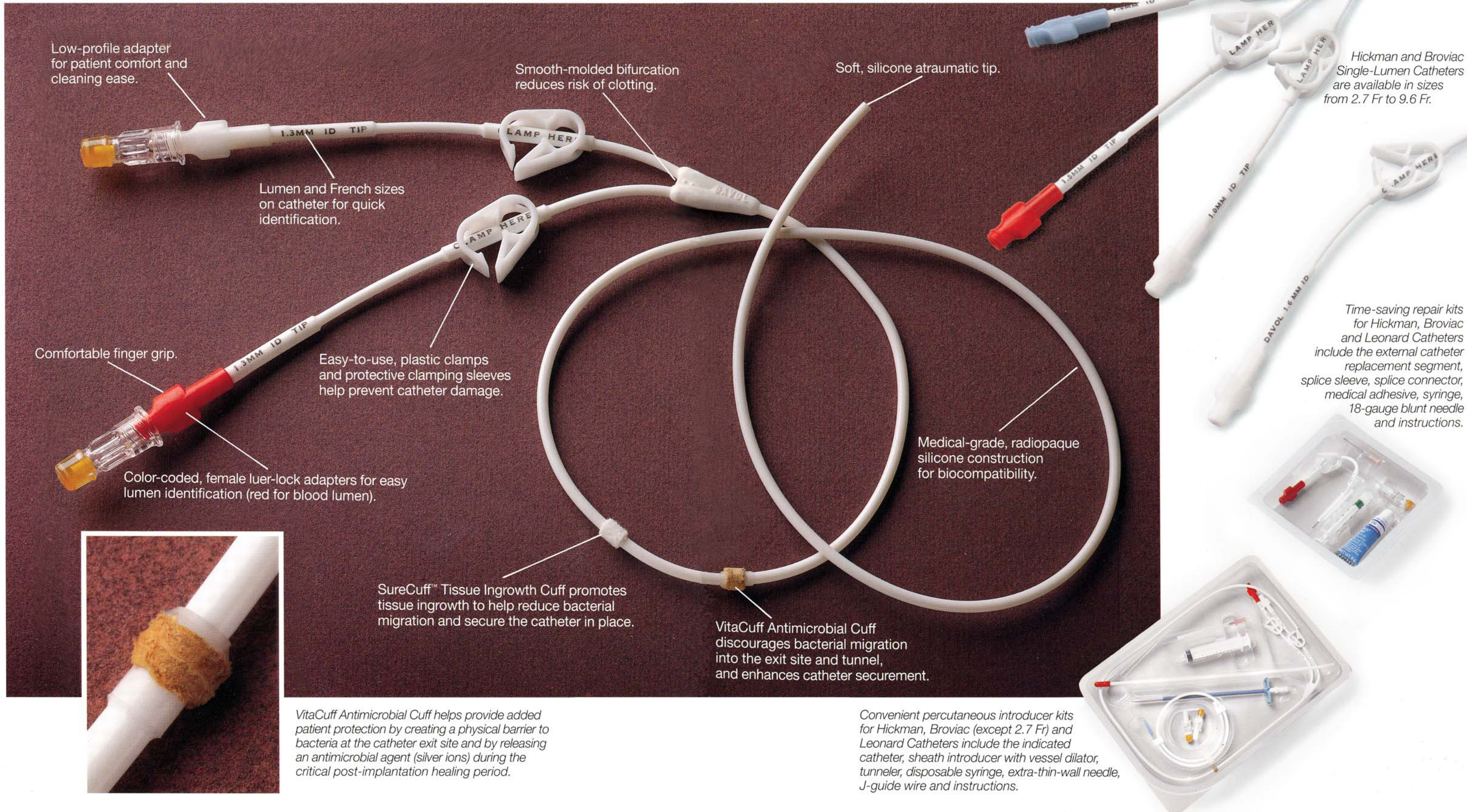

Low-profile adapter for patient comfort and cleaning ease.

Smooth-molded bifurcation reduces risk of clotting.

Soft, silicone atraumatic tip.

Lumen and French sizes on catheter for quick identification.

Comfortable finger grip.

Easy-to-use, plastic clamps and protective clamping sleeves help prevent catheter damage.

Medical-grade, radiopaque silicone construction for biocompatibility.

Color-coded, female luer-lock adapters for easy lumen identification (red for blood lumen).

SureCuff™ Tissue Ingrowth Cuff promotes tissue ingrowth to help reduce bacterial migration and secure the catheter in place.


VitaCuff Antimicrobial Cuff discourages bacterial migration into the exit site and tunnel, and enhances catheter securement.

VitaCuff Antimicrobial Cuff helps provide added patient protection by creating a physical barrier to bacteria at the catheter exit site and by releasing an antimicrobial agent (silver ions) during the critical post-implantation healing period.

Convenient percutaneous introducer kits for Hickman, Broviac (except 2.7 Fr) and Leonard Catheters include the indicated catheter, sheath introducer with vessel dilator, tunneler, disposable syringe, extra-thin-wall needle, J-guide wire and instructions.

Hickman®, Broviac® and Leonard® Catheters with VitaCuff® Antimicrobial Cuff

Catheter code numbers ending in "4" have VitaCuff Antimicrobial Cuff positioned immediately adjacent to the SureCuff™ Tissue Ingrowth Cuff; catheter code numbers ending in "2" have VitaCuff positioned approximately 4 cm from SureCuff. Catheter code numbers ending in "0" have SureCuff only.

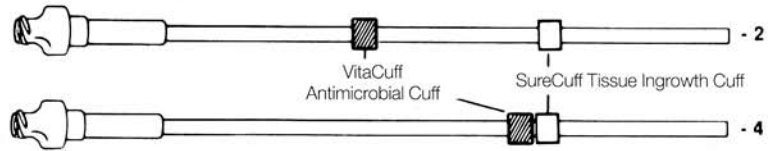

Ordering Information

Single-Lumen

	Intro-Eze™ Introducer			Peel-Apart Introducer			Cutdown			Total Length	Volume	Lumen Size (OD/ID)	Repair Kit
2.7 Fr Broviac	NA	NA	NA	NA	NA	NA	0600044	0600042	0600040	71 cm	0.15 ml	0.9 mm/ 0.5 mm	0601600
4.2 Fr Broviac	0603674	0603672	0603670	0600524	0600522	0600520	0600064	0600062	0600060	71 cm	0.3 ml	1.4 mm/ 0.7 mm	0601610
6.6 Fr Broviac	0603654	0603652	0603650	0600544	0600542	0600540	0600104	0600102	0600100	90 cm	0.7 ml	2.2 mm/ 1.0 mm	0601620
6.6 Fr Broviac Short Length	NA	NA	NA	NA	NA	NA	0600124	NA	0600120	90 cm	0.7 ml	2.2 mm/ 1.0 mm	0601620
9.6 Fr Hickman	0603114	NA	0603110	0600564	0600562	0600560	0600164	0600162	0600160	90 cm	1.8 ml	3.2 mm/ 1.6 mm	0601630

Dual-Lumen

	Intro-Eze™ Introducer			Peel-Apart Introducer			Cutdown			Total Length	Volume	Lumen Size (OD/ID)	Repair Kit
7 Fr Hickman	NA	0603712	0603710	0600574	0600572	0600570	NA	0600312	0600310	65 cm	0.8 ml red, 0.6 ml white	2.3 mm/ 1.0 mm red/ 0.8 mm white	0601690 red leg, 0601680 white leg
9 Fr Hickman Pediatric	NA	NA	0600010	NA	NA	0600580	0600324	NA	0600320	65 cm	1.3 ml red, 0.6 ml white	3.0 mm/ 1.3 mm red/ 0.7 mm white	0601700 body, 0601690 red leg, 0601680 white leg
9 Fr Hickman	0600024	NA	0600020	0600604	NA	0600600	0600334	NA	0600330	90 cm	1.3 ml red, 0.6 ml white	3.0 mm/ 1.3 mm red/ 0.7 mm white	0601700 body, 0601690 red leg, 0601680 white leg
10 Fr Leonard	0603104	0603102	0603100	0600634	0600632	0600630	0600344	0600342	0600340	90 cm	1.3 ml red, 1.3 ml white	3.3 mm/ 1.3 mm red/ 1.3 mm white	0601750 body, 0601690 red leg, 0601680 white leg
12 Fr Hickman	0600034	0600032	0600030	0600624	0600622	0600620	NA	0600352	0600350	90 cm	1.8 ml red, 1.8 ml white	4.0 mm/ 1.6 mm red/ 1.6 mm white	0601710 body, 0601690 red leg, 0601680 white leg

Triple-Lumen

	Intro-Eze™ Introducer			Peel-Apart Introducer			Cutdown			Total Length	Volume	Lumen Size (OD/ID)	Repair Kit
10 Fr Hickman	0606664	NA	0606660	0606564	NA	0606560	0606464	NA	0606460	97 cm	1.4 ml red, 0.8 ml white, 0.8 ml blue	3.3 mm/ 1.5 mm red/ 0.8 mm white/ 0.8 mm blue	0601790 body, 0601690 red leg, 0601680 white leg, 0601730 blue leg
12.5 Fr Hickman	0600054	0600052	0600050	0600654	0600652	0600650	NA	0600362	0600360	90 cm	1.6 ml red, 0.7 ml white, 0.7 ml blue	4.2 mm/ 1.5 mm red/ 1.0 mm white/ 1.0 mm blue	0601740 body, 0601690 red leg, 0601680 white leg, 0601730 blue leg

Adhesive Kit

Product Code 60172

Injection Caps

Product Code 60186

Bard Access Systems

Hickman. Groshong. Designs for Life.™

In the United States:

Bard Access Systems
5425 West Amelia Earhart Drive, Salt Lake City, Utah 84116
801-595-0700

For Customer Service: 800-545-0890
For Clinical Information: 800-443-3385

In Canada:

Bard Canada
2345 Stanfield Road, Mississauga, Ontario L4Y3Y3

For Customer Service: Western Canada: 800-268-2862
Eastern Canada: 800-387-7851
Ontario: 800-387-9473

BARD, Hickman, Broviac, Leonard and Groshong are registered trademarks of C.R. Bard, Inc. or an affiliate.
VitaCuff is a registered trademark of Vitaphore Corporation.
SureCuff and Intro-Eze are trademarks of C.R. Bard Inc. or an affiliate.

©1994 Bard Access Systems

BAS29094
129410M
Printed in U.S.A.