



2015 Product Reference Guide

Advancing Wound Care



TO OUR CUSTOMERS

Tissue Regeneration in Wound Care Across the Spectrum

Derma Sciences is a Tissue Regeneration Company at the forefront of research and innovation for the management of acute and chronic wounds, and burns. The products we bring to the market are in the best tradition of evidence-based medicine, helping clinicians around the world improve clinical outcomes. We are proud to bring you these innovations and are committed to your cause.

With you in mind, our portfolio has grown to provide clinicians with solutions to assist with tissue regeneration - from the start of a wound through to wound closure. From MEDIHONEY[®], to XTRASORB[®], to BIOGUARD[®], ALGICELL[®] Ag, TCC-EZ[®] and now, AMNIOEXCEL[®] and AMNIOMATRIX[®], our intelligent, evidence-based products help to advance healing for patients suffering from complex chronic or acute wounds such as diabetic foot ulcers, venous leg ulcers, pressure ulcers and burns.

MEDIHONEY®

MEDIHONEY®, with Active *Leptospermum* Honey (ALH), is the global leading medical-grade honey product for the management of wounds and burns. Derived from the Active *Leptospermum scoparium* species of plant in New Zealand, these unique dressings contain plant-derived properties beneficial throughout all phases of wound healing and is the most studied species of honey for the management of wounds and burns. Due to its multiple mechanisms of action (low pH and high osmolarity), MEDIHONEY® is a first-line choice to help in the management of chronic and stalled wounds and burns, and to assist with fast, effective autolytic debridement of necrotic tissue.

XTRASORB®

Derma Sciences' XTRASORB® dressings, with a patented Super Absorbent Polymer (SAP) technology, provide superior absorption, fluid handling and moisture management capabilities ideal moist wound healing environment in the wound bed and surrounding periwound tissue.

Absorbs 2x to 4x as much fluid as conventional dressings; Converts fluid to a gel – locking it within the dressing and away from the wound bed – even under compression; Keeps harmful components, like MMPs and bacteria locked in the gel; Helps to reduce the risk of maceration; Increases time between dressing changes, minimizes wound bed disruption and patient discomfort, and decreases costs; Multiple configurations to handle any type of wound exudate.

BIOGUARD®

BIOGUARD® Barrier Dressings, with a cationic biocide, PolyDADMAC, intrinsically bound to the substrate dressing, provide a >5 log-kill of bacteria within the dressing and without toxicity to healthy tissue. Due to its patented manufacturing process, BIOGUARD® Barrier Dressings are the only dressing on the market that are non-leaching, non-resistant and non-toxic.

Prevents bacteria from entering the wound area and inhibits its growth; Its active molecule, PolyDADMAC, is a cationic biocide that attracts bacteria and binds with it, disrupting and destroying the cell wall; PolyDADMAC is intrinsically bound to the substrate, so it doesn't leach away from the dressing (leaching antimicrobials can delay wound healing and pose resistance issues; PolyDADMAC is up to 100 times larger than PHMB (used in competitive dressings)

(the larger the molecule the lower the chances of creating resistance even over long exposure); BIOGUARD is non-toxic to healthy cells.

ALGICELL® Ag

Silver has been used as an antimicrobial since ancient times, and has undergone a resurgence in the last decade. The effectiveness of a silver dressing is influenced by many factors: the type of silver, the length of its antimicrobial activity, its absorptive capabilities and more. ALGICELL® Ag addresses these key issues powerfully and is designed for superiority in the management of critically colonized or infected wounds.

TCC-EZ® AND MEDE-KAST®

Derma Sciences is the only company to give you the flexibility of the quick and easy TCC-EZ® casting system and a traditional Total Contact Cast as well.

Derma Sciences TCC-EZ® offers a one-piece, roll-on, woven design that simplifies the application process while reducing the potential for causing additional tissue damage. The result is a process that takes about half the amount of time of traditional casting systems and a product that gives greater patient comfort. All while providing evidence-based off-loading for optimal clinical outcomes.

The traditional MedE-Kast® and MedE-Kast® Ultra Total Contact Cast Systems provide a healthcare provider with the necessary components to apply a total contact cast in one simple kit. The MedE-Kast® Ultra Total Contact Cast System contains two pre-made splints to help make the application process easier and faster.

AMNIOEXCEL® AND AMNIOMATRIX®

AMNIOEXCEL® is an Amniotic Allograft Membrane composed of an epithelial layer, basement membrane and compact layer. Dehydrated using our patent pending DryFlex® processing technology for a highly malleable and conformable product providing intimate contact with the wound bed and ease of application. This extracellular matrix is fully resorbable and provides a natural environment to enhance soft tissue reconstruction & regeneration.

AMNIOMATRIX® is an Amniotic Allograft Suspension derived from the components of the amniotic fluid and morselized amnion. It is cryopreserved and the proprietary processing method preserves the viable cellular materials, cytokines, and growth factors.

Table of Contents

Advanced Wound Care

MEDIHONEY® Wound & Burn Dressings 5
AMNIOEXCEL® 8
AMNIOMATRIX® 8
TCC-EZ® Total Contact Casting System 9
MEDE-KAST® AND MEDE-KAST® ULTRA 10
XTRASORB® Super Absorbent Dressings 11
BIOGUARD® Barrier Dressings 13
ALGICELL® Ag 14
Calcium Alginate Dressings 15
Silver Antimicrobial Wound Care Dressings 16
Dermagran® Advanced Wound Care Products 17
Hydrocolloid Dressings 19
Hydrogel Dressings 21
Foam Dressings 23
Wound Cleansers 25

Traditional Wound Care

Sponges 29
Wound Dressings 32
Packing Strips 35
Island and Barrier Border Gauze Dressings 36
Non-Adherent Impregnated Dressings 37
Conforming Bandages 39
Elastic Net Dressings (Retainer / Tube Gauze) 41
Compression Bandages and Bandaging Systems 49
Protective Products 55
Tapes 57
Laparotomy Sponges 58

Burn Care

Silver Antimicrobial Burn Care Dressings 61
Gauze Burn Dressings 63
Other Burn Care Products 64

Skin Care and Bathing

Incontinence Products 67
Hand Sanitizers / Lotions 68
Shampoos / Body Washes 69
Bath Additives 70

Specialty Fixation and Closure Devices

Fixation Devices 74
Wound Closure Strips 76



Advanced Wound Care

- **MEDIHONEY® Wound & Burn Dressings**
- **AMNIOEXCEL® Amniotic Allograft Membrane**
- **AMNIOMATRIX® Amniotic Allograft Suspension**
- **TCC-EZ® Total Contact Casting System**
- **XTRASORB® Super Absorbent Dressings**
- **BIOGUARD® Barrier Dressings**
- **ALGICELL® Ag Dressings**
- **Calcium Alginate Dressings**
- **Silver Antimicrobial Wound Care Dressings**
- **DERMAGRAN® Advanced Wound Care Products**
- **Hydrocolloid Dressings**
- **Hydrogel Dressings**
- **Foam Dressings**
- **Wound Cleansers**

MEDIHONEY® Dressings

MEDIHONEY® Gel Dressing

MEDIHONEY® Gel, a patented formulation with 80% Active *Leptospermum* Honey and 20% natural gelling agents, has improved physical stability, allowing the product to remain at the site of the wound – even in the presence of wound fluid and body heat – for an extended time. This results in easier application for caregivers, and offers patients an enhanced user experience with improved clinical outcomes. **Does not contain natural rubber latex.**



Product	Description	Pkg unit/Case	HCPCS
31805	0.5 oz tube	10/box, 4 boxes/case	A6240
31815	1.5 oz tube	1/box, 12 boxes/case	A6240

MEDIHONEY® Paste Dressing

MEDIHONEY® is a special type of honey indigenous to New Zealand, called Active *Leptospermum* (or Manuka) Honey. This honey has been shown to have unique properties that are advantageous to wound healing, including a low pH and high osmolarity. MEDIHONEY® is indicated for use on lightly to moderately exuding wounds, or wounds that are difficult to manage due to their odd shape or placement on the body. This non-toxic dressing can be used as a cornerstone for wound bed preparation, and throughout all the phases of wound healing. **Does not contain natural rubber latex.**



Product	Description	Pkg unit/Case	HCPCS
31505	0.5 oz tube	10/box, 4/case	N/A
31515	1.5 oz tube	1/box, 12/case	N/A
31535	3.5 oz tube	1/box, 12/case	N/A

MEDIHONEY® HCS (Hydrogel Colloidal Sheet) Dressing

MEDIHONEY® HCS may be our most versatile dressing yet. Combining the best of two novel technologies, MEDIHONEY® HCS is one new and unique dressing with an increased broad spectrum of application and use. *Does not contain natural rubber latex.*



Product	Description	Pkg unit/Case	HCPCS
NON-ADHESIVE			
31622	2.4" x 2.4"	10/box, 5 boxes/case	A6234
31644	4.3" x 4.3"	10/box, 5 boxes/case	A6235
31688	8" x 8"	5/box, 4 boxes/case	N/A
31612	8" x 12"	2/box, 5 boxes/case	N/A
FENESTRATED - NON-ADHESIVE			
31618	1.8" x 1.8"	10/box, 5 boxes/case	A6234
ADHESIVE			
31722	2.8" x 2.8" (4.3" x 4.3" with adhesive border)	10/box, 5 boxes/case	A6237
31744	4½" x 4½" (6" x 6" with adhesive border)	10/box, 5 boxes/case	A6238

MEDIHONEY® Dressings

MEDIHONEY® Honeycolloid™ Dressing

Dressing is a novel honey-based dressing for lightly to moderately exuding wounds. The dressing contains 80% Active *Leptospermum* Honey, can be used to cover a shallow wound or fill a deep wound (non-adhesive version only), and provides a moist environment conducive to wound healing. The adhesive version is specifically to be used as an all-in-one dressing, with no need for a secondary dressing to hold it in place. This non-toxic dressing can be used as a cornerstone for wound bed preparation, and throughout all the phases of wound healing. **Does not contain natural rubber latex.**



Product	Description	Pkg unit/Case	HCPCS
NON-ADHESIVE			
31222	2" x 2"	10/box, 10 boxes/case	A6234
31245	4" x 5"	10/box, 5 boxes/case e	A6235
ADHESIVE			
31422	2" x 2" (3½" x 3½" with adhesive border)	10/box, 10 boxes/case	A6237
31445	4½" x 4½" (6" x 6" with adhesive border)	10/box, 5 boxes/case	A6238

MEDIHONEY® Calcium Alginate Dressing


MEDIHONEY® Calcium Alginate Dressing is a novel honey-based dressing for moderately to heavily exuding wounds. The dressing contains 95% Active *Leptospermum* Honey and can be used to cover a shallow wound or fill a deep wound, while providing a moist environment conducive to wound healing. *Leptospermum* Honey has been shown to have several properties which assist in promoting an optimal healing environment. It has a low pH, and has been shown to help lower the overall wound pH, which is beneficial to chronic wounds. **Does not contain natural rubber latex.**



Product	Description	Pkg unit/Case	HCPCS
31012	¾" x 12"	5/box, 4 boxes/case	A6199
31022	2" x 2"	10/box, 10 boxes/case	A6196
31045	4" x 5"	10/box, 5 boxes/case	A6197

AMNIOEXCEL®

AMNIOEXCEL® is an Amniotic Allograft Membrane composed of an epithelial layer, basement membrane and compact layer. Dehydrated using our patent pending DryFlex® processing technology for a highly malleable and conformable product providing intimate contact with the wound bed and ease of application. This extracellular matrix is fully resorbable and provides a natural environment to enhance soft tissue reconstruction & regeneration. *Does not contain natural rubber latex.*

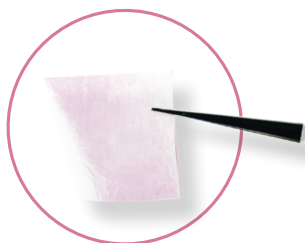

Product	Description	Pkg unit/Case	HCPCS
25015	1.5cm x 1.5cm (2.25cm ²)	1/box	Q4137
25023	2cm x 3cm (6.0cm ²)	1/box	Q4137
25044	4cm x 4cm (16.0cm ²)	1/box	Q4137
25048	4cm x 8cm (32.0cm ²)	1/box	Q4137

AMNIOMATRIX®

AMNIOMATRIX® is an Amniotic Allograft Suspension derived from the components of the amniotic fluid and morselized amnion. It is cryopreserved and the proprietary processing method preserves the viable cellular materials, cytokines, and growth factors. *Does not contain natural rubber latex.*



Product	Description	Pkg unit/Case	HCPCS
25510	0.25ml	1 vial	Q4139
25520	0.50ml	1 vial	Q4139
25530	1.0ml	1 vial	Q4139
25540	2.0ml	1 vial	Q4139
25550	3.0ml	1 vial	Q4139

TCC-EZ® TOTAL CONTACT CASTING SYSTEM

TCC-EZ® is a novel total contact casting system that replaces several rolls with one simple, easy to apply casting component. With an open weave, this fiberglass component can be applied in under 5 minutes, significantly reducing the time required to apply a total contact cast, while preserving the off-loading characteristics of a traditional casting system. When in place, TCC-EZ® effectively reduces (or "off loads") pressure exerted onto diabetic foot ulcers as patients stand and walk. Known as the Gold Standard for off-loading pressure, use of total contact casts as a key component in diabetic foot ulcer care has been shown in several randomized controlled studies to result in an ulcer heal rate of up to 89%, although healing rates and times varied.



Product	Description	Pkg Unit/Case	HCPCS
TCC-EZ CASTING SYSTEMS			
TCC23000	10 Casting Systems 3" with 2 Regular Boots	1/case	Q4038
TCC23001	10 Casting Systems 3" with No Boot	1/case	Q4038
TCC23002	Casting System 3" Single Application	1/case	Q4038
TCC23005	5 Casting System 3" with 1 Regular Boot	1/case	Q4038
TCC23051	5 Casting System 3" with No Boot	1/case	Q4038
TCC24000	10 Casting Systems 4" with 2 Regular Boots	1/case	Q4038
TCC24001	10 Casting Systems 4" with No Boot	1/case	Q4038
TCC24002	Casting System 4" Single Application	1/case	Q4038
TCC24005	5 Casting System 4" with 1 Regular Boot	1/case	Q4038
TCC24014	10 Casting Systems 4" with 1 Large Boot	1/case	Q4038
TCC24016	10 Casting Systems 4" with 1 Large Charcot Boot	1/case	Q4038
TCC24024	10 Casting Systems 4" with 1 Extra Large Boot	1/case	Q4038
TCC24026	10 Casting Systems 4" with 1 Extra Large Charcot Boot	1/case	Q4038
TCC24051	5 Casting System 4" with No Boot	1/case	Q4038
TCC24214	10 Casting Systems 4" with 2 Large Boots	1/case	Q4038
TCC24224	10 Casting Systems 4" with 2 Extra Large Boots	1/case	Q4038
TCC25050	5 – 3" and 5 – 4" with 2 Regular Boots	1/case	Q4038
TCC25051	5 – 3" and 5 – 4" with No Boot	1/case	Q4038
TCC-EZ BOOTS			
TCC21100	Regular Boot	1/case	N/A
TCC21114	Large Boot	1/case	N/A
TCC21116	Large Charcot Boot	1/case	N/A
TCC21124	Extra Large Boot	1/case	N/A
TCC21126	Extra Large Charcot Boot	1/case	N/A
TCC-EZ TOOLS			
TCC2SAW	Cast Removal Saw	1/case	N/A
TCC2SAWSSB	Saw Replacement Stainless Steel Blades (4 Blades)	1/case	N/A
TCC2VACFLT	Dust Vacuum Disposable Filter Cartridge	1/case	N/A
TCC2VAC	Cast Removal Dust Vacuum	1/case	N/A
TCC2SPRL	Cast Removal Spreader Large	1/case	N/A
TCC ADDITIONAL ITEMS			
TCC2FCT03	3" Rolls of Fiberglass (10 rolls)		N/A
TCC2FCT04	4" Rolls of Fiberglass (10 rolls)		N/A
TCC2FC045	4" x 5 yds Rolls of Plaster (12 rolls)		N/A

MEDE-KAST® AND MEDE-KAST® ULTRA TOTAL CONTACT CASTING SYSTEMS

The traditional MedE-Kast® and MedE-Kast® Ultra Total Contact Cast Systems provide a healthcare provider with the necessary components to apply a total contact cast in one simple kit. The MedE-Kast® Ultra Total Contact Cast System contains two pre-made splints to help make the application process easier and faster.



Product	Description	Pkg Unit/Case	HCPCS
TCC2MDKK	10 Casting Systems	1/case	Q4038
TCC2MDKKS	Casting System - Single Application	1/case	Q4038
TCC2ULTR	Ultra – 10 Casting Systems	1/case	Q4038
TCC2ULTRS	Ultra Casting System – Single Application	1/case	Q4038

XTRASORB® DRESSINGS

XTRASORB®

Super Absorbent Dressing

XTRASORB® is a novel super-absorbent wound dressing. At its core is a super-absorbent polymer that allows the dressing to absorb significantly increased amounts of fluid as compared to that of typical highly absorbent dressings. Whereas other absorbent dressings (including foams, ABD pads, and other cellulose based dressings) rely on spaces or cavities within the dressings to fill with fluid, XTRASORB®'s super absorbent polymer converts fluid within the dressing into a gel. This keeps the fluid bound within the dressing, and away from the patient's skin. As such, the dressing is designed to decrease dressing change frequency, to reduce the risk of maceration, and to protect the patient's skin from harmful components within wound fluid that can impede wound healing. *Does not contain natural rubber latex.*



Product	Description	Pkg unit/Case	HCPCS
89545	4" x 5" Non-Adhesive	10/box, 10 boxes/case	A6252
89569	6" x 9" Non-Adhesive	10/box, 5 boxes/case	A6253

XTRASORB® Foam

Super Absorbent Dressing

XTRASORB® Foam is a novel foam-based super-absorbent wound dressing. The dressing is comprised of two layers: The wound contact layer is a standard polyurethane hydrophilic foam. Onto this is laminated the second layer, a patented super-absorbent copolymer sheet. When fluid enters the dressing via the foam wound contact layer, it is quickly absorbed up into the super-absorbent layer. This converts the fluid into a gel, allowing it to be bound into the dressing, which can help to reduce the risk of maceration while also keeping the harmful components of wound fluid away from the patient's skin. The dressings can absorb a significantly increased amount of wound fluid as compared to that of standard foam dressings. XTRASORB® Foam comes in both adhesive and non-adhesive formats. *Does not contain natural rubber latex.*



Product	Description	Pkg unit/Case	HCPCS
NON-ADHESIVE			
86122	2" x 2"	10/box, 4 boxes/case	A6209
86144	4" x 4¾"	10/box, 4 boxes/case	A6210
86188	8" x 8"	5/box, 8 boxes/case	A6211
ADHESIVE			
86233	3.2" x 3.2"	10/box, 4 boxes/case	A6212
86244	4½" x 4½"	10/box, 4 boxes/case	A6212
86266	6" x 6"	10/box, 4 boxes/case	A6213

XTRASORB® HCS (Hydrogel Colloidal Sheet) Super Absorbent Dressing

XTRASORB® HCS is a novel hydrogel-hydrocolloid combination super-absorbent wound dressing. The clear dressing takes a low water content hydrogel and combines it with sodium carboxymethyl cellulose (sodium CMC – a common component in hydrocolloids and the base component in hydrofibers), along with other gelling/absorbing agents. The result is a dressing that can be used to help hydrate non-exuding wounds, while also being able to absorb fluid from lightly to moderately exuding wounds. The fluid that is absorbed converts to a gel, allowing it to be bound into the dressing, which can help to reduce the risk of maceration while also keeping the harmful components of wound fluid away from the patient's skin. The dressing can absorb significantly more wound fluid as compared to that of standard hydrocolloid dressings. XTRASORB® HCS comes in both adhesive and non-adhesive formats. **Does not contain natural rubber latex.**



Product	Description	Pkg unit/Case	HCPCS
NON-ADHESIVE			
86322	2.3" x 2.3"	10/box, 4 boxes/case	A6234
86344	4.3" x 4.3"	10/box, 4 boxes/case	A6235
86388	8" x 8"	5/box, 4 boxes/case	A6236
ADHESIVE			
86433	3" x 3"	10/box, 4 boxes/case	A6237
86466	6" x 6"	10/box, 4 boxes/case	A6238

BIOGUARD® BARRIER DRESSINGS

(Compare to Covidien AMD® dressings)

BIOGUARD® Barrier Dressings, with a cationic biocide, PolyDADMAC, intrinsically bound to the substrate dressing, provide a >5 log-kill of bacteria within the dressing and without toxicity to healthy tissue. Due to its patented manufacturing process, BIOGUARD® Barrier Dressings are the only dressing on the market that are non-leaching, non-resistant and non-toxic. **Does not contain natural rubber latex.**

Product	Description	Pkg Unit/Case	HCPCS
CONFORMING BANDAGE			
97241	2" x 4.1 yds	1 roll/pkg, 100 pkgs/case	A6266
97341	3" x 4.1 yds	1 roll/pkg, 100 pkgs/case	A6266
97441	4" x 4.1 yds	1 roll/pkg, 100 pkgs/case	A6266
GAUZE SPONGES			
97412	12-Ply, 4" x 4", 2's	50 dressings/tray, 24 trays/case (1200 dressings)	A6222
97208	8-Ply, 2" x 2", 2's	100 dressings/tray, 30 trays/case (3000 dressings)	A6222
GAUZE ROLL			
97322	6-Ply, 4.5" x 4.1 yds	1 roll/pkg, 100 pkgs/case	A6266
PACKING STRIPS			
97831	¼" x 5 yds	12/cs	A6266
97832	½" x 5 yds	12/cs	A6266
97833	1" x 5 yds	12/cs	A6266
NON-ADHERENT DRESSINGS			
97334	3" x 4"	50/box, 18/cs	A6222
97338	3" x 8"	50/box, 12/cs	A6223
ISLAND DRESSINGS			
97041	4" x 10"	25/box, 4 boxes/cs	A6223
97042	4" x 14"	25/box, 2 boxes/cs	A6223
97045	4" x 5"	25/box, 8 boxes/cs	A6223
97048	4" x 8"	25/box, 4 boxes/cs	A6223
READY – CUT GAUZE BURN DRESSING			
97118	100-Ply, 18" x 18"	1/pk, 10 pkgs/cs	N/A



Conforming Bandage



Sponges



Gauze Roll



Ready-Cut Gauze Burn Dressing



Non-Adherent Dressings



Packing Strips



Island Dressings

ALGICELL® Ag CALCIUM ALGINATE DRESSING WITH ANTIMICROBIAL SILVER

(Comparable to Aquacel® Ag, Silvercel®, and Maxorb® Extra Ag)

ALGICELL® Ag is a primary wound dressing made of Calcium Alginate containing 1.4% silver. In the presence of wound exudate, the sodium ions from the exudate take the place of the silver ions, releasing the silver ions. As wound exudate is absorbed, the alginate forms a gel, which assists in maintaining a moist environment for optimal wound healing, and allows intact removal. The

silver ions released in the presence of wound fluid protect the dressing from bacterial colonization, and provides an effective barrier to bacterial penetration. ALGICELL® Ag dressing provides an antimicrobial effect lasting up to seven days. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
88512	¾" x 12"	5/box, 4 boxes/case	A6199
88522	2" x 2"	10/box, 10 boxes/case	A6196
88544	4 ¼" x 4¼"	10/box, 5 boxes/case	A6197
88545	4" x 5"	10/box, 5 boxes/case	A6197
88548	4" x 8"	5/box, 4 boxes/case	A6197
88582	8" x 12"	5/box, 4 boxes/case	A6198

CALCIUM ALGINATE DRESSINGS

Algicell® Calcium Alginate Dressing – Sterile

(Comparable to Kaltostat®)

These soft, white sterile dressings absorb moderate to large amounts of wound drainage and may help control minor bleeding. Alginate ropes can be used to fill or pack the dead space

in a wound. Algicell® dressings are designed to keep their form, allowing for easy and trauma-free removal, while maintaining a moist wound environment through its ability to gel upon contact with wound fluids. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
88022	2" x 2"	10/box, 4 boxes/case, 40/case	A6196
88044	4" x 4"	10/box, 4 boxes/case, 40/case	A6196
88048	4" x 8"	10/box, 4 boxes/case, 40/case	A6197
88112	¾" x 12", Rope	10/box, 4 boxes/case, 40/case	A6199
88100	¾" x 36", Rope	1/box, 12 boxes/case, 12/case	A6199

Calcium Alginate Dressing – Sterile

These soft, white sterile dressings absorb moderate to large amounts of wound drainage and may help control minor bleeding. Calcium alginate ropes can be used to fill or pack the dead space in a wound. Calcium Alginate dressings are designed to keep their form, allowing for easy and trauma-free removal, while maintaining a moist wound environment through its ability to gel upon contact with wound fluids. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
88222	2" x 2"	10/box, 4 boxes/case, 40/case	A6196
88244	4" x 4"	10/box, 4 boxes/case, 40/case	A6196
88245	4" x 5"	10/box, 4 boxes/case, 40/case	A6197
88248	4" x 8"	10/box, 4 boxes/case, 40/case	A6197
88212	¾" x 12", Rope	10/box, 4 boxes/case, 40/case	A6199

SILVER ANTIMICROBIAL WOUND CARE DRESSINGS

Silverseal® Wound Contact Dressing – Sterile

Silverseal® Wound Contact Dressings are designed for placement directly on the wound. They provide a topical antimicrobial barrier that reduces the bacteria and fungi counts and within 2 to 4 millimeters of the wound surface. The dressing can be left in place for up to seven days. Silverseal® Dressings are composed of knitted fabrics that are flexible, minimally adhering, and antimicrobial. The knitted fabric is 100% covered, uniformly and circumferentially with a thin coat of silver approximately (4/100,000 inch). The circumferential plating of the fibers results in a three-dimensional surface area of pure metallic silver in contact with the wound fluids. Analysis of the surface of Silverseal® reveals 99.1% pure metallic silver and 0.9% silver oxide. These dressings are versatile and user friendly, conforming to hard-to-dress areas and providing mild compression when needed. Silverseal® Wound Contact Dressings are indicated for management of incisions, skin grafts, donor sites, lacerations, abrasions, pressure sores, chronic wounds, chronic ulcers, dermal ulcers, vascular ulcers, and diabetic ulcers.



Product	Description	Pkg Unit/Case	HCPCS
SS-WCD-1012	10" x 12"	5/box	A6208
SS-WCD-0412	4" x 12"	10/box	A6208
SS-WCD-0404	4¼" x 4¼"	10/box	A6207
SS-WCD-0466	4" x 66"	1/box	A6208

DERMAGRAN® ADVANCED WOUND CARE PRODUCTS

Dermagran® Ointment

Dermagran® Ointment is a pH balanced skin protectant. Formulated with Derma Sciences' proprietary zinc-nutrient technology, usage of Dermagran® Ointment provides a moist, mildly acidic environment conducive to wound healing. In addition, this product contains Vitamin A, Calcium, and Magnesium. Indicated for the temporary protection, comfort, and lubrication of minor skin irritations such as intertrigo, chafing, galling, rubbing or friction. *Does not contain natural rubber latex.*



Product	Description	Pkg Unit/Case	HCPCS
DT-4	4 oz tube	12/case	A6250
DG-4	4 oz jar	12/case	A6250

Dermagran®-B Hydrophilic Wound Dressings

Dermagran®-B Hydrophilic Wound Dressings provides a primary cover, filler or both for wound deficiencies, which absorbs wound exudate and creates a moist, mildly acidic environment to facilitate wound healing. Available as an impregnated gauze or tube, both contain Derma Sciences' zinc-nutrient dressing formulation and balanced pH technology. The dressing provides a moist, mildly acidic environment conducive to wound healing. In addition these products contain Vitamin A, Vitamin B₆, Calcium and Magnesium. Indicated for use on skin ulcers (diabetic, venous stasis), stage I-IV pressure sores, surgical incisions and superficial injuries such as minor partial thickness thermal burns, superficial lacerations, cuts or abrasions. *Does not contain natural rubber latex.*



Product	Description	Pkg Unit/Case	HCPCS
SPD-03	3 oz tube	12/case	A6261
IMPREGNATED GAUZE – STERILE			
SPD-20	2" x 2"	25/box, 8 boxes/case, 200/case	A6231
SPD-21	4" x 4"	15/box, 4 boxes/case, 60/case	A6231
SPD-24	8" x 4"	15/box, 4 boxes/case, 60/case	A6232

Dermagran® Moisturizing Spray

Dermagran® moisturizing spray is a pH balanced skin protectant. It conditions the skin before application of Dermagran® Ointment or another appropriate dressing. Indicated for on stage I pressure ulcers, abrasions, and preventive skin care. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
DM-4	4 oz spray	12/case	A6250

Dermagran® Amorphous Hydrogel Dressing

(Comparable to Intrasite®)

Our Amorphous Hydrogel is part of our Dermagran® family, and contains the same zinc-nutrient technology as the other products in the line. Hydrogel Wound Dressing provides a primary cover, filler or both for wound deficiencies. It absorbs wound exudate and creates a moist, mildly acidic environment to facilitate wound healing. Indicated for use on stage IV pressure ulcers, venous stasis ulcers, partial thickness thermal burns, surgical incisions, skin irritations, abrasions, cuts and skin conditions associated with peristomal care. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
DH-03	3 oz tube	12/case	N/A

HYDROCOLLOID DRESSINGS

Primacol® Hydrocolloid Dressing – Sterile

(Comparable to Duoderm CGF®)

Primacol® dressings mix with wound exudate to form a soft mass that allows for easy removal while protecting the wound from outside contamination from bacteria, fecal matter or urine. A moist wound environment is created, conducive to wound healing. The foam backing provides an extra layer of cushioning and comfort for patients, along with tapered edges. The tapered edges are designed to minimize bunching and rolling up, leading to improved wear times and cost-effectiveness. Indicated for use for the management of lightly exuding Pressure Ulcers (Stage I-IV), Diabetic foot ulcers, venous stasis ulcers, other chronic and acute wounds. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
85144	4" x 4"	5/box, 20 boxes/case, 100/case	A6234
85166	6" x 6"	5/box, 16 boxes/case, 80/case	A6235
85188	8" x 8"	5/box, 16 boxes/case, 80/case	A6236

Primacol® Thin Hydrocolloid Dressing – Sterile

(Comparable to Duoderm CGF Thin®)

Our version of this popular dressing has a lower coefficient of friction than the leading competitor, helping to reduce friction and shear forces. It's thin profile is designed to minimize bunching and rolling up, leading to improved wear times and cost-effectiveness. Additionally, ours has exceptional stretch, allowing for them to be removed with minimal trauma to the wound bed and patient. Indicated for use for the management of lightly exuding Stage I Pressure Ulcers, post-operative wounds, other partial thickness wounds, and for protection of at risk skin sites (such as over the bony prominences). **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
85322	2" x 2"	20/box, 20 boxes/case, 400/case	A6234
85344	4" x 4"	10/box, 20 boxes/case, 200/case	A6234
85366	6" x 6"	5/box, 16 boxes/case, 80/case	A6235

Primacol® Bordered Hydrocolloid Dressing – Sterile

(Comparable to DuoDERM® Signal)

Primacol® hydrocolloid dressings are similar to our regular hydrocolloids, except a film-backing replaces the foam backing. The smooth film backing provides additional transparency for wound inspection, and helps reduce friction shear forces. Our dressing has very thin tapered edges, designed to minimize bunching and rolling up, leading to improved wear times and cost-effectiveness. Indicated for use for the management of lightly exuding Stage I-IV Pressure Ulcers, Diabetic foot ulcers, venous stasis ulcers, other chronic and acute wounds. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
85222	2" x 2"	20/box, 20 boxes/case, 400/case	A6234
85244	4" x 4"	10/box, 20 boxes/case, 200/case	A6234
85266	6" x 6"	5/box, 16 boxes/case, 80/case	A6235
85288	8" x 8"	5/box, 16 boxes/case, 80/case	A6236
85192	Sacral Shape	5/box, 16 boxes/case, 80/case	A6235

HYDROGEL DRESSINGS

Aquasite[®] Amorphous Hydrogel Dressing

(Comparable to Intrasite[®])

Amorphous Hydrogels are useful for rehydrating sloughy or necrotic tissue and enhancing autolytic debridement. On minimally draining wounds, they help to provide a moist wound environment, conducive to wound healing. Amorphous Hydrogel is used to fill uneven and difficult to reach spaces. The dressing is in a novel Bellow's tube, making precise application easier. Indicated for the management of acute and chronic wounds.



Product	Description	Pkg Unit/Case	HCPCS
81110	1 fl oz Bellows	12/case	A6248

Aquasite® Sheet Hydrogel – Sterile

(Comparable to Vigilon®)

As with amorphous hydrogels, sheet hydrogels are useful for the hydration of wound site. Our Aquasite® Hydrogel Sheet Dressings are used as a cover over deep wounds, or as a primary dressing over partial thickness wounds. Patients like hydrogel sheets due to their cooling and soothing properties upon application. Ours is designed to maintain its physical integrity, even after prolonged wound contact, for ease of removal. Indicated for the management of acute and chronic wounds. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
82023	2" x 3"	10/box, 4 boxes/case, 40/case	A6242
82044	4" x 4"	10/box, 4 boxes/case, 40/case	A6242

Aquasite® Impregnated Gauze – Sterile

(Comparable to TransiGel®)

Our impregnated dressings incorporate our Aquasite® hydrogel into 100% Cotton USP Type VII gauze pads of various sizes. The dressing are useful for; rehydrating sloughy or necrotic tissue, enhancing autolytic debridement, helping to provide a moist wound environment on minimally draining wounds, conducive to healing and a sterile, moist packing for full thickness wounds. These dressings are ideal for the management of acute and chronic wounds including Stage IV pressure ulcers. These dressings are not indicated for the management of third degree burns. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
53022	2" x 2"	10/box, 4 boxes/case, 40/case	A6231
53044	4" x 4"	10/box, 4 boxes/case, 40/case	A6231
53048	4" x 8"	10/box, 4 boxes/case, 40/case	A6232

Hydrocell® Thin Adhesive Foam Dressing – Sterile

(Comparable to Flexan®)

Our Hydrocell® Thin Adhesive Foam Dressing is a very thin and conformable dressing, with a skin-friendly adhesive layer on the wound contact side. These dressing are easy for clinicians to apply and to remove, with minimal trauma to the patient and to the wound bed. The polyurethane foam sheets are hydrophilic in nature, so the dressing can absorb a large amount of fluid, and provide a soft and cushioning layer for added patient comfort. It also has an outer film layer providing a barrier to external contaminants. For use on moderately exuding wounds. **Does not contain natural rubber latex.**




Product	Description	Pkg Unit/Case	HCPCS
84423	2" x 3"	10/box, 4 boxes/case, 40/case	A6209
84444	4" x 4"	10/box, 4 boxes/case, 40/case	A6209

Hydrocell® Adhesive Foam Dressing – Sterile

(Comparable to Tielle®)

Our Hydrocell® Adhesive Foam Dressing adds a foam center island to our thin adhesive foam for additional absorption. It also has an outer film layer providing a barrier to external contaminants. These dressings are easy to apply and to remove. For use on moderately to heavily exuding wounds. **Does not contain natural rubber latex.**

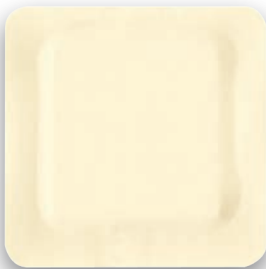

Product	Description	Pkg Unit/Case	HCPCS
84244	4" x 4"	10/box, 4 boxes/case, 40/case	A6212
84266	6" x 6"	10/box, 4 boxes/case, 40/case	A6213

Hydrocell® Foam Dressing – Sterile

(Comparable to Allevyn®)

Our Hydrocell® Foam Dressing is non-adhesive, and has an outer film layer providing a barrier to external contaminants. Like our other foams, they can absorb a large amount of fluid and are easy for clinicians to apply and remove with minimal trauma to the patient and to the wound bed. A secondary dressing or medical tape is needed to hold this dressing in place. For use on moderately to heavily exuding wounds. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
84544	4" x 4"	10/box, 4 boxes/case, 40/case	A6209
84566	6" x 6"	10/box, 4 boxes/case, 40/case	A6210

Sorbacell® Foam Dressing (Without Film Backing) – Sterile

(Comparable to Hydrasorb®)

Our Sorbacell® Foam Dressing provides maximum absorption of wound fluid. It is extremely flexible and conformable, and also can be used as a filler for deep cavities. Sorbacell is non-adhesive, therefore a secondary dressing or medical tape is needed to hold the dressing in place. This foam does not have an outer film layer and is ideal for use under compression. For use on moderately to heavily exuding wounds. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
84322	2" x 2"	10/box, 4 boxes/case, 40/case	A6209
84344	4" x 4"	10/box, 4 boxes/case, 40/case	A6209
84318	1" x 8" Foam Strips	10/box, 4 boxes/case, 40/case	A6209

WOUND CLEANSERS

Primaderm® Wound Cleanser – Non-Sterile

(Comparable to Saf-Cleans® AF)

This non-sterile, preserved dermal cleanser is intended for multi-use. It has a non-cytotoxic formula and is for use when extra cleansing power is needed. Available in a spray bottle that delivers solution at safe pressures. The nozzle can be set to either stream or spray to accommodate the desired pressure, within AHCPR guidelines, for removing debris or for cleansing the wound. The cleanser contains Poloxamer 188, an effective surfactant known to be gentle on tissue. Indicated for use for pressure ulcers, stasis ulcers, diabetic foot ulcers, post-surgical wounds, first and second degree burns, cuts, abrasions, and minor irritations of the skin.



Product	Description	Pkg Unit/Case	HCPCS
69201	500 ml/17.3 oz Spray Bottle	12/case	A6260
69202	118 ml/4.15 oz Squir Top Bottle	12/case	A6260

Primaderm® Dermal Wound Cleanser – Sterile

(Comparable to Shur-Cleans®)

This sterile, non-cytotoxic dermal cleanser is intended for single-use, making it the cleanser of choice for traumatic acute wounds. Ideally, the cleanser should be used when extra cleansing power is needed. It is packaged in a container with a squirt top that delivers the solution at safe pressures, within AHCPR guidelines. The cleanser contains Poloxamer 188, an effective surfactant known to be gentle on tissue. Can also be used for pressure ulcers, stasis ulcers, diabetic foot ulcers, post-surgical wounds, first and second degree burns, cuts, abrasions, and minor irritations of the skin.



Product	Description	Pkg Unit/Case	HCPCS
69101	118 ml/4.15 oz Squirt Top Bottle	12/case	A6260

Dermagran[®] Wound Cleanser With Zinc

Dermagran[®] Wound Cleanser is a gentle skin and wound cleanser that contains both Zinc and Vitamin B₆. As such, the wound cleanser is designed to not just cleanse, but to also prepare the wound bed with an optimal environment for wound healing. The multi-use spray bottle is designed so that the protective cap does not need to be fully removed in order to operate, increasing ease of use. This 360° delivery system meets AHCPH guidelines for irrigation pressure. The cleanser contains one of the lowest cytotoxic formulations on the market. When used along with Dermagran[®] B, these synergistic products can help to optimize your clinical outcomes.



Product	Description	Pkg Unit/Case	HCPCS
WC-04	4 oz bottle	12/case	A6260

TRADITIONAL WOUND CARE

- Sponges
- Wound Dressings
- Packing
- Island and Barrier Border Gauze Dressings
- Non-Adherent Impregnated Dressings
- Conforming Bandages
- Elastic Net Dressing (Retainer / Tube Gauze)
- Compression Bandages and Bandaging Systems
- Protective Products
- Tapes
- Laparotomy Sponges

SPONGES

Ducare[®] Woven Gauze Sponges

(Comparable to Curity[®] and J&J and Medline's Woven Gauze Sponges)

These woven sponges are made from 100% cotton. They are gold standard for general use or for debriding, covering and packing wounds. These dressings come in a variety of sizes and are available Sterile in packages of 2 and Non-Sterile in sleeves of 200. 4 x 4's also available in packages of 10, packed in trays. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
STERILE		
82208	2" x 2", 8 ply	2's, 50 pkgs/box, 30 boxes/case, 3000/case
82312	3" x 3", 12 ply	2's, 25 pkgs/box, 48 boxes/case, 2400/case
82408	4" x 4", 8 ply	2's, 25 pkgs/box, 24 boxes/case, 1200/case
82410	4" x 4", 12 ply	10's, 32 pkgs/tray, 4 trays/case, 1280/case
82412	4" x 4", 12 ply	2's, 25 pkgs/box, 24 boxes/case, 1200/case
82812	8" x 4", 12 ply	2's, 25 pkgs/box, 20 boxes/case, 1000/case
NON-STERILE		
90208	2" x 2", 8 ply	200/sleeve, 25 sleeves/case, 5000/case
90212	2" x 2", 12 ply	200/sleeve, 40 sleeves/case, 8000/case
90312	3" x 3", 12 ply	200/sleeve, 20 sleeves/case, 4000/case
90408	4" x 4", 8 ply	200/sleeve, 20 sleeves/case, 4000/case
90412	4" x 4", 12 ply	200/sleeve, 10 sleeves/case, 2000/case
90812	8" x 4", 12 ply (H Fold)	200/sleeve, 10 sleeves/case, 2000/case

Also offered in USP Type VII Gauze — See your Derma Sciences Representative for more information

Dupaque® X-Ray Detectable Woven Gauze Sponges

(Comparable to Vistec®, Ray-Tec®, and Medline's X-Ray Detectable Woven Gauze Sponges)



These woven sponges are made from 100% cotton and are low-linting, fluffy, and highly absorbent. They feature x-ray detectable barium sulfate strands for easy detection. Sterile sponges come in trays of 5 or 10, and the non-sterile sponges come in sleeves of 200. **Does not contain natural rubber latex.**

Product	Description	Pkg Unit/Case
STERILE		
9464S	4" x 4", 16 ply	10/tray, 128 trays/case, 1280/case
9467S	8" x 4", 12 ply	10/tray, 80 trays/case, 800/case
NON-STERILE		
94416	4" x 4", 16 ply	200/sleeve x 20 sleeves/case
94812	8" x 4", 12 ply	200/sleeve, 10 sleeves/case, 2000/case

Also offered in USP Type VII Gauze

Dusoft® Non-Woven Sponges

(Comparable to Versalon®, NuGauze® and Avant Gauze®)



These non-woven sponges are ideal for wound cleansing and also can be used as a secondary dressing for added absorbency. They produce less lint and are less adherent than woven gauze sponges, so cause less disruption to newly formed granulation tissue. The rayon/polyester blend provides fast wicking properties while being highly absorbent, minimizing skin maceration. **Does not contain natural rubber latex.**

Product	Description	Pkg Unit/Case
STERILE		
84122	2" x 2", 4 ply	2's, 25 pkgs/box, 60 boxes/case, 3000/case
84133	3" x 3", 4 ply	2's, 25 pkgs/box, 48 boxes/case, 2400/case
84144	4" x 4", 4 ply	2's, 25 pkgs/box, 24 boxes/case, 1200/case
84644	4" x 4", 6 ply	2's, 25 pkgs/box, 12 boxes/case, 600/case
NON-STERILE		
94122	2" x 2" 4 ply	200/sleeve, 40 sleeves/case, 8000/case
94133	3" x 3" 4 ply	200/sleeve, 20 sleeves/case, 4000/case
94144	4" x 4" 4 ply	200/sleeve, 10 sleeves/case, 2000/case

Dutop[®] Cover Sponges

(Comparable to Topper[®], Curity[®] Cover, and Medline's Post-Op Sponges)

These sponges are made with extra-absorbent cellulose fibers, and therefore they are more absorbent than traditional gauze sponges. The outer layer is made up of our non-woven fibers, making the dressing less adherent to newly formed tissue. Due to their high level of absorbency and non-adherent quality, these sponges are often used post-operatively. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
STERILE		
84943	4" x 3"	2's, 25 pkgs/box, 24 boxes/case, 1200/case
84944	4" x 4"	2's, 25 pkgs/box, 24 boxes/case, 1200/case
NON-STERILE		
94943	4" x 3"	100/sleeve, 20 sleeves/case, 2000/case (bulk)
94944	4" x 4"	100/sleeve, 20 sleeves/case, 2000/case (bulk)

Tracheotomy/Drain Sponges - Sterile

(Comparable to Exclon[®], Sof-Wick[®], and Medline's Woven and Avant Gauze Non-Woven Drain Sponges)

An added fenestration makes these sponges ideal for dressing around a tracheostomy tube or any type of tube on the body. Available in a soft 6 ply non-woven or in standard 12 ply gauze. The non-woven sponges provide superior absorption and are low-linting, fast-wicking, and non-adherent to wound beds. The woven sponges are ideal for when natural breathability is desired. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
84916	Non-Woven 4" x 4", 6-Ply	2's, 25 pkgs/box, 12 boxes/case, 600/case
84918	Non-Woven 2" x 2", 6-Ply	2's, 35 pkgs/box, 20 boxes/case, 1400/case
82910	Woven Gauze 4" x 4", 12-Ply	1's 25 pkgs/box, 24 boxes/case, 600/case

WOUND DRESSINGS

Gazetex[®] Bandage Rolls

(Comparable to Kerlix[®],
Kling-Fluff Rolls[®], and Bulkee II[®])



Gazetex bandage rolls are washed, low-linting woven gauze rolls. These bandages are used for wrapping or packing large and deep wounds. When used as a primary dressing, the dressing provides fast wicking, absorbency, and aeration. When used as a secondary dressing, the dressing provides additional cushioning and comfort. **Does not contain natural rubber latex.**

Product	Description	Pkg Unit/Case
STERILE		
9322	4½" x 147", 6 ply	1/pkg, 100 pkgs/case, 100/case
9321	2½" x 108", 6 ply	1/pkg, 100 pkgs/case, 100/case
NON-STERILE		
9320	4½" x 147", 6 ply	1/pkg, 100 pkgs/case, 100/case
9320-48	4½" x 147", 6 ply	12/sleeve, 4 sleeves/case, 48/case

Gazetex[®] 100% Cotton, 6-Ply, Fluff Sponge

(Comparable to Kerlix[®], Kling[®] and
Bulkee[®] Super Fluff Sponges)



These washed, fluffed gauze sponges are low-linting and soft. They are excellent for absorption, packing wounds, and managing exudate in and around the wound. The dressings are ideal for pre-op prepping and cleansing, as well as debriding. **Does not contain natural rubber latex.**

Product	Description	Pkg Unit/Case
STERILE		
9340	6" x 6¾", diagonally	2's, 240 pkgs/case, 480/case
9345	7¾" x 8¾", diagonally	2's, 80 pkgs/case, 160/case
NON-STERILE		
9333	6" x 6¾", diagonally	100/sleeve, 6 sleeves/case, 600/case
9334	7¾" x 8¾", diagonally	100/sleeve, 4 sleeves/case, 400/case

Dupad® Abdominal Pad

(Comparable to Curity® Abdominal Pads, SURGIPAD® Combine Dressings and Medline's ABD Pads)

Our abdominal pads are filled with cellulose fibers, making them super-absorbent. The fibers quickly absorb and disperse fluids laterally, to prevent pooling and maceration. The backing is constructed of a hydrophobic material, preventing strike-through. All edges are sealed to prevent linting. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
STERILE		
87059	5" x 9"	25/box, 8 boxes/case, 200/case
87088	8" x 7½"	25/box, 8 boxes/case, 200/case
87810	8" x 10"	25/box, 8 boxes/case, 200/case
NON-STERILE		
97159	5" x 9"	576/case (Bulk packaged)
97188	8" x 7½"	576/case (Bulk packaged)
97110	8" x 10"	432/case (Bulk packaged)

Primapad® Non-Adherent Absorbent Pad – Sterile

(Comparable to Telfa®, Release® Non-Adhering Dressings and Medline's Non-Adherent Pads)

Compares to the industry standard, this non adherent gauze pad is often chosen as the primary dressing when extra wound bed protection is needed. Air circulates through the soft, perforated film, while fluid is absorbed by the soft cotton pad. These dressings can be easily cut to fit any shape. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
83323	2" x 3"	100/box, 24 boxes/case, 2400/case
83334	3" x 4"	50/box, 18 boxes/case, 900/case
83338	3" x 8"	50/box, 12 boxes/case, 600/case

Surecare® Eye Pads – Sterile

(Comparable to Medline's Eye Pads)

Constructed of absorbent cotton filler covered with fine mesh gauze on both sides. Sterile and ready to use, this pad is ideal for applying medication or to cover and protect the eye. ***Does not contain natural rubber latex.***



Product	Description	Pkg Unit/Case
82911	2.625" x 1.625"	1's, 25/tray, 24 trays/case, 600/case

PACKING

Pak-Its® Gauze Packing Strips – Plain and Iodoform

(Comparable to Curity®, NuGauze® and Medline's Packing Strips)

These sterile, Plain and Iodoform Gauze Strips can be used dry to absorb drainage in exuding wounds or soaked with a medication of choice. Pre-cut sizes suit different wound types, saving time and money. A secondary dressing is needed to hold the strips in place and to protect the wound from outside contaminants. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
PLAIN PACKING STRIPS			
59121	¼" x 5 yd	12/case	A6407
59221	½" x 5 yd	12/case	A6407
59321	1" x 5 yd	12/case	A6407
59421	2" x 5 yd	12/case	A6407
IODOFORM PACKING STRIPS			
59146	¼" x 5 yd	12/case	A6266
59246	½" x 5 yd	12/case	A6266
59346	1" x 5 yd	12/case	A6266
59446	2" x 5 yd	12/case	A6266

Pak-Its® Gauze Packing Strips – Plain and Iodoform with Strip Delivery System

(Comparable to Curity®, NuGauze® and Medline's Packing Strips)

These packing strips are identical to our other strips except the bottles each have a special strip delivery system that ensures packing material is kept clean and tangle free, with the end always easy to locate at the top of the bottle. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
PLAIN PACKING STRIPS			
59120	¼" x 5 yd	12/case	A6407
59220	½" x 5 yd	12/case	A6407
59320	1" x 5 yd	12/case	A6407
59420	2" x 5 yd	12/case	A6407
IODOFORM PACKING STRIPS			
59145	¼" x 5 yd	12/case	A6266
59245	½" x 5 yd	12/case	A6266
59345	1" x 5 yd	12/case	A6266
59445	2" x 5 yd	12/case	A6266

ISLAND AND BARRIER BORDER GAUZE DRESSINGS

Compdress® Island Dressing – Sterile

(Comparable to Telfa® Island and Medline's Border Gauze Dressings)

These border gauze dressings can be used as either primary or secondary dressings. They are considered a lower cost-per-unit alternative for when absorptive dressings are required. The product is comprised of a non-adherent, absorptive central pad, with a soft conformable layer of gauze as a backing. There is a skin-friendly adhesive border around the perimeter.

Does not contain natural rubber latex.



Product	Description	Pkg Unit/Case	HCPCS
89022	2" x 2"	25/box, 4 boxes/case, 100/case	A6219
89044	4" x 4"	25/box, 4 boxes/case, 100/case	A6219
89066	6" x 6"	25/box, 4 boxes/case, 100/case	A6219
89046	4" x 6"	25/box, 4 boxes/case, 100/case	A6219
89048	4" x 8"	25/box, 4 boxes/case, 100/case	A6219
89041	4" x 10"	25/box, 4 boxes/case, 100/case	A6219
89068	6" x 8"	25/box, 4 boxes/case, 100/case	A6220
89042	4" x 14"	25/box, 4 boxes/case, 100/case	A6220

Dudress® Film Top Barrier Dressing – Sterile

(Comparable to AllDress® and Medline's Stratasorb® Composite Dressings)

These all-in-one type dressings are identical to our CompDress® Island Gauze dressings, with the addition of a polyurethane film barrier, which keeps out external contaminants, as well as an effective water barrier to permit showering and bathing. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
89144	4" x 4"	25/box, 4 boxes/case, 100/case	A6203
89166	6" x 6"	25/box, 4 boxes/case, 100/case	A6203
89168	6" x 8"	25/box, 4 boxes/case, 100/case	A6204

NON-ADHERENT IMPREGNATED DRESSINGS

Petrolatum Gauze Non-Adhering Dressing – Sterile

(Comparable to Curity® and Adaptic® PG and Medline's Petrolatum Gauze Dressings)

Petrolatum impregnated dressings are designed to provide non-adherent packing for full-thickness wounds, providing a moist environment that is conducive to wound healing. They are made of fine, soft, conformable gauze impregnated with white petrolatum. Being non-adherent, removal of these dressings causes minimal trauma to the wound bed and patient. Indicated for use of management of full-thickness chronic wounds such as stage II-IV pressure ulcers, tunneling wounds and non-infected wounds. The "overwrap" version of these dressings provides an additional layer of sterility for use in environments such as operating rooms. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
DKC20049	½" x 72"	12/box, 6 boxes/case, 72/case	A6266
DKC20053	1" x 36"	12/box, 6 boxes/case, 72/case	A6266
DKC20056	3" x 9"	12/box, 6 boxes/case, 72/case	A6223
DKC20059	3" x 18"	12/box, 6 boxes/case, 72/case	A6266
DKC20062	3" x 36"	12/box, 6 boxes/case, 72/case	A6266
DKC20065	6" x 36"	12/box, 6 boxes/case, 72/case	A6266
DKC77030	1" x 8"	50/box, 4 boxes/case, 200/case	A6266

Xeroform Petrolatum Dressing – Sterile

(Comparable to Curity[®] and Adaptic[®] X and Medline's Xeroform Petrolatum Gauze Dressings)

Xeroform impregnated dressings are similar to our Petrolatum Gauze dressings, with the addition of 3% bismuth tribromophenate. They are the impregnated dressing of choice for when mild medication or deodorizing are required, such as post-operative applications, 1st and 2nd degree burns, and skin grafts. The "overwrap" version of these dressings provides an additional layer of sterility for use in environments such as operating rooms. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
DKC77034	1" x 8"	50/box, 4 boxes/case, 200/case	A6266
DKC20068	5" x 9"	50/box, 4 boxes/case, 200/case	A6223
DKC71164	2" x 2"	25/box, 6 boxes/case, 150/case	A6222
DKC71165	4" x 4"	25/box, 6 boxes/case, 150/case	A6222

Shur-Conform[®] Oil Emulsion Non-Adhering Dressing – Sterile

(Comparable to Curity[®] Non-Adhering Dressings and Adaptic[®] and Medline's Oil Emulsion Dressings)

Oil Emulsion impregnated dressings are designed to provide non-adherent packing for wounds, providing a moist environment conducive to wound healing. They are made of a knitted cellulose acetate fabric coated with a formulated petrolatum emulsion. The knitted fabric allows wound fluid to move through the dressing and into a secondary absorbent dressing. Being non-adherent, removal of these dressings causes minimal trauma to the wound bed and patient. Used for the management of skin grafts, surgical sites, abrasions, lacerations, dermal ulcers, and cosmetic surgery sites. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
DKC77041	3" x 3"	50/box, 12 boxes/case, 600/case	A6222
DKC77042	3" x 8"	24/box, 6 boxes/case, 144/case	A6223
DKC77043	3" x 16"	36/box, 6 boxes/case, 216/case	A6223
DKC77044	3" x 8", 3's	36/box, 6 boxes/case, 216/case	A6223

CONFORMING BANDAGES

Dutex[®] 100% Cotton 2 Ply Conforming Bandage

(Comparable to Kerlix Lite[®], Kling[®]
and Medline's Bulkee[®] II Lite)

This 100% cotton 2-ply woven bandage is ideal as a secondary dressing for wrapping legs and arms and to hold dressings in place. The dressing provides light stretch which enables it to conform to difficult body contours, while also providing light compression. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case	HCPCS
STERILE			
77781	2" x 4.1 yards	12/sleeve, 8 sleeves/case, 96/case	A6445
77782	3" x 4.1 yards	12/sleeve, 8 sleeves/case, 96/case	A6446
77783	4" x 4.1 yards	12/sleeve, 8 sleeves/case, 96/case	A6446
77784	6" x 4 ½" yards	6/sleeve, 8 sleeves/case, 48/case	A6447
NON-STERILE			
76781	2" x 4.1 yards	12/sleeve, 8 sleeves/case, 96/case	A6442
76782	3" x 4.1 yards	12/sleeve, 8 sleeves/case, 96/case	A6443
76783	4" x 4.1 yards	12/sleeve, 8 sleeves/case, 96/case	A6443
76784	6" x 4 ½" yards	6/sleeve, 8 sleeves/case, 48/case	A6444

Duform[®] Synthetic Conforming Bandage, Knitted

(Comparable to Conform[®],
Sof-Kling[®], Medline's Sof-Form)



This knitted stretch gauze roll is ideal as a secondary dressing for wrapping legs and arms and to hold dressings in place. The dressings gently secure to body contours without restricting movement. *Does not contain natural rubber latex.*

Product	Description	Pkg Unit/Case	HCPCS
STERILE			
75202	2" x 4.1 yds stretched	12/sleeve, 8 sleeves/case, 96/case	A6445
75203	3" x 4.1 yds stretched	12/sleeve, 8 sleeves/case, 96/case	A6446
75204	4" x 4.1 yds stretched	12/sleeve, 8 sleeves/case, 96/case	A6446
75206	6" x 4½" yds stretched	6/sleeve, 8 sleeves/case, 48/case	A6447
NON-STERILE			
75102	2" x 4.1 yds stretched	12/sleeve, 8 sleeves/case, 96/case	A6442
75103	3" x 4.1 yds, stretched	12/sleeve, 8 sleeves/case, 96/case	A6443
75104	4" x 4.1 yds, stretched	12/sleeve, 8 sleeves/case, 96/case	A6443
75106	6" x 4½" yds, stretched	6/sleeve, 8 sleeves/case, 48/case	A6444

ELASTIC NET DRESSING (RETAINER / TUBE GAUZE)

Surgilast® Tubular Elastic Dressing Retainer

(Comparable to Medline's Elastic Net Retention Dressings)

A tubular elastic stretch net designed to serve as a secondary dressing, applying gentle pressure to keep bandages or other medical devices securely in place without adhesive tape. Ideal for securing burn or post-op dressings, hot or cold packs, IV tubes, incontinence pads and decubitus ulcer dressings. Easy to customize. Allows maximum air-flow. Available in a wide range of sizes.



Product	Size	Working Stretch	Applications	Pkg Unit/ Case
GL701	1	6 7/8"	Fingers, toes, wrist	25 yards, 1/box
GL501	1	6 7/8"	Fingers, toes, wrist	50 yards, 1/box
GL702	2	8"	Small: hand, arm, leg, foot	25 yards, 1/box
GL502	2	8"	Small: hand, arm, leg, foot	50 yards, 1/box
GL720	3	10 1/8"	Medium: hand, arm, leg, foot	10 yards, 1/box
GL703	3	10 1/8"	Medium: hand, arm, leg, foot	25 yards, 1/box
GL503	3	10 1/8"	Medium: hand, arm, leg, foot	50 yards, 1/box
GL704	4	12 1/4"	Large: hand, arm, leg, foot	25 yards, 1/box
GL504	4	12 1/4"	Large: hand, arm, leg, foot	50 yards, 1/box
GL722	5	15"	Small: head, shoulder, thigh	10 yards, 1/box
GL705	5	15"	Small: head, shoulder, thigh	25 yards, 1/box
GL505	5	15"	Small: head, shoulder, thigh	50 yards, 1/box
GL706	5 1/2"	19 3/4"	Medium: head, shoulder, thigh	25 yards, 1/box
GL506	5 1/2"	19 3/4"	Medium: head, shoulder, thigh	50 yards, 1/box
GL724	6	25 1/2"	Large: head, shoulder, thigh	10 yards, 1/box
GL707	6	25 1/2"	Large: head, shoulder, thigh	25 yards, 1/box
GL507	6	25 1/2"	Large: head, shoulder, thigh	50 yards, 1/box
GL708	7	29"	Small: chest, back, perineum, axilla	25 yards, 1/box
GL508	7	29"	Small: chest, back, perineum, axilla	50 yards, 1/box
GL709	8	32 1/2"	Medium: chest, back, perineum, axilla	25 yards, 1/box
GL509	8	32 1/2"	Medium: chest, back, perineum, axilla	50 yards, 1/box
GL710	9	36"	Large: chest, back, perineum, axilla	25 yards, 1/box
GL510	9	36"	Large: chest, back, perineum, axilla	50 yards, 1/box
GL726	10	38"	X-Large: chest, back, perineum, axilla	10 yards, 1/box
GL711	10	38"	X-Large: chest, back, perineum, axilla	25 yards, 1/box
GL511	10	38"	X-Large: chest, back, perineum, axilla	50 yards, 1/box
GL712	11	40"	Special Sizing	25 yards, 1/box
GL512	11	40"	Special Sizing	50 yards, 1/box
GL713	12	42"	Special Sizing	25 yards, 1/box
GL714	13	44"	Special Sizing	25 yards, 1/box
GL715	14	46"	Special Sizing	25 yards, 1/box
GL2522	22	72"	Special Sizing	25 yards, 1/box

Surgilast® Tubular Elastic Dressing Retainer

(Comparable to Medline's Elastic Net Retention Dressings)

Provides the same features and benefits as Surgilast® while eliminating the latex sensitivity.

Does not contain natural rubber latex.



Product	Size	Working Stretch	Applications	Pkg Unit/ Case
GLLF2501	1	5 3/8"	Fingers, toes, wrist	25 yards, 1/box
GLLF2502	2	7"	Small: hand, arm, leg, foot	25 yards, 1/box
GLLF2503	3	9"	Medium: hand, arm, leg, foot	25 yards, 1/box
GLLF2504	4	11"	Large: hand, arm, leg, foot	25 yards, 1/box
GLLF2505	5	13 3/4"	Small: head, shoulder, thigh	25 yards, 1/box
GLLF2505.5	5 1/2"	18"	Medium: head, shoulder, thigh	25 yards, 1/box
GLLF2506	6	23"	Large: head, shoulder, thigh	25 yards, 1/box
GLLF2507	7	28"	Small: chest, back, perineum, axilla	25 yards, 1/box
GLLF2508	8	31 1/2"	Medium: chest, back, perineum, axilla	25 yards, 1/box
GLLF2509	9	33"	Large: chest, back, perineum, axilla	25 yards, 1/box
GLLF2510	10	37 1/2"	X-Large: chest, back, perineum, axilla	25 yards, 1/box
GLLF2511	11	40"	Special Sizing	25 yards, 1/box
GLLF2512	12	41 1/2"	Special Sizing	25 yards, 1/box
GLLF2514	14	45"	Special Sizing	25 yards, 1/box

Surgilast® Pre-Cut Tubular Elastic Dressing Retainer

A pre-cut tubular elastic dressing retainer that is pre-shaped to accommodate the more difficult-to-fit parts of the body. Ideal for securing hot or cold packs, ear dressings, electrodes, etc.



Product	Width	Application	Pkg Unit/ Case
GL750	Stress Vest	Small/Medium	20/box
GL752	Stress Vest	Large/X-Large	20/box
GL600	Head	Universal Fit	20/box
GL620	Perineum Panty	Small/Medium	20/box
GL622	Perineum Panty	Large/X-Large	20/box
GL641	Knee, Foot, Elbow, Hand	Universal Fit	20/box

Surgitube® Tubular Gauze

50 yard rolls

(For use WITH applicator)

A tubular gauze bandage designed to protect and secure primary dressings. Made from 100% white or beige cotton yarn knitted in a continuous seamless tube. Easy to apply. Available in a range of sizes for use on every area of the body. Manufactured in two styles for use with or without a special applicator. **Does not contain natural rubber latex.**



Product	Size	Working Stretch	Width	Pkg Unit/ Case
GL219	1	Small fingers, toes (WHITE)	5/8"	1/box
GL209	1	Small fingers, toes (BEIGE)	5/8"	1/box
GL220	2	Large fingers, toes (WHITE)	1"	1/box
GL210	2	Large fingers, toes (BEIGE)	1"	1/box
GL221	3	Hands, wrists, feet (WHITE)	1 1/2"	1/box
GL211	3	Hands, wrists, feet (BEIGE)	1 1/2"	1/box
GL222	4	Arms, lower leg (WHITE)	2 5/8"	1/box
GL212	4	Arms, lower leg (BEIGE)	2 5/8"	1/box
GL223	5	Adult legs, thighs, head (WHITE)	3 5/8"	1/box
GL224	T-1	Torso (WHITE)	5"	1/box
GL225	T-2	Torso (WHITE)	7"	1/box

Traditional Wound Care

Surgitube® Tubular Gauze

(For use Without applicator)

Tube gauze that secures dressings without tape. Allows maximum airflow and will not unravel when cut. Made of 100% Cotton. Also available in other sizes upon request. **Does not contain natural rubber latex.**



Product	Size	Description	Width	Pkg Case	Unit/
50 YARD ROLLS					
GL241	1P	Small fingers, toes (WHITE)	5/8"		1/box
GL246	1P	Small fingers, toes (BEIGE)	5/8"		1/box
GL242	2P	Large fingers, toes (WHITE)	7/8"		1/box
GL247	2P	Large fingers, toes (BEIGE)	7/8"		1/box
GL242A	2AP	Infant legs, bulky fingers (WHITE)	1 1/8"		1/box
GL243	3P	Hands, feet, lower legs, shoulders (WHITE)	1 1/2"		1/box
GL248	3P	Hands, feet, lower legs, shoulders (BEIGE)	1 1/2"		1/box
GL244	4P	Wrists, elbows (tight weave) (WHITE)	1 1/2"		1/box
GL245	5P	Head, shoulders, breasts, legs (WHITE)	2"		1/box
GL260	5P	Head, shoulders, breasts, legs (BEIGE)	2"		1/box
15 YARD ROLLS					
GL251	1	Small fingers, toes (WHITE)	5/8"		1/box
GL252	2	Large fingers, toes (WHITE)	7/8"		1/box
GL252A	2A	Infant legs, bulky fingers (WHITE)	1 1/8"		1/box
10 YARD ROLLS					
GL110W	1	Small fingers, toes (WHITE)	5/8"		1/box
GL210W	2	Large fingers, toes (WHITE)	7/8"		1/box
GL210F	2	Large fingers, toes (BEIGE)	7/8"		1/box
GL310W	3	Hands, feet, lower legs, shoulders (WHITE)	1 1/2"		1/box
5 YARD ROLLS					
GL105W	1	Small fingers, toes (WHITE)	5/8"		1/box
GL205W	2	Large fingers, toes (WHITE)	7/8"		1/box
GL205F	2	Large fingers, toes (BEIGE)	7/8"		1/box
GL305W	3	Hands, feet, lower legs, shoulders (WHITE)	1 1/2"		1/box
GL305F	3	Hands, feet, lower legs, shoulders (WHITE)	1 1/2"		1/box
GL405W	4	Wrists, elbows (tight weave) (WHITE)	1 1/2"		1/box

SurgiGrip® Tubular Elastic Support Bandage

(Comparable to Tubigrip®)

A Does not contain natural rubber latex tubular elastic bandage designed to provide support and compression. Made of a soft cotton knit interlaced with elastic threads. Perfect for sprains and strains, support following cast removal, soft-tissue injuries, joints and muscles, post-graft pressure for burn patients, and edema management/ control. Convenient alternative to flat wrap bandages. Can be applied to sensitive areas using an applicator and stays in one place with no clips or tape. Machine washable. 11 yard (10 meter) length. 8-12 mm Hg compression at 50% stretch. **Does not contain natural rubber latex.**



Product	Size	Application	Pkg Unit/ Case
GLB10	2 ½" (6.25 cm)	Small hands and limbs	1/box
GLC10	2 ¾" (6.75 cm)	Adult hands, arms and legs	1/box
GLD10	3" (7.5 cm)	Large arms and legs	1/box
GLE10	3 ½" (8.75 cm)	Legs or small thighs	1/box
GLF10	4" (10.00 cm)	Large knees or thighs	1/box
GLG10	4 ½" (12.00 cm)	Large thighs	1/box
GLJ10	6 ¾" 17.5 cm)	Small trunk	1/box
GLK10	8" (20.32 cm)	Medium trunk	1/box
GLL10	12¾" (32.38 cm)	Large trunk	1/box

SurgiGrip® Tubular Elastic Support Bandage

(Comparable to Tubigrip®)

A tubular elastic bandage containing latex, designed to provide support and compression. Perfect for sprains and strains, support following cast removal, soft-tissue injuries, joints and muscles, post-graft pressure for burn patients, and edema management/control. Convenient alternative to flat wrap bandages. Can be applied to sensitive areas using an applicator and stays in one place with no clips or tape. Machine washable. 11 yard (10 meter) length.



Product	Size	Application	Pkg Unit/Case
WMB10	2 ½" (6.25 cm)	Small hands and limbs	1/box
WMC10	2 ¾" (6.75 cm)	Adult hands, arms and legs	1/box
WME10	3 ½" (8.75 cm)	Legs or small thighs	1/box
WMJ10	6 ¾" (17.5 cm)	Small trunk	1/box

Surgitube® and SurgiGrip® Tubular Gauze Applicators

Used with certain styles of Surgitube® Tubular Gauze, so primary bandage is not dislodged. Also used with certain styles of SurgiGrip® Tubular Elastic Support Bandages to eliminate possible irritation to sensitive areas.



Product	Size	Used With	Pkg Unit/Case
METAL CAGE APPLICATORS			
GL230	0M	GL209, GL219	1/box
GL231	1M	GL209, GL219	1/box
GL232	2M	GL210, GL220, GLB10, WMB10 GL211, GL221, GLC10, GLD10, GLE10,	1/box
GL234	3M	WMC10, WMD10, WME10, GL212, GL222, GLF10, GLG10, WMF10	1/box
GL236	4M	WMG10	1/box
LUCITE TUBE APPLICATORS			
GL231P	1L	GL209, GL219	1/box
GL232P	2L	GL210, GL220, GL211, GLC10, GLD10, GLE10	1/box
GL234P	3L	WMC10, WMD10, WME10	1/box
ALUMINUM SPLINT TYPE APPLICATORS			
GLALUAPP		GL209, GL219, GL110W, GL210W, GL210F, GL105W, GL105F, GL205W, GL205F	12/case
APPLICATOR SETS			
GL227	1, 2, 3	Metal Cage Applicators	1 set/box

COMPRESSION BANDAGES AND BANDAGING SYSTEMS

Dusor® Elastic Bandage – Economy Quality

(Comparable to Tensor® and Medline's Sure Wrap®)



Quality, cotton-wrapped heat resistant rubber strands provides firm compression and comfortable support. The woven fabric is ideal for pressure dressings, vascular and muscle support. Enclosed clips provide fast convenient application. The bandage is durable for repeated laundering.

Product	Description	Pkg Unit/Case
71412	2" x 5 yds	12/case
71413	3" x 5 yds	12/case
71414	4" x 5 yds	12/case
71416	6" x 5 yds	12/case

PrimaCare® Elastic Bandage with Velcro Closure

(Comparable to Dyna Flex® and Medline's Matrix Elastic Bandages)



This is a high-quality rubber elastic bandage ideal for high compression treatment on various conditions. The Velcro self-closure eliminates the need for clips or tape. **Does not contain natural rubber latex.**

Product	Description	Pkg Unit/Case
71422	2" x 5 yds	12/case
71423	3" x 5 yds	12/case
71424	4" x 5 yds	12/case
71426	6" x 5 yds	12/case

Flexoplast® Elastic Adhesive Bandage

(Comparable to Elasticon®, Elastoplast® and Co-Flex™)



An all elastic adhesive bandage that provides compression. Adapts easily to contours of the body. Used in combination with Primer® Unna Boot to effectively manage venous leg ulcers. Also provides support for fractures and muscle injuries. Sticks to itself and will not slip, bunch, or loosen.

Product	Description	Pkg Unit/Case
GL1001	1" x 5 yds, stretched	12/case
GL1002	2" x 5 yds, stretched	6/case
GL1003	3" x 5 yds, stretched	4/case
GL1004	4" x 5 yds, stretched	6/case
GL1003R	3" x 5 yds, stretched	12/case

CompriBand® Short-Stretch Compression Bandage

(Comparable to Comprilan®)

This bandage is specifically designed to manage the effects of venous insufficiency. This 100% Cotton elastic bandage can be worn by both ambulatory and non-ambulatory patients and is machine washable. It may remain in place for up to 7 days and can be removed for sleeping or bathing. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
71201	1½" x 5.4 yds	12/case
71202	2.3" x 5.4 yds	12/case

Duban[®] Cohesive Bandage

(Comparable to Coban[®], Dyna Flex[®] Cohesive Compression Bandages and Medline's Soft Wrap[®] Elastic Bandage)

This cohesive elastic bandage is a high quality elastic wrap that offers protection and support. It adheres to itself, but never to the skin or hair. It conforms to body contours while offering controlled compression. It is lightweight and water resistant. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
73201	1" x 5 yds	12/case
73202	2" x 5 yds	12/case
73203	3" x 5 yds	12/case
73204	4" x 5 yds	12/case

Sorban[®] Undercast Padding Bandage

(Comparable to Webri[®] and Medline's Wytex Undercast Padding)

This bandage is a lightweight 100% Cotton bandage that offers protection of bony prominences and absorption of wound exudate. It is primarily used under other forms of compression bandaging and hard casts. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
79504	4" x 5 yds	12/case

Dufore® Four-Layer Bandaging System

(Comparable to Profore®)

This system is designed to sustain effective levels of graduated compression for up to seven days. Sustained compression has been recognized as the most important factor in managing venous leg ulcers. The system is comprised of four layers. The first layer, an absorptive padding layer, provides padding for bony prominences, distributes compression evenly over the leg and absorbs exudate. The second layer provides short-stretch compression by resisting expansion of the calf muscle during limb movement. It works with the third and fourth layers elastic compression bandages to effectively minimize edema and increase venous return. A wound contact dressing is included with the system to help maintain a moist wound healing environment and to protect the wound from lint or outside contaminants. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
72414	4-Layer System	8/case

Primer[®] Modified Unna Boot

For the ambulatory management of venous stasis leg ulcers. A gauze bandage that has been evenly impregnated with a non-hardening paste of zinc oxide. Available with or without calamine. Ideal for managing: venous stasis ulcers, lymphatic edema, minor sprains, strains, fractures, thrombophlebitis and skin grafts on



burn wounds. Contains no gelatin which prevents the hardening or caking of the bandage. Made without preservatives and with non-raveling thread locked gauze. Packaged in a unique vacuum-sealed pouch.

Product	Description	Pkg Unit/Case
GL3001	3" x 10 yds	12/case
GL4001	4" x 10 yds	12/case
GL3001C	3" x 10 yds (with Calamine)	12/case
GL4001C	4" x 10 yds (with Calamine)	12/case

Unna-Pak[®]: Primer[®] Modified Unna Boot and Duban[®] Self Adherent Bandage

(Comparable to Unna-Flex[®] Plus)

For the ambulatory management of venous stasis leg ulcers. A convenient system containing two high performance bandages that together form a two-layer Unna Boot. The Unna-Pak Primer has all the advantage of the Primer[®] Modified Unna Boot with the additional benefit of a flesh colored secondary dressing that is lightweight, strong and easy to apply. The Unna Boot is 10 yards and the self-adherent bandage is 5 yards.



Product	Description	Pkg Unit/Case
GL2003	3"	12/case
GL2004	4"	12/case

Primer® Flex Elastic Unna Boot

(Comparable to Unna-Flex)

For the ambulatory management of venous stasis leg ulcers. A stretchable and non-rigid gauze bandage that has been evenly impregnated with a non-hardening paste of zinc oxide. Available with or without calamine. Ideal for managing: venous stasis ulcers, lymphatic edema, minor sprains, strains, fractures, thrombophlebitis and skin grafts on burn wounds. Contains no gelatin which prevents the hardening or caking of the bandage. Made without preservatives and remains flexible. Packaged in a unique vacuum-sealed pouch.



Product	Description	Pkg Unit/Case
GL3001F	3" x 10 yds	12/case
GL4001F	4" x 10 yds	12/case

PROTECTIVE PRODUCTS

Glen-Sleeve® II Arm and Leg Protectors

(Comparable to Skinsleeve®)

Knitted tubular sleeves designed to provide mild compression and protect the fragile skin of long-term care residents against everyday scrapes, shearing and bumps. Cotton and lycra construction with a polypropylene fiber that creates a wicking action to remove excess moisture from the skin. Durable sewn lockstitch construction with finished hems on both ends. Unique shaped heel (in leg protector). **Does not contain natural rubber latex.**



Product	Description	Size	Pkg Unit/Case
GL1000	Hand, Wrist, Arm (WHITE)	18" (L) x 3"(W)	12 pairs/case
GL1000B	Hand, Wrist, Arm (BEIGE)	18" (L) x 3"(W)	12 pairs/case
GL1000WP	Padded Hand, Wrist, Arm (WHITE)	13 ½" (L) x 3"(W)	12 pairs/case
GL1100B	Petite Hand, Wrist, Arm (BEIGE)	16 ½" (L) x 2"(W)	12 pairs/case
GL2000	Hand, Wrist, Thumb, Arm (WHITE)	18" (L) x 3"(W)	12 pairs/case
GL3000	Below Knee (WHITE)	15" (L) x 3 ½" (W)	12 pairs/case
GL3000B	Below Knee (BEIGE)	15" (L) x 3 ½" (W)	12 pairs/case

Glen-Sleeve® Arm Protectors

(Comparable to Skinsleeve®)

Knitted tubular sleeves designed to provide mild compression and protect the fragile skin of long-term care residents against everyday scrapes, shearing and bumps. The original Glen-Sleeve® designed with a color-coded strip for easy identification. Made from polyester, cotton, and nylon spandex.



Product	Description	Size	Pkg Unit/Case
GL899	Small/Medium (yellow band)	18" – 20" (L) x 3 ¼" (W)	12 pairs/case
GL900	Large/X-Large (blue band)	18" – 20" (L) x 3 ½" (W)	12 pairs/case

Heel & Elbow Protectors

(Comparable to Heelbo®)

Designed to protect against decubitus ulcers and other skin damage common among bed or wheelchair-bound patients. Consists of a knitted tube of soft fabric with a tricot-covered urethane foam pad stitched inside. Outer fabric moves independently of foam pad, foam pad stays in contact with the skin. Easy to apply, machine washable and dryable.



Product	Description	Size	Pkg Unit/Case
GL800	Small (yellow band)	9 ½"	12 pairs/case
GL801	Medium/Large (blue band)	11"	12 pairs/case
GL802	X-Large (green band)	11 ¾" – 12 ¼"	12 pairs/case

TAPES

Ultrafix[®] Self-Adhesive Dressing Retention Tape

(Comparable to Mefix[®], Hypafix[®], and Medline's Medfix[®])

A non-sterile and non-woven fixative tape that is breathable and especially gentle to the skin. Available in multiple widths that can be cut to size. Used to secure primary dressings such as gauze and foams. *Does not contain natural rubber latex.*



Product	Description	Pkg Unit/Case
62036	2" x 11 yards	10/case
62037	4" x 11 yards	10/case
62038	6" x 11 yards	10/case

LAPAROTOMY SPONGES

Laparotomy Sponges, Prewashed, X-Ray Detectable

(Comparable to Curity[®] and Medical Action Laparotomy Sponges)



Our Laparotomy Sponges are made of 100% cotton gauze. The sponges are soft, have high absorbency, and are packaged pre-washed. All edges are folded in and stitched with cotton thread, preventing the possibility of any loose fibers. **Does not contain natural rubber latex.**

Product	Description	Pkg Unit/Case
STERILE		
1205	12" x 12", 4 ply, (molded pouch)	10/pkg, 40 pkgs/case, 400/case
1206	12" x 12", 4 ply, (flat pouch)	5/pkg, 60 pkgs/case, 300/case
1213	18" x 18", 4 ply, (flat pouch)	5/pkg, 40 pkgs/case, 200/case
1214	18" x 18", 4 ply, (molded pouch)	5/pkg, 40 pkgs/case, 200/case
1223	4" x 18", 4 ply, (molded pouch)	5/pkg, 60 pkgs/case, 300/case
1236	8" x 36", 4 ply, (molded pouch)	2/pkg, 50 pkgs/case, 100/case
12731	8" x 24", 8 ply, (molded pouch)	2/pkg, 50 pkgs/case, 100/case
12751	8" x 24", 8 ply, (molded pouch)	1/pkg, 100 pkgs/case, 100/case
12931	8" x 36", 4 ply, (molded pouch)	2/pkg, 50 pkgs/case, 100/case
NON-STERILE		
1201	12" x 12", 4 ply, folded	5/pkg, 100 pkgs/case, 500/case
1211	18" x 18", 4 ply folded	5/pkg, 50 pkgs/case, 250/case
1221	4" x 18", 4 ply folded	5/pkg, 128 pkgs/case, 640/case
1231	8" x 36", 4 ply	1/pkg, 100 pkgs/case, 100/case
1271	18" x 24", 8 ply	1/pkg, 100 pkgs/case, 100/case

A decorative graphic at the top of the page featuring overlapping, semi-transparent geometric shapes in shades of yellow and orange, creating a layered, abstract effect.

BURN CARE

- Silver Antimicrobial Burn Care Dressings
- Gauze Burn Dressings
- Other Burn Care Products

SILVER ANTIMICROBIAL BURN CARE DRESSINGS

Silverseal® Burn Contact Dressing

Silverseal® Burn Contact Dressings are designed for placement directly on the burn. The Silverseal® Burn Contact Dressing provides a topical antimicrobial barrier that reduces the bacteria and fungi counts on the wound surface and within 2 to 4 millimeters of the burn surface. The dressing can be left in place for up to seven days depending on the amount of drainage and necrotic tissue present. Silverseal® Burn Contact Dressings are composed of knitted fabrics that are flexible, minimally adhering, and antimicrobial. The knitted fabric is 100% covered, uniformly and circumferentially with a thin coat of pure metallic silver approximately (4/100,000 inch). The fabric allows for these dressings to be extremely versatile and user friendly, conforming to hard-to-dress areas and providing mild compression when needed. Indicated for use as an antimicrobial barrier over 1st and 2nd degree burns, full and partial thickness wounds as well as graft and donor sites.



Product	Description	Pkg Unit/Case
SS-BWCD-0404	4 ¼" x 4 ¼"	10/box
SS-BWCD-0408	4" x 8"	10/box
SS-BWCD-0816	8" x 16"	5/box
SS-BWCD-1616	16" x 16"	5/box

Silverseal® Burn Wrap

Silverseal® Burn Wrap is identical to the Silverseal® Burn Contact Dressing with the exception that the product is cut in rolls. These large sizes allow for multiple uses over hard-to-dress areas and large body surface areas, while also providing slight compression when needed.

May be used as an antimicrobial barrier over 1st and 2nd degree burns, full and partial thickness wounds as well as graft and donor sites.



Product	Description	Pkg Unit/Case
SS-EBW-0466	4" x 66"	1/pouch
SS-EBW-06108	6" x 108"	1/pouch

Silverseal® Tubular Stretch Knit

Silverseal® Tubular Stretch Knit Products are designed to bring antimicrobial protection to digits and certain hard to dress extremities. These easy to use products are gently stretched to apply, and can be used as primary dressings. The Silverseal® dressings should then be covered or secured with conventional technique. Indicated for use as an antimicrobial barrier over 1st and 2nd degree burns, full and partial thickness wounds as well as graft and donor sites.



Product	Description	Pkg Unit/Case
SS-BKWCD-5015	½" x 1 ½"	5/box
SS-BKWCD-0510	½" x 10"	5/box
SS-BKWCD-1035	1" x 3 ½"	5/box
SS-BKWCD-115	1" x 15"	5/box
SS-BKWCD-34	3" x 4"	5/box

Silverseal® Burn Gloves

Silverseal® Burn Gloves are made from the same material used for the Burn Contact Dressings. The gloves are open finger and oversized 20% for use immediately after acute injury or post surgery for both partial to full thickness burns. Use of these gloves can save time and increase patient compliance. They offer broad antimicrobial protection to often hard to dress parts of the hand. May be used as an antimicrobial barrier over 1st and 2nd degree burns, full and partial thickness wounds as well as graft and donor sites.



Product	Description	Pkg Unit/Case
SS-BWCD-M1	Silver Acute Burn Gloves – Medium	1/pouch
SS-BWCD-L1	Silver Acute Burn Gloves – Large	1/pouch

GAUZE BURN DRESSINGS

Dupress® Burn Dressings

These lint free pads help minimize trauma and adherence to the wound. A wide-array of sizes and designs are available to meet your individual needs. Our burn pads are made of a fine mesh gauze cover, cellulose fill to absorb large amounts of fluid and stitched edges.



Product	Description	Pkg Unit/Case
STERILE		
11903	4" x 4", 32 ply, Burn Gauze, 20 x 12 mesh wrapped	10/pkg, 50 pkgs/case 500/case
11904	16" x 16", 16 ply, Burn Gauze 20 x 12 mesh wrapped	2/pkg, 40 pkgs/case, 80/case
11907	18" x 22" compress, rayon / cellulose wrapped in Gauze	1/pkg, 30 pkgs/case, 30/case
11908	24" x 36" compress, rayon / cellulose wrapped in Gauze	1/pkg, 24 pkgs/case, 24/case
8002	4" x 12", 12 ply, Burn Gauze, 20 x 12 mesh	20/pkg, 40 pkgs/case, 800/case
8007	9" x 5yds, Durlix, fluff gauze	1/pkg, 50 pkgs/case, 50/case
8118S	18" x 36" folded gauze, 20 x 12 mesh 200/case	2/pkg, 100 pkgs/case,
8136S	18" x 72" folded gauze, 20 x 12 mesh	2/pkg, 50 pkgs/case, 100/case
9825	9" x 9", 8 ply, gauze pad (10's)	10/pkg, 40 pkgs/case, 400/case

OTHER BURN CARE PRODUCTS

IN ADVANCED WOUND CARE

MEDIHONEY® Wound & Burn Dressings	5-7
ALGICELL® Ag Silver Antimicrobial Dressings	14
Aquasite® Amorphous Hydrogel Dressing	21
Aquasite® Sheet Hydrogel Dressing	22
Primaderm® Dermal Cleanser	25
Primaderm® Dermal Wound Cleanser	25

IN TRADITIONAL WOUND CARE

Xeroform Petrolatum Dressing	38
Shur-Conform® Oil Emulsion Non-adhering Dressing	38
Surgilast® Tubular Elastic Net Retainer	41-43
Ultrafix Self-Adhesive Dressing Retention Tape	57



Skin Care and Bathing

- Incontinence Products
- Hand Sanitizers / Lotions
- Shampoos / Body Washes
- Bath Additives

INCONTINENCE PRODUCTS

In Between® Perineal Spray Skin Cleanser

(Comparable to Aloe Vesta® Perineal/Skin Cleanser and Secura® Personal Cleanser)

In Between® is a versatile product for the quick and effective clean up of all body exudates, particularly those from incontinent episodes. It eliminates odor from urine and feces while gently and thoroughly cleansing the entire perineal skin area. In Between® may be sprayed directly on the skin, soil mass, clothing, bed linens, floors, walls, furniture, carpets, and any other surface that needs to be clean and free of odors. In Between® cleans, soothes, and freshens the skin. It leaves no soap film and does not require rinsing. In Between® contains a special antimicrobial preserving system that aids in eliminating cross infection concerns.



Product	Description	Pkg Unit/Case
S201-1G	1 gallon	4/case
S201-8Z	8 oz	24/case

HAND SANITIZERS / LOTIONS

Mysotrol® No Rinse Hand Sanitizer

Mysotrol® is a convenient and highly practical healthcare personnel hand sanitizer that offers safety and effectiveness by providing germicidal and virucidal action in seconds. It meets the OSHA protocol for a healthcare personnel handwash, reducing the risk of nosocomial infections. As a waterless, no-rinse hand sanitizer, Mysotrol® also evaporates quickly to save time. Lastly, it is formulated without isopropyl alcohol, eliminating the strong alcohol scent common to other hand sanitizers.



Product	Description	Pkg Unit/Case
S505-800ML	800 ml bag	12/case
S505-16Z	16 oz	12/case
S505-4Z	4 oz	24/case

Skin Care Lotion

Derma Sciences' Skin Care Lotion is a greaseless, urea-based hand and body lotion for dry, overworked skin. It is available with or without a pleasant honey almond fragrance. The lotion has a light, greaseless feel with the consistency suitable to healthcare use. Upon initial rub-in, it is quickly absorbed to produce a moisturized, supple and soothing feel without leaving a greasy, tacky, or unpleasant oily after-feel. It can be used in many capacities; as part of an overall skin care program, in Physical Therapy as a lubricant and for massage; in Extended Care Nursing for preventative skin care, and in General Nursing areas for daily care and therapeutic body rubs. The product contains a special antimicrobial preserving system to aid in eliminating cross infection concerns.



Product	Description	Pkg Unit/Case
S401-1G	1 gallon	4/case
S401-8Z	8 oz	24/case
S401-4Z	4 oz	24/case

SHAMPOOS / BODY WASHES

ApriVera®

ApriVera® is an effective odor reducing, non-alkaline body and hair cleanser designed to make the bathing process a positive experience for staff and resident/patient. Its special combination of conditioners and emollients, in a pearlescent lotion base, provides this gentle product with its unique rich lather. ApriVera® is both safe and effective for the tub, shower, bedside basin, and whirlpool bathing. With its proprietary apricot fragrance, ApriVera® leaves both user and facility smelling clean and odor free. The product contains a special antimicrobial preserving system, designed to aid in eliminating cross infection concerns.



Product	Description	Pkg Unit/Case
S107-1G	1 gallon	4/case
S107-800ML	800 ml bag	12/case
S107-16Z	16 oz	12/case
S107-8Z	8 oz	24/case

Bathe Away®

Bathe Away® is a complete body, hair, and scalp cleansing product specially formulated to gently bathe away soil and odor. Bathe Away® contains both glycerin and coconut oil products to minimize the depletion of the skin's own natural oils during cleansing, and chamomile to soothe the skin. It is pH balanced to protect the acid mantle of normal, healthy skin and hair. Its fresh fragrance helps eliminate odors and makes the bathing experience a pleasant one. The product contains a special antimicrobial preserving system to help aid in eliminating cross infection concerns.



Product	Description	Pkg Unit/Case
S101-1G	1 gallon	4/case
S101-16Z	16 oz	12/case
S101-8Z	8 oz	24/case

Primaderm® Shampoo & Body Wash (No Rinse Formula)

Whenever traditional shampooing or bathing is inconvenient or unavailable, Primaderm® Shampoo & Body Wash is an excellent option. Caregivers can shampoo and bathe patients with absolutely no rinsing. Simply apply the product and pat patient dry with a towel. Use of this product can help to reduce laundry and linen costs. The product is pH balanced. *Does not contain natural rubber latex.*



Product	Description	Pkg Unit/Case
S106-1G	1 gallon	4/case
S106-8Z	8 oz	24/case
S106-4Z	4 oz	24/case

Primaderm® Anti-Dandruff Shampoo

Primaderm® Anti-Dandruff Shampoo starts to work immediately to give you fast, effective dandruff control. It contains mild cleansers that leaves your hair looking and feeling clean. It is pH balanced and gentle enough for everyday use. Helps prevent recurrence of flaking and itching associated with dandruff. Does not contain natural rubber latex.



Product	Description	Pkg Unit/Case
S105-8Z	8 oz	24/case

BATH ADDITIVES

Hydro-Soft®

Hydro-Soft® is a concentrated blend of skin emollients and gentle hair cleansers. When used in whirlpool units, Hydro-Soft® functions in two ways. The first function is it reduces the surface tension between the skin and water. This helps the whirlpool unit to effectively soften, loosen, and lift soil, sebum, and flaking tissue from the skin. The second function is it also aids in relieving the drying effect of whirlpool bathing and hydrotherapy treatments by controlling the skin's moisture loss due to prolonged water contact. Hydro-Soft® is completely safe to use in any whirlpool unit and with any type of dispensing system. It is water soluble, and will not leave an oil residue on tub surfaces.



Product	Description	Pkg Unit/Case
---------	-------------	---------------

S301-8Z	8 oz	24/case
---------	------	---------

Disinfectant and Cleanser Hardsurface/ Whirlpool Detergent/Disinfectant

Derma Sciences' Whirlpool/Hard Surface Detergent/Disinfectant is designed specifically for cleaning hard surfaces and whirlpool units in nursing homes, hospitals, and other institutions where controlling the hazards of cross contamination is of prime importance. One ounce of product per gallon of water has been proven effective as a bactericide, mildewstat, sanitizer, virucide, and fungicide. It effectively destroys some of the most problematic organisms in health care facilities including MRSA and *Pseudomonas aeruginosa*.



Product	Description	Pkg Unit/Case
---------	-------------	---------------

S511-1G	1 gallon	4/case
---------	----------	--------



SPECIALTY FIXATION AND CLOSURE DEVICES

- Fixation Devices
- Wound Closure Strips

FIXATION DEVICES

Cath-Strip® Reclosable Catheter Fastener



Cath-Strip® is a Does not contain natural rubber latex, multi-use, reclosable catheter fastener. The exclusive material features flexibility and dynamic adherence which provide excellent tape-to-skin adhesion while minimizing the tension that may result in irritation, blistering, or skin shear. Cath-Strip® is designed to accommodate a variety of catheter fastening requirements including urinary catheters, gastrostomy and jejunostomy tubes, wound drainage systems, central line catheters, and multi-port IVs. *Does not contain natural rubber latex.*

Product	Description	Pkg Unit/Case
CS50	1/pouch, 50/box	50/box

UC Strip® Catheter Tubing Fastener



UC Strip® is a Does not contain natural rubber latex, economical, one-piece catheter / tubing fastener. This flexible adhesive fastener secures urinary and gastrostomy catheter tubing to the patient without the irritation or tissue trauma caused by other tapes. UC Strip® is made of a flexible material with a skin-friendly adhesive. Comfortable for the patient, the fastener allows limited side-to-side movement while it prevents the tubing from being pulled out. *Does not contain natural rubber latex.*

Product	Description	Pkg Unit/Case
UC20	1/pouch, 20/box	20/box

NG Strip® Nasal Tube Fastener



NG Strip® is a Does not contain natural rubber latex securement device for any nasogastric or nasal feeding tube. NG Strip® eliminates the time consuming ritual of securing nasal tubing with strips of irritating tape. NG Strip® features flexibility and dynamic adherence which provide excellent tape-to-skin adhesion while minimizing the tension that may result in irritation, blistering, or skin shear. NG Strip® also minimizes the chance of irritation and tissue necrosis that are common with conventional tapes. *Does not contain natural rubber latex.*

Product	Description	Pkg Unit/Case
NG25	Adult - 1/pouch, 25/box	25/box
NG40	Pediatric - 10/pouch, 4/box	10/pkg, 4 pkgs/box, 40/box
NG50	Adult - 10/pouch, 5/box	10/pkg, 5 pkgs/box, 50/box
NG60	Infant - 10/pouch, 6/box	10/pkg, 6 pkgs/box, 60/box

The Marsupial® Pouch Wound Drainage Container Holder



The Marsupial® pouch is an adjustable belt that fits comfortably around the waste. The pouch is designed to comfortably hold cumbersome wound drainage containers. It allows free movement of the arms and hands. Ideal for patients undergoing breast, chest, or abdominal surgery. The belt and pouch are made of terry cloth. This is a great product designed to keep patients comfortable while helping them to maintain dignity.

Product	Description	Pkg Unit/Case
MAR-01	1 per pouch	1/pkg

WOUND CLOSURE STRIPS

Suture Strip® Plus Flexible Wound Closure Strips

Suture Strip® Plus with maximum dynamic adherence is a Does not contain natural rubber latex premium wound closure strip that provides superior wound security. It is made of macroporous non-woven polyamide with a pressure sensitive adhesive. The proprietary materials allow Suture Strip® Plus to have a perfect balance between elasticity and adherence, closely mimicking the elastic modulus of human skin. The inherent elasticity of the product allows the strips to elongate and yield to the movement of skin, preventing the concentration of forces at the proximal ends and reducing the chances of traction blisters or skin shear. The strips are water resistant and have a skin-friendly adhesive. Suture Strip® Plus is ideal for primary wound closure of surface tissue as well as traumatic wound closure such as skin tears. It is also ideal for instances when aesthetic results are of great importance, such as in plastic surgery. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
TP1100	1/8" x 3"	5/envelope, 50 envelopes/box
TP1101	1/4" x 3"	3/envelope, 50 envelope/box
TP1102	1/4" x 4"	10/envelope, 50 envelope/box
TP1103	1/2" x 4"	6/envelope, 50 envelope/box
TP1104	1/4" x 1 1/2"	6/envelope, 50 envelope/box
TP1105	1" x 5"	4/envelope, 25 envelopes/box

Suture Strip® Wound Closure Strips

Suture Strip® provides superior wound security, but with slightly less flexibility than Suture Strip® Plus, making it an economical alternative. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
TS3100	1/8" x 3"	5/envelope, 50 envelopes/box
TS3101	1/4" x 3"	3/envelope, 50 envelopes/box
TS3102	1/4" x 4"	10/envelope, 50 envelopes/box
TS3103	1/2" x 4"	6/envelope, 50 envelopes/box
TS3104	1/4" x 1 1/2"	6/envelope, 50 envelopes/box

Shur Strip® Wound Closure Strips Sterile

(Comparable to Steri Strips®)

Shur Strip® sterile skin closures provide a wound closure alternative more similar to the market leading brand. They are made of a porous, non-woven backing coated with a pressure-sensitive, hypoallergenic adhesive. These strips are not elastic, keeping wound edges securely together to maximize wound healing. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
DKC81118	¼" x 3"	3/envelope, 50 envelopes/box
DKC81120	¼" x 4"	10/envelope, 50 envelopes/box
DKC81121	½" x 4"	6/envelope, 50 envelopes/box
DKC81119	¼" x 1 ½"	6/envelope, 50 envelopes/box

Shur Strip® Wound Closure Strips Non-Sterile Bulk

(Comparable to Steri Strips®)

These wound closure strips are bulk packaged and are a non-sterile version of our Shur Strip®. **Does not contain natural rubber latex.**



Product	Description	Pkg Unit/Case
DKC71152	¼" x 4"	5/card, 1100 cards, 5500 total
DKC71153	½" x 4"	6/card, 1100 cards, 6600 total

Index: Product Name

Product Name	Page
ALGICELL® Ag CALCIUM ALGINATE DRESSING WITH ANTIMICROBIAL SILVER	11
ALGICELL® CALCIUM ALGINATE DRESSING - STERILE	13
AMNIOEXCEL®	8
AMNIOMATRIX®	8
APRIVERA®	67
AQUASITE® AMORPHOUS HYDROGEL DRESSING	21
AQUASITE® SHEET HYDROGEL - STERILE	20
AQUASITE® IMPREGNATED GAUZE - STERILE	20
BATHE AWAY®	67
BIOGUARD® BARRIER DRESSINGS	13
CATH-STRIP® RECLOSABLE CATHETER FASTENER	72
COMPRESS® ISLAND DRESSING - STERILE	34
COMPRIBAND® SHORT-STRETCH COMPRESSION BANDAGE - DOES NOT CONTAIN NATURAL RUBBER LATEX	48
DERMAGRAN® AMORPHOUS HYDROGEL DRESSING	16
DERMAGRAN®-B HYDROPHILIC WOUND DRESSINGS	15
DERMAGRAN® MOISTURIZING SPRAY	16
DERMAGRAN® OINTMENT	15
DERMAGRAN® WOUND CLEANSER WITH ZINC	24
DISINFECTANT AND CLEANSER HARDSURFACE/ WHIRLPOOL DETERGENT/DISINFECTANT	69
DUBAN® COHESIVE BANDAGE - DOES NOT CONTAIN NATURAL RUBBER LATEX	49
DUCARE® WOVEN GAUZE SPONGES	27
DUDRESS® FILM TOP BARRIER DRESSING - STERILE	34
DUFORE® FOUR-LAYER BANDAGING SYSTEM - DOES NOT CONTAIN NATURAL RUBBER LATEX	50
DUFORM® SYNTHETIC CONFORMING BANDAGE, KNITTED	38
DUPAD® ABDOMINAL PAD	31
DUPAQUE® X-RAY DETECTABLE WOVEN GAUZE SPONGES	28
DUPRESS® BURN DRESSINGS	61
DUSOFT® NON-WOVEN SPONGES	28
DUSOR® ELASTIC BANDAGE - ECONOMY QUALITY - DOES NOT CONTAIN NATURAL RUBBER LATEX	47
DUTEX® 100% COTTON 2 PLY CONFORMING BANDAGE	37
DUTOP® COVER SPONGES	29
FLEXOPLAST® ELASTIC ADHESIVE BANDAGE	48
GAZETEX® 100% COTTON, 6-PLY, FLUFF SPONGE	30
GAZETEX® BANDAGE ROLLS	30
GLEN-SLEEVE® II DOES NOT CONTAIN NATURAL RUBBER LATEX ARM AND LEG PROTECTORS	53
GLEN-SLEEVE® ARM PROTECTORS	53
HEEL & ELBOW PROTECTORS	54
HYDROCELL® ADHESIVE FOAM DRESSING - STERILE	21
HYDROCELL® FOAM DRESSING - STERILE	22
HYDROCELL® THIN ADHESIVE FOAM DRESSING - STERILE	21
HYDRO-SOFT®	69
IN BETWEEN® PERINEAL SPRAY SKIN CLEANSER	65
LAPAROTOMY SPONGES, PREWASHED, X-RAY DETECTABLE	56
THE MARSUPIAL® POUCH WOUND DRAINAGE CONTAINER HOLDER	73
MEDIHONEY® CALCIUM ALGINATE DRESSING	8
MEDIHONEY® GEL DRESSING	7
MEDIHONEY® HCS (HYDROGEL COLLOIDAL SHEET) DRESSING	6
MEDIHONEY® HONEYCOLLOID™ DRESSING	7
MEDIHONEY® PASTE DRESSING	5

MYSOTROL® NO RINSE HAND SANITIZER	67
NG STRIP® NASAL TUBE FASTENER	73
PAK-ITS® GAUZE PACKING STRIPS - PLAIN AND IODOFORM.	33
PAK-ITS® GAUZE PACKING STRIPS - PLAIN AND IODOFORM WITH STRIP DELIVERY SYSTEM	33
PETROLATUM GAUZE NON-ADHERING DRESSING - STERILE	35
PRIMACOL® BORDERED HYDROCOLLOID DRESSING - STERILE	18
PRIMADERM® ANTI-DANDRUFF SHAMPOO.	68
PRIMADERM® DERMAL WOUND CLEANSER - STERILE.	23
PRIMACOL® HYDROCOLLOID DRESSING - STERILE	17
PRIMADERM® SHAMPOO & BODY WASH (NO RINSE FORMULA)	68
PRIMADERM® SKIN PROTECTANT LOTION WITH 3% DIMETHICONE	68
PRIMACOL® THIN HYDROCOLLOID DRESSING - STERILE.	17
PRIMADERM® WOUND CLEANSER - NON-STERILE.	23
PRIMACARE® ELASTIC BANDAGE WITH VELCRO CLOSURE - DOES NOT CONTAIN NATURAL RUBBER LATEX	47
PRIMAPAD® NON-ADHERENT ABSORBENT PAD - STERILE	31
PRIMER® FLEX ELASTIC UNNA BOOT	52
PRIMER® DOES NOT CONTAIN NATURAL RUBBER LATEX MODIFIED UNNA BOOT	51
SHUR-CONFORM® OIL EMULSION NON-ADHERING DRESSING - STERILE	36
SHUR STRIP® WOUND CLOSURE STRIPS, NON-STERILE BULK	75
SHUR STRIP® WOUND CLOSURE STRIPS, STERILE	75
SILVERSEAL® BURN CONTACT DRESSING.	59
SILVERSEAL® BURN GLOVES	60
SILVERSEAL® BURN WRAP	59
SILVERSEAL® HYDROGEL DRESSINGS - STERILE.	12
SILVERSEAL® TUBULAR STRETCH KNIT	60
SILVERSEAL® WOUND CONTACT DRESSING - STERILE.	11
SKIN CARE LOTION	66
SORBACELL® FOAM DRESSING (WITHOUT FILM BACKING) - STERILE	22
SORBAN® UNDERCAST PADDING BANDAGE - DOES NOT CONTAIN NATURAL RUBBER LATEX	49
SURECARE® EYE PADS - STERILE	32
SURGIGRIP® DOES NOT CONTAIN NATURAL RUBBER LATEX TUBULAR ELASTIC SUPPORT BANDAGE	44
SURGIGRIP® TUBULAR ELASTIC SUPPORT BANDAGE	45
SURGILAST® DOES NOT CONTAIN NATURAL RUBBER LATEX TUBULAR ELASTIC DRESSING RETAINER	40
SURGILAST® PRE-CUT TUBULAR ELASTIC DRESSING RETAINER.	41
SURGILAST® TUBULAR ELASTIC DRESSING RETAINER	39
SURGITUBE® DOES NOT CONTAIN NATURAL RUBBER LATEX TUBULAR GAUZE.	43
SURGITUBE® DOES NOT CONTAIN NATURAL RUBBER LATEX TUBULAR GAUZE, 50 YARD ROLLS.	42
SURGITUBE® & SURGIGRIP® TUBULAR GAUZE APPLICATORS	46
SUTURE STRIP® PLUS FLEXIBLE WOUND CLOSURE STRIPS	74
SUTURE STRIP® WOUND CLOSURE STRIPS	74
TCC-EZ® TOTAL CONTACT CASTING SYSTEM	9
TRACHEOTOMY/DRAIN SPONGES - STERILE	29
UC STRIP® CATHETER TUBING FASTENER	72
ULTRAFIX® SELF-ADHESIVE DRESSING RETENTION TAPE.	55
UNNA-PAK®: PRIMER® DOES NOT CONTAIN NATURAL RUBBER LATEX MODIFIED UNNA BOOT AND DUBAN® SELF ADHERENT BANDAGE	51
XEROFORM PETROLATUM DRESSING - STERILE	36
XTRASORB® FOAM SUPER ABSORBENT DRESSING	9
XTRASORB® HCS (HYDROGEL COLLOIDAL SHEET) SUPER ABSORBENT DRESSING.	8
XTRASORB® SUPER ABSORBENT DRESSING	9

Index: Product Name

Product Number	Page	Product Number	Page
1201	58	59246	35
1205	58	59320	35
1206	58	59321	35
1211	58	59345	35
1213	58	59346	35
1214	58	59420	35
1221	58	59421	35
1223	58	59446	35
1231	58	59445	33
1236	58	62036	57
1271	58	62037	57
8002	63	62038	57
8007	63	69101	25
8118S	63	69201	25
8136S	63	69202	25
9320	32	71201	50
9320-48	32	71202	50
9321	32	71412	49
9322	32	71413	49
9333	32	71414	49
9334	32	71416	49
9340	32	71422	49
9345	32	71423	49
9464S	28	71424	49
9467S	28	71426	49
9825	63	72414	52
11903	63	73201	51
11904	63	73202	51
11907	63	73203	51
11908	63	73204	51
12731	58	75102	40
12751	58	75103	40
12931	58	75104	40
31012	7	75106	40
31022	7	75202	40
31045	7	75203	40
31222	7	75204	40
31245	7	75206	40
31422	7	76781	39
31445	7	76782	39
31505	5	76783	39
31515	5	76784	39
31535	5	77781	39
31622	6	77782	39
31644	6	77783	39
31722	6	77784	39
31744	6	79504	51
31805	5	81110	21
31815	5	82023	22
53022	22	82044	22
53044	22	82208	29
53048	22	82312	29
59120	35	82408	29
59121	35	82410	29
59145	35	82412	29
59146	35	82812	29
59220	35	82910	31
59221	35	82911	34
59245	35	83323	33

Product Number	Page
83334	33
83338	33
84122	30
84133	30
84144	30
84244	23
84266	23
84318	24
84322	24
84344	24
84423	23
84444	23
84544	24
84566	24
84644	30
84916	31
84918	31
84943	31
84944	31
85144	19
85166	19
85188	19
85192	20
85222	20
85244	20
85266	20
85288	20
85322	19
85344	19
85366	19
86122	11
86144	11
86188	11
86233	12
86244	12
86266	12
86322	12
86344	12
86388	12
86433	12
86466	12
87059	33
87088	33
87810	33
88022	15
88044	15
88048	15
88100	15
88112	15
88512	14
88522	14
88544	14
88545	14
88548	14
88582	14
89022	36
89041	36
89042	36
89044	36

Product Number	Page
89046	36
89048	36
89066	36
89068	36
89144	36
89166	36
89168	36
89545	11
89569	11
90208	27
90212	27
90312	27
90408	27
90412	27
90812	27
94122	30
94133	30
94144	30
94416	30
94812	30
94943	30
94944	30
97041	10
97042	10
97045	10
97048	10
97110	33
97159	33
97188	33
97208	13
97322	13
97334	13
97338	13
97412	13
97441	13
97831	13
97832	13
97833	13
CS50	74
DG-4	17
DH-03	18
DKC8119	75
DKC20049	37
DKC20053	37
DKC20056	37
DKC20059	37
DKC20062	37
DKC20065	37
DKC20068	38
DKC71149	75
DKC71150	75
DKC71151	75
DKC71152	77
DKC71153	77
DKC71154	75
DKC71164	38
DKC71165	38
DKC77030	37
DKC77034	38

Index: Product Name

Product Number	Page	Product Number	Page
DKC77041	38	GL507	41
DKC77042	38	GL508	41
DKC77043	38	GL509	41
DKC77044	38	GL510	41
DKC81117	75	GL511	41
DKC81118	77	GL512	41
DKC81120	77	GL600	43
DKC81121	77	GL620	43
DM-4	18	GL622	43
DT-4	17	GL641	43
GL105W	45	GL701	41
GL110W	45	GL702	41
GL205F	45	GL703	41
GL205W	45	GL704	41
GL209	44	GL705	41
GL210	44	GL706	41
GL210F	45	GL707	41
GL210W	45	GL708	41
GL211	44	GL709	41
GL212	44	GL710	41
GL219	44	GL711	41
GL220	44	GL712	41
GL221	44	GL713	41
GL222	44	GL714	41
GL223	44	GL715	41
GL224	44	GL720	41
GL225	44	GL722	41
GL227	48	GL724	41
GL230	48	GL726	41
GL231	48	GL750	43
GL231P	48	GL752	43
GL232	48	GL800	56
GL232P	48	GL801	56
GL234	48	GL802	56
GL234P	48	GL899	55
GL236	48	GL900	55
GL241	45	GL1000	55
GL242	45	GL1000B	55
GL242A	45	GL1000WP	55
GL243	45	GL1001	50
GL244	45	GL1002	50
GL245	45	GL1003	50
GL246	45	GL1003R	50
GL247	45	GL1004	50
GL248	45	GL1100B	55
GL251	45	GL2000	55
GL252	45	GL2003	53
GL252A	45	GL2004	53
GL260	45	GL2522	43
GL305F	45	GL3000	55
GL305W	45	GL3000B	55
GL310W	45	GL3001	53
GL405W	45	GL3001C	53
GL501	41	GL3001F	54
GL502	41	GL4001	53
GL503	41	GL4001C	53
GL504	41	GL4001F	54
GL505	41	GLB10	46
GL506	41	GLC10	46

Product Number	Page
GLD10.....	46
GLE10.....	46
GLF10.....	46
GLG10.....	46
GLJ10.....	46
GLK10.....	46
GLL10.....	46
GLLF2501.....	42
GLLF2502.....	42
GLLF2503.....	42
GLLF2504.....	42
GLLF2505.....	42
GLLF2505.5.....	42
GLLF2506.....	42
GLLF2507.....	42
GLLF2508.....	42
GLLF2509.....	42
GLLF2510.....	42
GLLF2511.....	42
GLLF2512.....	42
GLLF2514.....	42
GLALUAPP.....	46
MAR-01.....	75
NG25.....	75
NG40.....	75
NG50.....	75
NG60.....	75
S101-1G.....	69
S101-8Z.....	69
S101-16Z.....	69
S105-8Z.....	70
S106-1G.....	70
S106-4Z.....	70
S106-8Z.....	70
S107-1G.....	69
S107-8Z.....	69
S107-16Z.....	69
S107-800ML.....	69
S120-4Z.....	67
S120-8Z.....	67
S201-1G.....	67
S201-8Z.....	67
S301-8Z.....	71
S401-1G.....	68
S401-4Z.....	68
S401-8Z.....	68
S505-4Z.....	68
S505-800ML.....	68
S505-16Z.....	68
S511-1G.....	71
SPD-03.....	17
SPD-20.....	17
SPD-21.....	17
SPD-24.....	17
SS-BWCD-0404.....	61
SS-BWCD-0408.....	61
SS-BWCD-0816.....	61
SS-BWCD-1616.....	61
SS-BWCD-L1.....	62
SS-BWCD-M1.....	62

Product Number	Page
SS-BKWCD-34.....	62
SS-BKWCD-115.....	62
SS-BKWCD-0510.....	62
SS-BKWCD-1035.....	62
SS-BKWCD-5015.....	62
SS-EBW-0466.....	61
SS-EBW-06108.....	61
SS-HYG-0404.....	14
SS-WCD-0404.....	16
SS-WCD-0412.....	16
SS-WCD-0466.....	16
SS-WCD-1012.....	16
TCC21100.....	9
TCC21114.....	9
TCC21116.....	9
TCC21124.....	9
TCC21126.....	9
TCC23000.....	9
TCC23001.....	9
TCC23002.....	9
TCC23005.....	9
TCC23051.....	9
TCC24000.....	9
TCC24001.....	9
TCC24002.....	9
TCC24005.....	9
TCC24014.....	9
TCC24016.....	9
TCC24024.....	9
TCC24026.....	9
TCC24051.....	9
TCC24214.....	9
TCC24224.....	9
TCC25050.....	9
TCC25051.....	9
TCCEQP.....	12
TCC2MDKK.....	10
TCC2MDKKS.....	10
TCCSAW.....	9
TCCSPR.....	9
TCCVAC.....	9
TCC2ULTR.....	10
TCC2ULTRS.....	10
TP1100.....	76
TP1101.....	76
TP1102.....	76
TP1103.....	76
TP1104.....	76
TP1105.....	76
TS3100.....	76
TS3101.....	76
TS3102.....	76
TS3103.....	76
TS3104.....	76
UC20.....	74
WC-04.....	26
WMB10.....	47
WMC10.....	47
WME10.....	47
WMJ10.....	47

CONTACT US

TELEPHONE ORDERS

U.S. Market

Please call Derma Sciences' customer service at:

(800) 825-4325

We are open Monday through Friday, 9:00-5:30 EST

Canadian Market

Please call Derma Sciences' customer service at:

(800) 387-5302

We are open Monday through Friday, 9:00-5:30 EST

Rest of World Distribution

Please call (609) 514-4744 for all your distribution needs.

WEBSITE

www.dermasciences.com