

Urinalysis

HENRY SCHEIN® Urispec® Plus

Urine test strip analyzer

- CLIA-waived
- Ease of use
- Speed
- Productivity



For use with HENRY SCHEIN® BRAND



URISPEC® 11-Way



URISPEC® 5-Way Plus



URISPEC® 4-Way



URISPEC® GP+A

Automated urinalysis at the point of care

URISPEC® Plus provides dependable urine status results to detect early stages of many diseases such as diabetes, kidney disease and urinary tract infections. Instrument-read results have long proven to be advantageous for both, busy health care professionals and patients. **URISPEC® Plus** readings eliminate the subjectivity of visual colour interpretation. The comprehensive interface options and the optimized printouts minimize risks associated with manual transcriptions. Reliable results can be obtained virtually immediate at the point of care.


URISPEC® Plus is for the use with Henry Schein **URISPEC®** 11-Way, 5-Way Plus, 4-Way and GP+A test strips for up to 11 parameters for the semi-quantitative detection of blood (occult), urobilinogen, bilirubin, protein, nitrite, ascorbic acid, ketones, glucose, pH, specific gravity and leukocytes.

7 convincing features for optimized usability



No-touch operation

After placing the strip on the instrument, the measurement is started automatically. Results can be obtained within a minute, making urinalysis faster and more hygienic. Additionally, the automatic self-test and calibration makes using the urinalysis reader even faster.

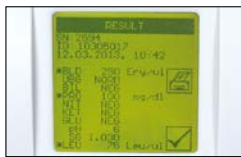

Graphic display

Graphic screen results for all 11 parameters are conveniently displayed at a glance. Pathological results are clearly marked. For further diagnosis results can be printed with the integrated printing unit.



Touchscreen

All functions can be accessed with easy touch-screen operation. Optionally, an external barcode reader or keyboard can also be connected for easy entry of patient information.



Small footprint for mobile operation

Small footprint allows for easy mobility and fit into operatory areas. Irrespective of the instrument size it is still capable of reading 50 test strips per hour.

Details: Depth: 7.9 inches · Width: 6.3 inches · Height: 3.0 inches · Weight: 1.90 lb



Power supply options

The instrument can be run using a 110/240 V autoswitch mains adaptor, or optional battery operation (6 AA batteries). Hence, the instrument is also ideal for near patient testing.



200 patient test results memory

Up to 200 measurements can be stored and retrieved conveniently for further diagnosis when necessary. The instrument can easily be connected via USB or RS 232 to existing laboratory information systems or to PCs.



Easy switch between different test strips

The new **URISPEC® Plus** allows you to read multiple **URISPEC®** test strips, including the proven 11-Way, 5-Way Plus, 4-Way and GP+A strips. With the intuitive selection function, you can effortlessly choose the perfect strip for your individual diagnosis needs.

Ordering information

570-1309	URISPEC® Plus UA Reader	1 unit w/12 boxes of URISPEC® 11-Way
102-2285	URISPEC® 11-Way	100 strips / box
900-7508	URISPEC® 5-Way Plus	100 strips / box
107-7157	URISPEC® 4-Way	100 strips / box
100-6567	URISPEC® GP+A	100 strips / box

