

HENRY SCHEIN®



Contents

<u>Essentials 776HS Pocket Aneroid</u>	3
<u>Pro Series 760HS Pocket Aneroid</u>	4
<u>Pro Plus Series 700HS Pocket Aneroid</u>	5
<u>Pro Plus Series Sphyg/Sprague Kit</u>	6
<u>Pro Plus Series 703HS Pocket Aneroid</u>	7
<u>Pro Plus Series 705GPK-BKHS MultiCuff</u>	8
<u>Pro Plus Series 750HS Clockface Aneroids</u>	9
<u>Inflation Sysytems</u>	10
<u>Institutional Stethoscopes</u>	11
<u>Sprague Rappaport Stethoscopes</u>	12
<u>Pro Plus 603HS Stethoscope</u>	13
<u>Pro Plus 601HS Stethoscope</u>	14
<u>Essentials Teaching Stethoscope</u>	15
<u>Pro 413HS Thermometer</u>	16
<u>351P-HS Disposable Penlight</u>	17



ESSENTIALS 776HS

An economical aneroid sphygmomanometer for today's cost conscious healthcare market.

The Essentials 776HS series features:

- Black enamel 300mmHg CE no-pin stop manometer
- Navy blue nylon cuff with range markings to facilitate selection of correct cuff size; artery label, and gauge holder
- Latex free inflation bladder and bulb
- Chrome-plated deflation valve
- Complete with leatherette zippered carrying case and operating instructions
- Assembled, inspected, and packaged in the U.S.A. from foreign components



PRIVATE LABELING	
COMPONENT	PACKAGING
N/A	1 Pc. Box/Label

Feature	Benefit
Gauge is 100% visually inspected for 0 point, 2% re-inspected for accuracy	To ensure conformance with National Standard ANSI SP 10
Cuff has a properly positioned artery mark at the bladder midpoint	Ensures proper compression of the brachial artery. Most low cost imports fail to meet this requirement from ANSI SP 10
Leak tested Valve	Valves are tested to ensure conformance with ANSI SP 10 for dump rate and leakage
100% inspected, assembled and packaged in the USA	Ensures unmatched consistency, quality, and performance in a value priced BP instrument
5 year calibration, 1 year system warranty	To protect the consumer's investment

ORDERING INFORMATION					
Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
1126076	776CHS	ESSENTIALS 776, Child	Latex Free	ea.	20
1126077	776HS	ESSENTIALS 776, Adult	Latex Free	ea.	20
1126080	776XHS	ESSENTIALS 776, Lrg. Adult	Latex Free	ea.	20



PRO 760HS SERIES

Our most popular aneroid, the 760HS is the model of choice for demanding institutions.

The 760HS series features:

- Contemporary black enamel 300mmHg CE no-pin stop manometer, provides day-in, day-out dependability
- Nylon cuff with proprietary Size Guide™ marking system, prevents mis-cuffing
- PVC latex free inflation bladder withstands heavy use
- Chrome-plated brass air release valve
- Leatherette carrying case with durable nylon zipper
- Assembled, inspected, and packaged in the U.S.A. from foreign components



PRIVATE LABELING	
COMPONENT	PACKAGING
Gauge & Cuff	1 Pc. Box/Label

Feature	Benefit
Nissei Made Aneroid Gauge EXCLUSIVELY - The MOST experienced and respected Asian Gauge manufacturer	Unsurpassed accuracy, quality, and reliability. Our gauges are randomly tested to 30,000 cycles, 3x MORE than required under the National Standard (ANSI SP 10)
Gauge is 100% visually inspected for 0 point, 5% re-inspected for accuracy.	To ensure conformance with National Standard ANSI SP 10
Accurate to +/-2mmHg	Exceeds industry requirements by 33%
Gauge features extended graduation marks to prevent parallax error	Makes it easier to read from oblique angles
Cuff has a properly positioned artery mark at the bladder midpoint	Ensures proper compression of the brachial artery. Most low cost imports fail to meet this requirement from ANSI SP 10
Size Guide marking system (points out when to choose larger or smaller cuffs, properly positioned artery mark)	Prevents mis-cuffing which can lead to measurement error
Leak tested Valve	Valves are tested to ensure conformance with ANSI SP 10 for dump rate and leakage.
Premium filter equipped valve and end valve	Large valve is more ergonomic, filters prevent damaging dust build up
100% inspected, assembled and packaged in the USA	Ensures unmatched consistency, quality, and performance in a value priced BP instrument
20 year calibration, 3 year system warranty	To protect the consumer's investment

ORDERING INFORMATION					
Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
1126071	760-10SABKHS	PRO 760, Small Adult, Black	Latex-Free	ea.	20
1126072	760-11ABKHS	PRO 760, Adult, Black	Latex-Free	ea.	20
1126073	760-9CBKHS	PRO 760, Child, Black	Latex-Free	ea.	20
1126074	760-13TBKHS	PRO 760, Thigh, Black	Latex-Free	ea.	12
1126075	760-12XBKHS	PRO 760, Lrg. Adult, Black	Latex-Free	ea.	20



PRO PLUS 700HS

Our finest pocket aneroid, and it is the instrument of choice for the discriminating professional.

The PRO PLUS 700HS series features:

- Chrome-plated precision crafted 300mmHg CE no-pin stop manometer
- Luminescent dial features bold numerals and extended graduation marks for easy reading in all light conditions and from virtually any angle
- Nylon cuff with proprietary Size Guide™ marking system, prevents mis-cuffing
- PVC Latex free inflation bladder withstands heavy use
- Latex free inflation bulb fitted with filter screen end valve to reduce dust build-up
- Chrome-plated brass air release valve with filter screen
- Leatherette carrying case with durable nylon zipper
- Assembled, inspected, and packaged in the U.S.A. from foreign components

PRIVATE LABELING	
COMPONENT	PACKAGING
Gauge & Cuff	2 Pc. Box/Label



Feature	Benefit
Nissei Made Aneroid Gauge EXCLUSIVELY - The MOST experienced and respected Asian Gauge manufacturer	Unsurpassed accuracy, quality, and reliability. Our gauges are randomly tested to 30,000 cycles, 3x MORE than required under the National Standard (ANSI SP 10)
Gauge is 100% visually inspected for 0 point, 5% re-inspected for accuracy.	To ensure conformance with National Standard ANSI SP 10
Accurate to +/-2mmHg	Exceeds industry requirements by 33%
Gauge features extended graduation marks to prevent parallax error	Makes it easier to read from oblique angles
Gauge is chrome plated	More durable finish
Gauge features Luminescent dial	Easier to read in low light conditions
Cuff has a properly positioned artery mark at the bladder midpoint	Ensures proper compression of the brachial artery. Most low cost imports fail to meet this requirement from ANSI SP 10
Size Guide marking system (points out when to choose larger or smaller cuffs, properly positioned artery mark)	Prevents mis-cuffing which can lead to measurement error
Unit features Latex Free bulb and bladder	Eliminates risk of latex reaction
Leak tested Valve	Valves are tested to ensure conformance with ANSI SP 10 for dump rate and leakage.
Premium filter equipped valve and end valve	Large valve is more ergonomic, filters prevent damaging dust build up
100% inspected, assembled and packaged in the USA	Ensures unmatched consistency, quality, and performance in a value priced BP instrument
Lifetime calibration, 3 year system warranty	To protect the consumer's investment

ORDERING INFORMATION					
Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
1126064	700-10SABKHS	PRO PLUS 700, Sm. Adult, Black	Latex-Free	ea.	20
1126065	700-11ABKHS	PRO PLUS 700, Adult, Black	Latex-Free	ea.	20
1126066	700-12XABKHS	PRO PLUS 700, Lg. Adult, Black	Latex-Free	ea.	20



PRO PLUS Sphyg/Sprague Kit

A professional kit that combines our Pro Pocket Aneroid with our Pro Sprague stethoscope in a color-coordinated premium nylon case.

The PRO PLUS Sphyg/Sprague Kit features:

- Contemporary black enamel 300mmHg manometer
- Nylon cuff with proprietary Size Guide™ marking system prevents mis-cuffing
- Chrome-plated deflation valve
- Color-matched Sprague stethoscope with 22" dual tubing, adjustable binaurals, and interchangeable chestpiece and accessory kit
- Extra-large, color-matched dual-pocket nylon carry case, measures 10" x 6" x 3"
- Inspected, Assembled and Packaged in the U.S.A. from foreign components



PRIVATE LABELING	
COMPONENT	PACKAGING
Gauge & Cuff	Label

Feature	Benefit
Nissei Made Aneroid Gauge EXCLUSIVELY - The MOST experienced and respected Asian Gauge manufacturer	Unsurpassed accuracy, quality, and reliability. Our gauges are randomly tested to 30,000 cycles, 3x MORE than required under the National Standard (ANSI SP 10)
Gauge is 100% visually inspected for 0 point, 5% re-inspected for accuracy.	To ensure conformance with National Standard ANSI SP 10
Accurate to +/-1mmHg	Exceeds industry requirements by 67%
Gauge features extended graduation marks to prevent parallax error	Makes it easier to read from oblique angles
Cuff has a properly positioned artery mark at the bladder midpoint	Ensures proper compression of the brachial artery. Most low cost imports fail to meet this requirement from ANSI SP 10
Size Guide marking system (points out when to choose larger or smaller cuffs, properly positioned artery mark)	Prevents mis-cuffing which can lead to measurement error
Unit features Latex Free bulb and bladder	Eliminates risk of latex reaction
Leak tested Valve	Valves are tested to ensure conformance with ANSI SP 10 for dump rate and leakage.
100% inspected, assembled and packaged in the USA	Ensures unmatched consistency, quality, and performance in a value priced BP instrument
Lifetime calibration, 3 year system warranty	To protect the consumer's investment

ORDERING INFORMATION					
Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
1126780	768-641-11ABKHS	PRO PLUS Sphyg/Sprague Kit, Adult, Black	Latex-Free	ea.	10
1126781	768-641-11ARBHS	PRO PLUS Sphyg/Sprague Kit, Adult, Royal Blue	Latex-Free	ea.	10
1126782	768-641-11ATLHS	PRO PLUS Sphyg/Sprague Kit, Adult, Teal	Latex-Free	ea.	10
1126784	768-641-11AVHS	PRO PLUS Sphyg/Sprague Kit, Adult, Purple	Latex-Free	ea.	10



PRO PLUS 703HS

A versatile palm style aneroid sphygmomanometer that combines precision crafted performance with rugged durability.

The Pro Plus 703HS series features:

- Revolutionary palm style manometer with integral latex-free bulb and trigger style air release valve for one-handed convenience
- Precision crafted 300mmHg CE no-pin stop movement tested to 200,000 cycles
- Movement housed in a zinc safety cage for unsurpassed durability
- Ambidextrous infinite control trigger release air valve
- Nylon cuff with proprietary Size Guide™ marking system, prevents mis-cuffing
- PVC latex free inflation bladder withstands heavy use
- Quick release luer slip connector
- Complete with leatherette carrying case with self repairing nylon zipper
- Assembled, inspected, and packaged in the U.S.A. from foreign components



PRIVATE LABELING	
COMPONENT	PACKAGING
Gauge & Cuff	2 Pc. Box/Label

Feature	Benefit
Nissei Made Aneroid Gauge EXCLUSIVELY - The MOST experienced and respected Asian Gauge manufacturer	Unsurpassed accuracy, quality, and reliability. Our gauges are randomly tested to 30,000 cycles, 3x MORE than required under the National Standard (ANSI SP 10)
Gauge is 100% visually inspected for 0 point, 5% re-inspected for accuracy.	To ensure conformance with National Standard ANSI SP 10
Accurate to +/-1mmHg	Exceeds industry requirements by 67%
Gauge features extended graduation marks to prevent parallax error	Makes it easier to read from oblique angles
Gauge features Luminescent dial	Easier to read in low light conditions
Gauge movement housed in zinc cage	Unsurpassed shock resistance
Cuff has a properly positioned artery mark at the bladder midpoint	Ensures proper compression of the brachial artery. Most low cost imports fail to meet this requirement from ANSI SP 10
Size Guide marking system (points out when to choose larger or smaller cuffs, properly positioned artery mark)	Prevents mis-cuffing which can lead to measurement error
100% inspected, assembled and packaged in the USA	Ensures unmatched consistency, quality, and performance in a value priced BP instrument
Lifetime calibration, 3 year system warranty	To protect the consumer's investment

ORDERING INFORMATION					
Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
1126067	703-11ABKHS	PRO PLUS 703, Adult, Black	Latex-Free	ea.	20

PRO PLUS 705GPK-BKHS

A versatile palm style aneroid 3-cuff sphygmomanometer system that combines precision crafted performance, rugged durability and affordability.



The Pro Plus 705GPK-BKHS series features:

- Revolutionary palm style manometer with integral latex-free bulb and trigger style air release valve for one-handed convenience
- Precision crafted 300mmHg CE no-pin stop movement tested to 50,000 cycles
- Movement housed in a zinc safety cage for unsurpassed durability
- Ambidextrous infinite control trigger release air valve
- 3 Nylon cuffs with proprietary Size Guide™ marking system, prevents mis-cuffing (Includes Small Adult, Adult, and Large Adult sizes)
- PVC latex free inflation bladder withstands heavy use
- Quick release luer slip connector for rapid selection of appropriate cuff
- Four section, smoke gray plexiglass tray features small 14" x 3.5" footprint
- Tray may be wall mounted with included hardware
- Assembled, inspected, and packaged in the U.S.A. from foreign components

PRIVATE LABELING	
COMPONENT	PACKAGING
Gauge & Cuff	Label



Feature

Benefit

Nissei Made Aneroid Gauge EXCLUSIVELY - The MOST experienced and respected Asian Gauge manufacturer	Unsurpassed accuracy, quality, and reliability. Our gauges are randomly tested to 30,000 cycles, 3x MORE than required under the National Standard (ANSI SP 10)
Gauge is 100% visually inspected for 0 point, 5% re-inspected for accuracy.	To ensure conformance with National Standard ANSI SP 10
Accurate to +/-1mmHg	Exceeds industry requirements by 67%
Gauge features extended graduation marks to prevent parallax error	Makes it easier to read from oblique angles
Gauge features Luminescent dial	Easier to read in low light conditions
Gauge movement housed in zinc cage	Unsurpassed shock resistance
Cuff has a properly positioned artery mark at the bladder midpoint	Ensures proper compression of the brachial artery. Most low cost imports fail to meet this requirement from ANSI SP 10
Size Guide marking system (points out when to choose larger or smaller cuffs, properly positioned artery mark)	Prevents mis-cuffing which can lead to measurement error
100% inspected, assembled and packaged in the USA	Ensures unmatched consistency, quality, and performance in a value priced BP instrument
Lifetime calibration, 3 year system warranty	To protect the consumer's investment

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M
1126778	705GPK-BKHS	PRO PLUS GP Multicuff Kit, Black	Latex-Free	ea.



750W-11ABKHS



752M-11ABKHS

Pro Plus 750HS Series Clockface Aneroids

A large face aneroid available in 2 configurations that offer the superb precision, styling, and dependability.

The Pro Plus 750HS series features:

- Large, easy-to-read 6" luminescent dial for high visibility in virtually any light
- 300mmHg CE no-pin stop movement withstands the heaviest use
- Inflation system includes black nylon cuff, latex free inflation bladder, deluxe bulb and valve, and eight foot length (when extended) of coiled latex free tubing
- Contemporary Euro-style ivory casing compliments any decor
- Complete with all necessary hardware and assembly instructions
- Available in Wall or Mobile configurations
- Assembled, inspected, and packaged in the U.S.A. from foreign components

The Pro Plus 750W Wall Aneroid includes:

- Wall bracket that swivels through 110° rotation for easy viewing with integral cuff storage

The Pro Plus 752M Mobile Aneroid includes:

- Height adjustable Spider-Leg™ caster stand with range of 45" to 54" (from bottom of base to top of unit)
- Integral ABS cuff storage compartment
- Manometer angled back 15° for easy viewing
- 18" base diameter

**Some assembly required*

PRIVATE LABELING	
COMPONENT	PACKAGING
Gauge & Cuff	Label



Feature	Benefit
Nissei Made Aneroid Gauge EXCLUSIVELY - The MOST experienced and respected Asian Gauge manufacturer	Unsurpassed accuracy, quality, and reliability. Our gauges are randomly tested to 30,000 cycles, 3x MORE than required under the National Standard (ANSI SP 10)
100% visually inspected for 0 point, 15% re-inspected for accuracy.	To ensure conformance with National Standard ANSI SP 10
Gauge features dampening mechanism	Reduces "needle bounce" and promotes longer life.
Gauge features Luminescent dial	Easier to read in low light conditions
Gauge crystal specially designed to seal out dust	Prolongs life of movement.
has a properly positioned artery mark at the bladder midpoint	Ensures proper compression of the brachial artery. Most low cost imports fail to meet this requirement from ANSI SP 10
Size Guide marking system (points out when to choose larger or smaller cuffs, properly positioned artery mark)	Prevents mis-cuffing which can lead to measurement error
Leak tested Valve	Valves are tested to ensure conformance with ANSI SP 10 for dump rate and leakage.
Premium filter equipped valve and end valve	Large valve is more ergonomic, filters prevent damaging dust build up
Unit features Latex Free bulb and bladder	Eliminates risk of latex reaction
100% inspected, assembled and packaged in the USA	Ensures unmatched consistency, quality, and performance in a value priced BP instrument
Lifetime calibration, 3 year system warranty	To protect the consumer's investment

ORDERING INFORMATION					
Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
1126069	750W-11ABKHS	PRO PLUS 750W, Adult, Black	Latex-Free	ea.	12
1126070	752M-11ABKHS	PRO PLUS 752M, Adult, Black	Latex-Free	ea.	4



840-11ABK-2HS



872NHS

PRIVATE LABELING	
COMPONENT	PACKAGING
Gauge & Cuff	1 Pc. Box/Label

CUFF & BLADDERS

- Durable nylon construction
- Exclusive Size Guide™ marking system
- Latex-free inflation bladder with 55cm tubing length
- For use with most palm, traditional pocket, clock, or mercury style aneroid instruments
- Individually boxed
- Three Year Warranty
- Luer connector not included

BULB & VALVE COMBINATIONS

- Completely latex-free
- For use with most aneroid or mercurial BP instruments
- End valve included
- Individually packaged
- Three Year Warranty

Feature

Benefit

Nissei Made Aneroid Gauge EXCLUSIVELY - The MOST experienced and respected Asian Gauge manufacturer	Unsurpassed accuracy, quality, and reliability. Our gauges are randomly tested to 30,000 cycles, 3x MORE than required under the National Standard (ANSI SP 10)
Gauge is 100% visually inspected for 0 point, 15% re-inspected for accuracy.	To ensure conformance with National Standard ANSI SP 10
Accurate to +/-2mmHg	Exceeds industry requirements by 33%
Pro Plus Gauge features extended graduation marks to prevent parallax error	Makes it easier to read from oblique angles
Cuff has a properly positioned artery mark at the bladder midpoint	Ensures proper compression of the brachial artery. Most low cost imports fail to meet this requirement from ANSI SP 10
Size Guide marking system (points out when to choose larger or smaller cuffs, properly positioned artery mark)	Prevents mis-cuffing which can lead to measurement error
Leak tested Valve	Valves are tested to ensure conformance with ANSI SP 10 for dump rate and leakage.
Premium filter equipped valve and end valve	Large valve is more ergonomic, filters prevent damaging dust build up
100% inspected, assembled and packaged in the USA	Ensures unmatched consistency, quality, and performance in a value priced BP instrument

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
9004841	800HS	Gauge for Pro Plus 700 Series	Latex Free	ea.	12
9004836	802HS	Gauge for Pro 760 Series	Latex Free	ea.	12
9004851	809EHS	Gauge for Essentials 776 Series	Latex Free	ea.	12
1126081	845-10SABK-2HS	Cuff & Bladder 2 Tube, Sm. Adult, Black	Latex Free	ea.	15
1126082	845-11ABK-2HS	Cuff & Bladder 2 Tube, Adult, Black	Latex Free	ea.	15
1126086	845-12XBK-2HS	Cuff & Bladder 2 Tube, Lg. Adult, Black	Latex Free	ea.	15
1126085	845-13TBK-2HS	Cuff & Bladder 2 Tube, Thigh, Black	Latex Free	ea.	16
1126084	845-9CBK-2HS	Cuff & Bladder 2 Tube, Child, Black	Latex Free	ea.	15
1126087	845-10SABK-1HS	Cuff & Bladder 1 Tube, Sm. Adult, Black	Latex Free	ea.	15
1126088	845-11ABK-1HS	Cuff & Bladder 1 Tube, Adult, Black	Latex Free	ea.	15

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
1126092	845-12XBK-1HS	Cuff & Bladder 1 Tube, Lg. Adult, Black	Latex Free	ea.	15
1126091	845-13TBK-1HS	Cuff & Bladder 1 Tube, Thigh, Black	Latex Free	ea.	16
1126090	845-9CBK-1HS	Cuff & Bladder 1 Tube, Child, Black	Latex Free	ea.	15
1126094	865-10SABKHS	Inflation System, Sm. Adult, Black	Latex Free	ea.	15
1126095	865-11ABKHS	Inflation System, Adult, Black	Latex Free	ea.	15
1126096	865-12XBKHS	Inflation System, Lg. Adult, Black	Latex Free	ea.	15
1126097	865-13TBKHS	Inflation System, Thigh, Black	Latex Free	ea.	16
1126093	865-9CBKHS	Inflation System, Child, Black	Latex Free	ea.	15
1126098	872NHS	Deluxe Bulb & Valve	Latex Free	ea.	100
1126099	875NHS	Standard Bulb & Valve	Latex Free	ea.	100

ESSENTIALS™

Henry Schein Essentials combine lightweight comfort, sensitive acoustics, and handsome styling, with truly affordable pricing.

The ESSENTIALS Series features:


- Choice of diaphragm only or combination chestpiece
- Lightweight aluminum chestpiece construction
- Non-chill diaphragm rim for patient comfort
- Proprietary diaphragm enhances acoustic response by 50% over traditional designs
- Adjustable aluminum binaurals with PVC eartips
- Flexible 22" PVC tubing
- Spare pair of black mushroom eartips included
- Inspected and packaged in the U.S.A.
- One Year Warranty

ESSENTIALS 660 series features:

- Non-chill diaphragm only nurse style chestpiece (1³/₄" dia.)
- Ideally suited for blood pressure measurement
- Weighs 2.5 oz. / Overall length 32"

ESSENTIALS 670 series features:

- Combination diaphragm/bell chest piece for wide frequency response
- Diaphragm diameter 1³/₄", bell diameter 3³/₄"
- Non-chill bell ring and diaphragm for patient comfort
- Weighs 3 oz. / Overall length 32.5"

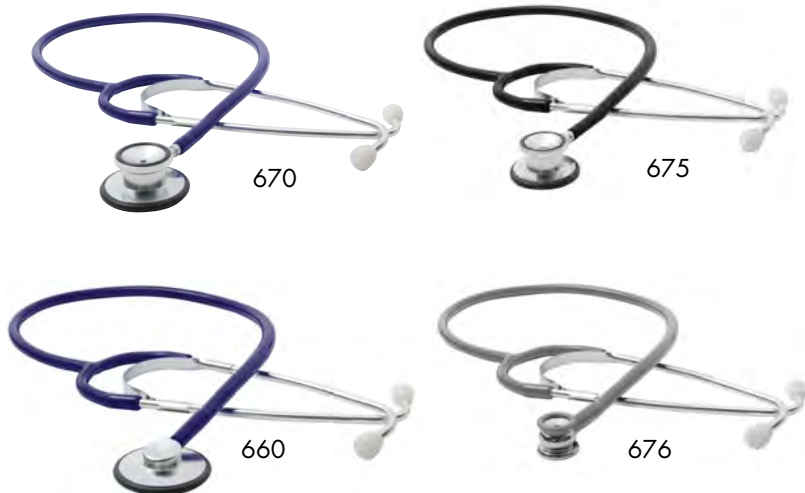

ESSENTIALS 675, Pediatric Dual-head
ESSENTIALS 676, Infant Dual-head

PRIVATE LABELING	
COMPONENT	PACKAGING
N/A	1 Pc. Box/Label

Feature	Benefit
660, 670, 675, 676 Exclusive component sourcing	Unsurpassed quality and consistency
660, 670, 675, 676 Aluminum chestpiece and binaurals	Outstanding acoustics in an ultra-lightweight design
660, 670 Proprietary diaphragm/rim design	Improves acoustic performance by up to 50% over competing brands with threaded rims
660, 670, 675, 676 Non chill rim	For patient comfort
660, 670, 675, 676 Binaural up to 1/2" longer than competing brands	More practitioner comfort
660, 670, 675, 676 Y tubing maintained at tight specs	For greater consistency
660, 670, 675, 676 Spare pair of mushroom eartips	For greater comfort
660, 670, 675, 676 Inspected and packaged in the USA	For more consistent quality
660, 670, 675, 676 1 year warranty	To protect the consumer's investment

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
9004818	660BKHS	ESSENTIALS Nursescope, Black	Latex-Free	ea.	50
9004819	660GHS	ESSENTIALS Nursescope, Gray	Latex-Free	ea.	50
9004820	660RBHS	ESSENTIALS Nursescope, Royal Blue	Latex-Free	ea.	50
9004821	660RHS	ESSENTIALS Nursescope, Red	Latex-Free	ea.	50
9004809	670BKHS	ESSENTIALS Dualhead, Black	Latex-Free	ea.	50

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
9004810	670GHS	ESSENTIALS Dualhead, Gray	Latex-Free	ea.	50
9004811	670RBHS	ESSENTIALS Dualhead, Royal Blue	Latex-Free	ea.	50
9004812	670RHS	ESSENTIALS Dualhead, Red	Latex-Free	ea.	50
9004813	675BKHS	ESSENTIALS Pediatric, Black	Latex-Free	ea.	50
9004814	676GHS	ESSENTIALS Infant, Gray	Latex-Free	ea.	50



641



647

PRIVATE LABELING	
COMPONENT	PACKAGING
N/A	2 Pc. Box/Label



PRO SPRAGUE RAPPAPORT

Our exclusive design features an improved valve mechanism that virtually eliminates acoustic leakage and a chestpiece key to facilitate after market service.

The PRO Spragues feature:

- Threaded chestpiece drum of chrome plated zinc permits the use of five interchangeable chestpiece fittings - an adult and pediatric diaphragm, and three bells (1 1/4", 1", 3/4") - for complete diagnostic capability in most patients
- Latex-free tubing
- Proprietary valve mechanism eliminates acoustic leakage
- Adjustable chrome-plated brass binaurals
- Complete with an accessory pouch containing the bell chestpieces, chestpiece key, 2 spare diaphragms, one pair of mushroom eartips, and one pair of silicon eartips
- Convenient 22" tubing length in choice of 29 fashion tubing colors
- Weighs 7 oz. (Traditional) 5.7 oz. (Single Tube)
- Overall length 30"
- 5 Year Warranty
- Inspected and packaged in the U.S.A.

Feature	Benefit
Exclusive Source	Unsurpassed quality and consistency
Chestpiece features patented valve mechanism	Virtually eliminates acoustic leakage
Convertible chestpiece can accommodate any combination of two – adult diaphragm, pediatric diaphragm, adult, pediatric, and infant bell	For the ultimate in versatility. Great scope for students, or EMS
Gasket equipped chestpiece drum	Improves leakage and improves security of chestpiece fittings
Chestpiece “key”	Facilitates removal or worn or broken diaphragms reducing the need for factory service
Dual tube design (Traditional Only)	Improves acoustic separation
Two tubing clips (Traditional Only)	Reduce incidental rubbing
Extra long binaurals	Improved practitioner comfort
Spare diaphragms	Permits self service
Includes 3 sets of eartips	For user comfort
Individually inspected and packaged in the USA	For more consistent quality
5 year warranty	To protect the consumer’s investment

ORDERING INFORMATION					
Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
9004822	641BKHS	PRO Sprague, 22", Black	Latex-Free	ea.	20
9004823	641NHS	PRO Sprague, 22", Navy	Latex-Free	ea.	20
9004824	641NOHS	PRO Sprague, 22", Neon Orange	Latex-Free	ea.	20
9004825	641RBHS	PRO Sprague, 22", Royal Blue	Latex-Free	ea.	20
9004826	647BKHS	PRO 1-Tube Sprague, 22", Black	Latex-Free	ea.	20

PRO PLUS 603HS

One of our most popular acoustic scopes, the 603 combines outstanding acoustic performance with rugged durability.



The Pro Plus 603HS features:

- Combination adult size chestpiece precisely machined to exacting tolerances from surgical stainless steel for outstanding performance and rugged durability
- Extra large bell (1³/₈" diameter) for unsurpassed low frequency response
- Ultra-sensitive diaphragm (1⁷/₈" diameter) for greater amplification and crisper high frequency response
- Color coordinated non-chill bell and diaphragm retaining rim for patient comfort
- Stainless steel binaurals with double leaf internal spring fixed at 15° angle for greater comfort
- Reinforcing yoke molded into flexible 22" PVC "Y" tubing
- Deluxe PVC eartips (in small and large) for the ultimate in wearing comfort and acoustic seal
- Spare diaphragm included
- Available in 3 striking colors for departmental coding or personal identification
- Scope ID Tag included
- Weighs 6.5 oz.
- Overall length 31"
- Lifetime Warranty

PRIVATE LABELING	
COMPONENT	PACKAGING
N/A	2 Pc. Box/Label



Feature

Benefit

Exclusive source has been producing premium scopes longer than any foreign mfg	Unsurpassed quality and consistency
CNC machined chestpieces	Manufacturing and acoustic consistency
Stainless steel chestpiece AND binaurals	The ultimate in durability and acoustic performance
Non chill bell and diaphragm rims color coordinated to the tubing	For patient comfort and striking style
Binaurals fixed at 15°	Ergonomically correct position for better acoustic seal.
Dual leaf internal binaural spring	Tested to 10,000 cycles for durability
Reinforcing yoke at tube Y	Reinforces stress point for greater durability
Combination diaphragm and bell chestpiece	Allows attenuating of desired frequencies. Diaphragm for high frequency sounds, bell for lower frequency sounds
Threaded eartips	Stay on securely – for the ultimate in wearing comfort and acoustic seal
100% inspected in the USA	To ensure quality, and performance of EVERY scope
Lifetime warranty	To protect the consumer's investment

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
9004816	603BDHS	PRO PLUS 603, Adult, Burgundy	Latex-Free	ea.	10
9004815	603BKHS	PRO PLUS 603, Adult, Black	Latex-Free	ea.	10
9004817	603NHS	PRO PLUS 603, Adult, Navy	Latex-Free	ea.	10

PRO PLUS 601HS

A convertible scope that combines the versatility of our Sprague with acoustic excellence.



PRIVATE LABELING

COMPONENT	PACKAGING
Chestpiece	2 Pc. Box/Label



The Pro Plus 601HS features:

- Threaded chestpiece drum allows rapid conversion from a pediatric diaphragm to an adult bell
- Chestpiece precisely machined to exacting tolerances from surgical stainless steel for outstanding performance and rugged durability
- Extra deep bell for unsurpassed low frequency response
- Ultra-sensitive diaphragms for greater amplification and crisper high frequency transmission
- Non-chill bell and diaphragm retaining rims for patient comfort
- Large internal bore stainless steel binaurals are fixed at a 15° angle for greater comfort
- Custom molded bi-lumen 19" PVC tubing stays soft and flexible
- Extra-thick tubing walls minimize extraneous noise
- Double leaf internal spring construction reduces bulk and prevents kinking when folded
- Complete with one pair of PVC eartips and deluxe PVC eartips in small and large sizes for the ultimate in wearing comfort and acoustic seal
- Complete with an accessory case containing the adult bell, 2 spare diaphragms, and two sets of eartips
- Latex-free
- Scope ID Tag included
- Weighs 7.6 oz. (adult bell) and 7.4 oz. (pediatric diaphragm)
- Overall length 28"
- Limited Lifetime Warranty on metal parts

Feature

Benefit

Exclusive source has been producing premium scopes longer than any foreign mfg	Unsurpassed quality and consistency
CNC machined chestpieces	Manufacturing and acoustic consistency
Stainless steel chestpiece AND binaurals	The ultimate in durability and acoustic performance
Non chill bell and diaphragm rims color coordinated to the tubing	For patient comfort and striking style
Binaurals fixed at 15°	Ergonomically correct position for better acoustic seal
Dual leaf internal binaural spring	Tested to 10,000 cycles for durability
Reinforcing yoke at tube Y	Reinforces stress point for greater durability
Combination diaphragm/ bell chestpiece	Allows attenuating of desired frequencies. Diaphragm for high frequency sounds, bell for lower frequency sounds
Threaded eartips	Stay on securely – for the ultimate in wearing comfort and acoustic seal
100% inspected in the USA	To ensure quality, and performance of EVERY scope
Lifetime warranty	To protect the consumer's investment

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
9007457	601BKHS	PRO PLUS 601, Convertible Cardiology, Black	Latex-Free	ea.	5

ESSENTIALS DUALHEAD TEACHING STETHOSCOPE

A training aid that allows student and instructor to listen simultaneously.



The Essentials Dualhead Scope features:

- Combination diaphragm/bell chestpiece
- Two sets of adjustable aluminum binaurals
- Two sets of flexible 22" and one straight 10" length of black PVC tubing
- Two spare pairs of black mushroom eartips included
- Weighs 6.5 oz. (184 gr.)
- Overall length 44.5"
- Inspected, and Packaged in the U.S.A.
- One-year warranty



PRIVATE LABELING

COMPONENT	PACKAGING
N/A	Label

Feature

Benefit

Exclusive component sourcing	Unsurpassed quality and consistency
Aluminum chestpiece and binaurals	Outstanding acoustics in an ultra-lightweight design
Proprietary diaphragm/rim design	Improves acoustic performance by up to 50% over competing brands with threaded rims
Non chill rim	For patient comfort
Binaural up to 1/2" longer than competing brands	More practitioner comfort
Y tubing maintained at tight specs	For greater consistency
Spare pair of mushroom eartips	For greater comfort
Inspected and packaged in the USA	For more consistent quality
1 year warranty	To protect the consumer's investment

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
1126777	671HS	ESSENTIALS, Teaching Dualhead Scope, Black	Latex-Free	ea.	30

PRO 413HS THERMOMETER



413HS

Features:

- Compact size - under 5" long, weighs just 11 grams
- Auto-off function conserves battery life
- Last reading memory
- Water resistant
- Audible tone
- Switchable C°/F° Scale
- Range 90°F-109.9°F ±.2°F /32°C - 43.9°C ±.1°C
- Replaceable 1.55v (LR41) battery provides 1,500 measurements
- Plastic case keeps instrument sanitary
- Boxed unit complete with 5 sheaths
- 2 Year Warranty



PRIVATE LABELING	
COMPONENT	PACKAGING
413 - Thermometer	1 Pc. Box/Label

Feature	Benefit
Voted best pen type thermometer mfg in 2006 in USA by Frost and Sullivan	Confident in our product quality
Meets ASTM E1112 standards	For accuracy and performance
C°/F° switchable	For any application
Water-resistant	Allows cleaning without risk of damage
Last reading memory	Recalls last reading in case unit powers off before practitioner has had a chance to record reading
Auto off	Conserves battery life
Beeps when measurement complete	Signals practitioner that measurement is completed

ORDERING INFORMATION					
Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
9004856	413HS	Pro Digital Thermometer, °F/°C	Latex-Free	ea.	144



Disposable Penlight

The 351P-HS series features:

- Sealed concentric beam, disposable diagnostic penlight
- Activated by depressing pocket clip
- Pupil gauge
- Measures 5" long, 1/2" diameter
- Sold in 6 pack

PRIVATE LABELING

COMPONENT	PACKAGING
Barrel	Sleeve (6)

ORDERING INFORMATION

Schein Item#	ADC Item#	Product Description	Latex Content	U/M	Cs Pk
9004970	351P-HS	Disposable Penlight 6/pack	Latex-Free	pkg.	60