

PROCARE[®]

REMEDY PRO[™]

Off Loading Shoe

Reduce forefoot pressure
& **promote** faster healing



The Remedy Pro™ Off Loading Shoe helps to reduce forefoot pressure and is indicated for: Forefoot Trauma; Post Surgical Healing; and Diabetic Ulcers or Wounds

The wedge design of the Remedy Pro Shoe promotes healing by transferring weight behind the MP joints and reducing pressure on the forefoot. The off-loading effect allows wounds under the metatarsal heads or toes and post-forefoot trauma to heal faster.



Easy Forefoot Closure
Wide, removable forefoot closure strap can be easily adjusted for maximum security without buckles.



15° Wedge design
Reduces forefoot pressure by shifting weight to midfoot and heel.




Square Toe Box
Protects the toes and provides improved fit for either left or right foot.



Durable Outer Sole
Mid foot & heel tread provide traction where patient needs it most for safety and support.

Durable, Padded Liner & Adjustable Ankle Strap
Secure ankle strap maintains heel in proper position, plus a soft padded liner help reduce friction and provide maximum patient comfort and protection.

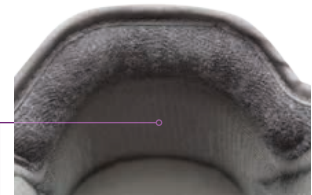

Seamless, Padded Heel
No seam at the heel eliminates friction that can cause wounds or ulcerations.

Part Number	Shoe Size Male	Shoe Size Female	Size
79-81722	2 - 4	3 - 5	X-Small
79-81723	4 - 6	5 - 7	Small
79-81725	6 - 8	7 - 9	Medium
79-81727	8 - 10	9 - 11	Large
79-81728	>10	11 - 13	X-Large

Contact your DJO sales representative to learn more about the **Remedy Pro** or call 800-793-6065

DJOGlobal.com/procaresoftgoods



T 800.793.6065 D 760.727.1280 F 800.580.9625
1430 Decision Street | Vista, CA 92081-8553 | U.S.A.
DJOGlobal.com

