




SCHMITZ

A SENSE OF CONFIDENCE

arco-matic® URO Urological examination chairs

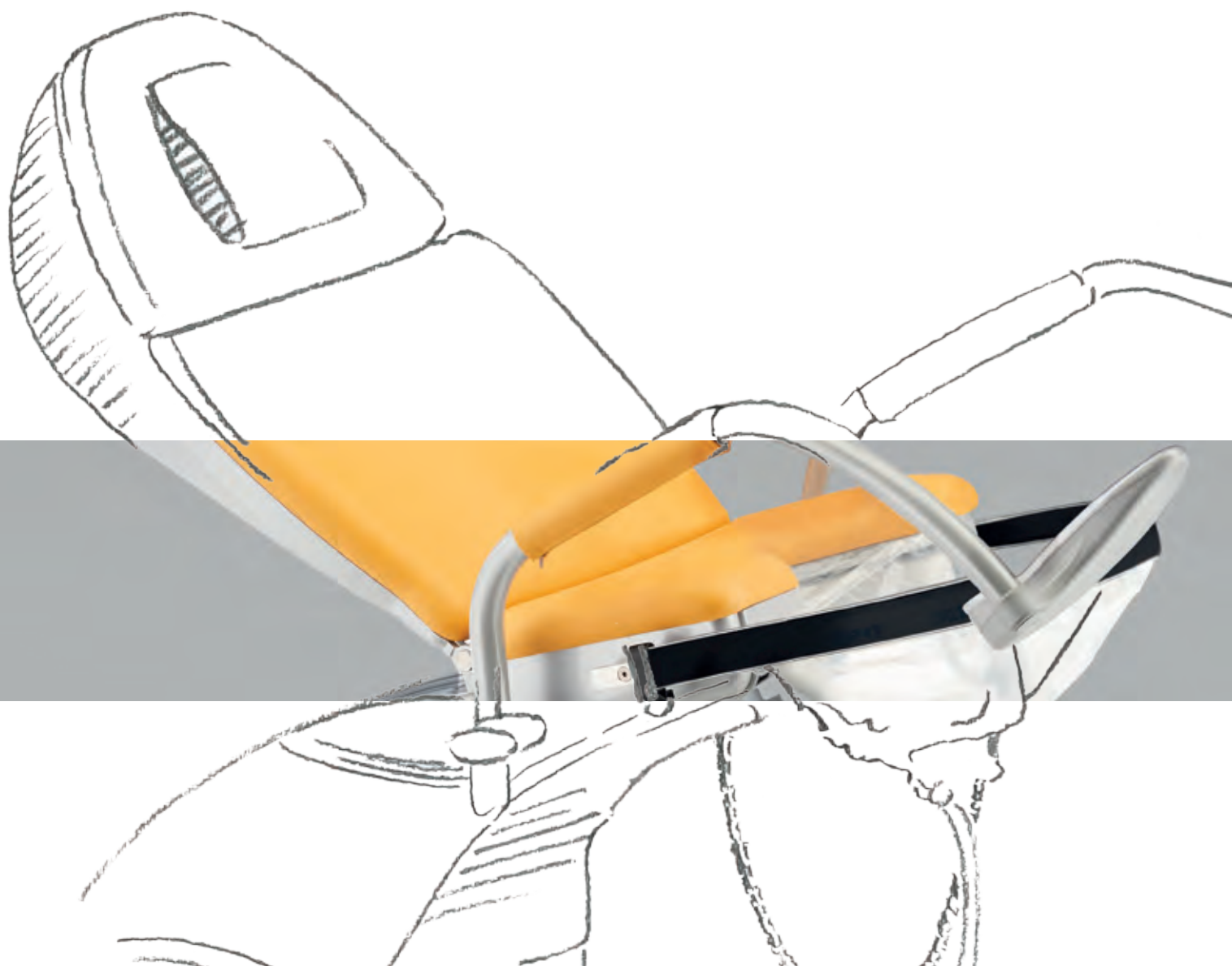

arco-matic® URO

Optimized use, high level of comfort, first-class design

arco-matic® URO, the series of urological examination chairs, combines excellent design with state-of-the-art technology and great value for money.

Ergonomic

A width of upholstery of up to 27" (685mm) at the back section and an individually adjustable head cushion attached magnetically to the chair offer patients pleasant and comfortable seating. With an access height of only 22.1" (561 mm), even older or physically challenged patients can easily take a seat on the chair, and even at its lowest level, the seat position remains horizontal.

Strong

Designed for a working load of 550 lb (250 kg), the chair offers great safety and many treatment options even for adipose patients.

Time-saving

Optimized adjustment speed and four freely selectable memory positions allow quick adjustment to the relevant position.

Functional

The backrest can be moved to a horizontal position; this means that patients can be positioned flat, e.g. for an ultrasound examination, or in a shock position in case of an emergency. The deep gyn. cut-out of 6" (155 mm) offers improved access to the examination area.

Mobile

Concealed castors make the examination chair a mobile unit, enabling flexible use of the treatment area and room layout*. The central locking mechanism allows fast and secure positioning.

Hygienic

Level, enclosed areas and the absence of corners and edges make for improved hygienic conditions and easy cleaning. Even the area under the arco-matic® can be cleaned thanks to its mobility.

Made in Germany





With its versatile, colorful options and an award-winning design, arco-matic® URO leaves nothing to be desired and adds a visual highlight to the examination room.

Design award:



**reddot award 2017
winner**

Illustration showing:
114.9702.0 with foot supports 101.0505.0

arco-matic® URO colors:
Upholstery gold yellow 86
Plastic elements silver gray

*Note: Must not be used for patient transport!

arco-matic® URO

Standard accessories

1

arco-matic® URO 114.9702.0

Examination and treatment chair for urology, with memory function for four positions, pelvic elevation, back section and foot supports (please choose!) adjustable by electromotor, seat height adjustable from 22.1" to 39.8" (561 mm to 1,011 mm) by pelvic elevation. Adjustment and memory functions activated by hand and foot control.

Including cushion for head section, one pair of side rails for seat section, LED Ambient illumination, 2 power sockets, an extractable leg plate and a fluid collection system (AllenAmatech) incl. 10 drapes. Mobile by means of 2 swivel-type double castors, Ø 2.9" (75 mm), 4 fixed-type castors Ø 2" (50 mm) mm, electrically conductive, with central lock.

Attention!

Please choose your favorite leg positioning system, see below!

You may choose between two different types of leg positioning:

2

Foot supports (pair) 101.0505.0

with integrated handles, powder coated RAL 9006 white aluminum, can be swiveled to the side, lockable. Plastic covers of transparent material as protection for the foot supports.

3

Leg supports, Göpel type, with integrated handle (pair) 101.0506.0

powder coated RAL 9006 white aluminum. Leg supports with integral foam upholstery in silver gray, including fastening strap, rod of stainless steel 18/10.



114.9702.0




101.0505.0



101.0506.0

When ordering:

Please choose your favourite color for upholstery, see color chart on page 7!

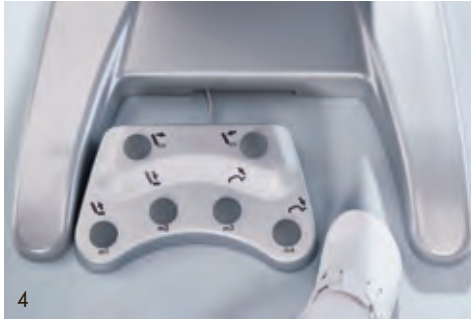

4



4 Hand and foot control

- Flexible positioning of the detachable foot control for a stress-free and comfortable working environment
- Saves and retrieves up to 4 memory positions
- Soft start-up of the chair movements
- Control of all motion sequences through foot and hand controls



5

5 Leg plate

integrated into the seat frame, extractable, padded surface of soft-touch plastics, gray. Dimensions of pad: 18.8" x 14.6" x 1.2" (480 x 373 x 30 mm)

6 LED Ambient illumination

LED illumination on both sides, tone-in-tone color shade adapted to the selected pad color (white LED for gray pads).

7 Side rail (pair)

for the seat section, stainless steel 18/10, 1.1" x 0.38" (28,6 mm x 9,54 mm), 5.9" (150 mm) long.

8 Power sockets

9 Mobility

2 swivel-type double castors, Ø 2.9" (75 mm), 4 fixed-type castors Ø 2" (50 mm), electrically conductive, with central lock. Note: Must not be used for patient transport!



6



7



8



9

arco-matic® URO for urology

Standard equipment:

Ergonomic upholstery on seat, back rest and head rest. Magazine for holding paper roll of 15.7" to 19.7" (400 to 500 mm) width, max. diam. 5.5" (140 mm), possibility to cover either the back and seat section or only the seat section with paper. With foot control and silver grey plastic cladding.

Upholstery according to color chart. Foot support mounts for quick removal or individual swiveling on a horizontal level of foot supports or Göpel-type leg supports with integrated handles.

Including cushion for head section, one pair of side rails for seat section, LED Ambient illumination, 2 power sockets, an extractable leg plate and a fluid collection system (AllenAmatech) incl. 10 drapes. Mobile by means of 2 swivel-type double castors, Ø 2.9" (75 mm), 4 fixed-type castors Ø 2" (50 mm), electrically conductive, with central lock.

Connection cable 118.1" (3 m) with mains plug. Connection: 100 - 240 V, 50/60 Hz, max. 0.94 kW.

Net weight approx. 220 lb (100 kg)
(Model no. 114.9702.0)

Safe working load: 550 lb (250 kg)


Powder coated finish:

High quality powder coating in RAL 9006 white aluminum. Gloss finish, extremely knock and scratch resistant. Easy to clean, resistant to disinfectants and environmental conditions.



Pelvic elevation

Positionings, Dimensions:

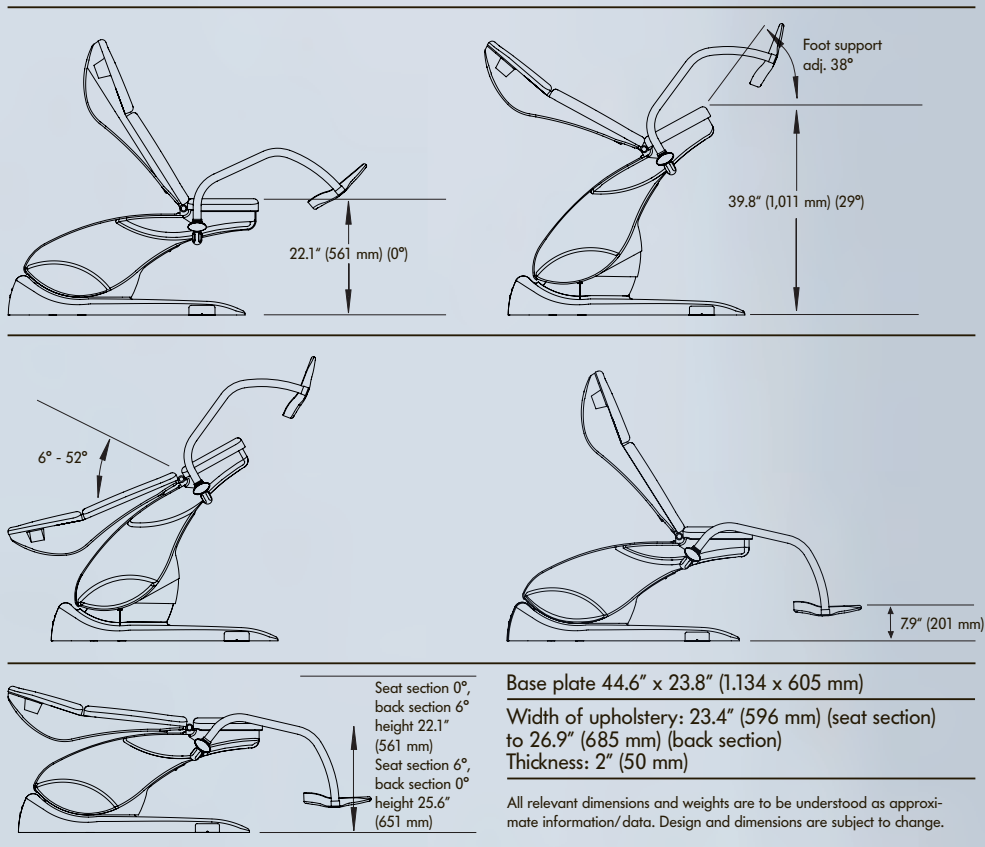

Upholstery:

Foam material on support plate, with leatherette coverings. Colors available as per SCHMITZ color chart in 6 colors.



Back section adjustment

When ordering:

Please state the following when ordering a chair from the arco-matic® URO series: Upholstery color see color chart on page 7.



Electromotive adjustment of foot supports

Colors for upholsteries and plastic elements

Leatherette



gold yellow 86 101.0532.0




red carmine 92 101.0533.0




caribbean blue 82 101.0536.0



lime 99 101.0535.0

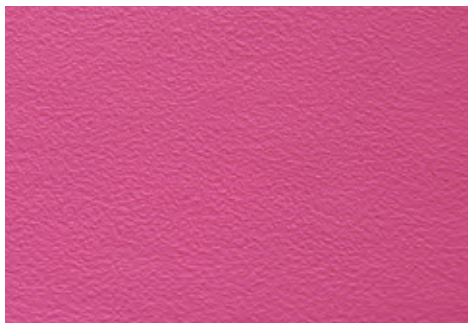

raspberry 103 101.0534.0

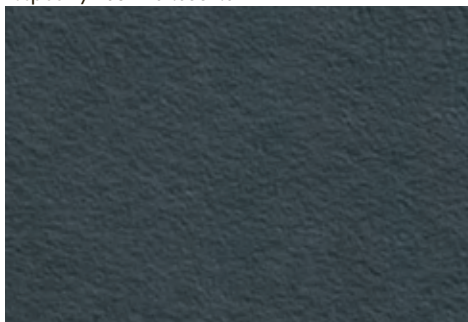

titan grey 91 101.0537.0

Plastic elements



silver-gray

Please choose your favorite color for upholstery!

Each set of pads in the color you choose consists of:

- Pads for head, back and seat section
- One cushion for the head section, with magnetic attachment
- One pair of leatherette covers for the leg positioning system



Legal manufacturer:

**SCHMITZ u. Söhne
GmbH & Co. KG**
P. O. Box 14 61
58734 Wickede (Ruhr)
Germany (postal address)
Zum Ostenfeld 29
58739 Wickede (Ruhr)
Germany (visiting address)
T +49 (0)2377 84 0
F +49 (0)2377 84 135
www.schmitz-soehne.com
export@schmitz-soehne.de

Technical service, hotline:
T +49 (0)2377 84 549
F +49 (0)2377 84 210
service@schmitz-soehne.de

The products may only be used
for their intended purpose of use.
Please also observe our operating
manual!

United States of America:

SCHMITZ USA LP
5126 South Royal Atlanta Drive
Tucker GA 30084-3052
U.S.A.
T +1 770-349-6308
F +1 770-934-3384
info@schmitzusa.com

**Information on
our establishments
worldwide can be
found on our website.**



 /schmitz.soehne

