

StatSpin Express 2 Centrifuge



High Productivity. Lower Cost.

Eliminate sample preparation bottlenecks in busy clinical laboratories with the StatSpin Express 2 Centrifuge. This compact and quiet centrifuge spins standard 5 or 3 mL blood collection tubes along with inserts for 1.5-2mL microtubes at three factory-set timed cycles 30, 120, and 180 seconds.

For peak work flow, several units can provide higher throughput than larger, more costly alternatives. Run up to 4 tubes per cycle with minimal risk of tube breakage due to the unique rotor design.

Key Product Features:

- Compact low-profile design
- No need to carry the rotor back and forth to the centrifuge.

Range of Applications

- Variety of selectable run times to precisely replicate most procedures developed on other centrifuge equipment. Run up to 4 tubes at a time.




StatSpin Express 2 Centrifuge Specifications¹

Dimensions:	Height 6 in (15.2 cm) Width 6.6 in (16.8 cm) Depth 8.6 in (21.9 cm) Weight 5 lbs (2.3 kg)
Cycle Times:	Quick Spin: 30 seconds Normal: 120 seconds Hard Spin: 180 seconds
Speeds (Forces):	8500 rpm (\pm 5%) / 4440 xg
Acceleration Time:	Approximately 20 seconds
Deceleration Time:	Approximately 25 seconds
Tube Accommodation:	Accepts up to 5 mL collection tubes (13 x 75 mm), Monovette syringe tubes and smaller tube formats. Adapters available for 1.5 – 5 mL tubes.
Power Requirements:	100-240V, 50/60 Hz, 1.7 amp (Remote power supply has a 50 cm/3.5 ft line to the centrifuge and a 5 ft North American power cord. Replacement AC power cord must have an IEC 320 connector and meet <HAR> standard).
Warranty:	One year

Ordering Information

SSX2	StatSpin Express 2 Centrifuge for 100-240 V, 50/60 Hz. Includes 4-place rotor (RTX4), universal power supply, removable chamber liner, four 5 mL inserts and North American power cord.
------	---

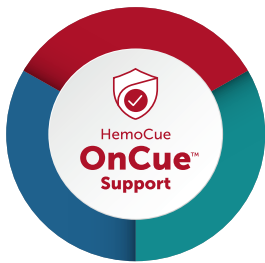

HemoCue OnCue™ Support

HemoCue OnCue™ Support is a comprehensive program designed to meet your needs as quickly and conveniently as possible.

Powered by HemoCue employees, you get answers and education from the experts who know our products best.

At your office, online, or on the phone, we're ready to do more for you.

¹StatSpin Express 2 Centrifuge Instructions For Use 55-003910-001FD

Specifications subject to change without prior notice. The StatSpin Express drive and suspension system are covered under U.S. Patent No. 5,566,919.

Certified by a Nationally Recognized Testing Laboratory (NRTL) for the US and Canada.

©2020 Beckman Coulter. All rights reserved.

StatSpin products are developed and manufactured by Beckman Coulter, Inc. and sold through HemoCue America, a division of Radiometer America Inc. Beckman Coulter, the stylized logo, and the Beckman Coulter product and service marks mentioned herein are trademarks or registered trademarks of Beckman Coulter, Inc. in the United States and other countries.

HemoCue America | California, USA | Phone 800.881.1611 | Fax 800.333.7043
web@hemocue.com hemocue.us

©2020 HemoCue America. All rights reserved. MMUS-0053[3]

